
SILVERGEAR

MASSAGE GUN PRO
RELAX AND RECOVER

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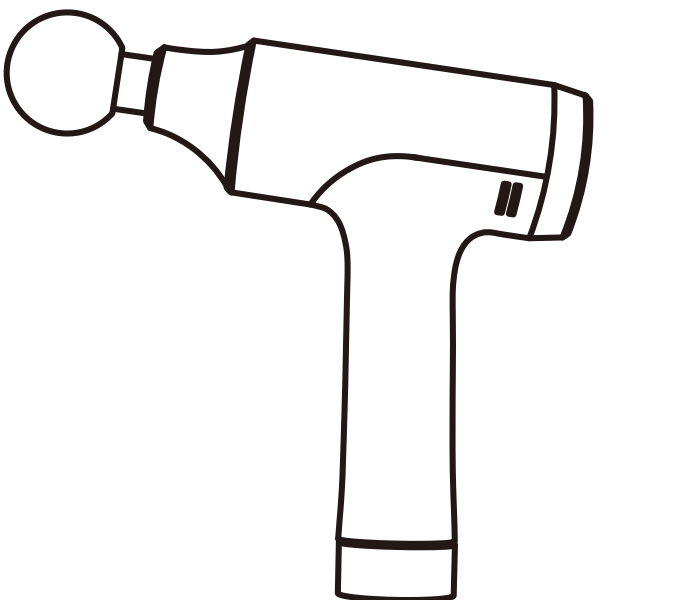160.00 mm

160.00 mm

EN
Thank you for purchasing our Silvergear Massage Gun Pro. This manual will provide instructions for effective and safe use of the product. Please read it carefully before use and store the manual in a safe place for future reference.
If you ever pass the product on to others, this manual should be included.

1. Package content
-Massage gun
-4 Massage heads
-Charging cable
-Manual

2. Technical data
Input: 5V 1A – 2A
Battery: Li-ion, 3.7V 2400 mAh
Rated Voltage: 100-240V, 50-60Hz
Charging time: approx. 2-3 hours
Operating time: approx. 3 hours
Dimensions: approx. 18 x 7 x 23 cm
Weight : approx. 907 g

3. Description of parts
1.LED display
2.Power button
3.Charging port
4.Massage head

Massage heads
A.Ball head for large muscle groups
B.Bullet head for body joints
C.Flat head for all body parts
D.U head for the vertebrae

4. Safety instructions
•For intended use only, as described in this manual.
•Do not let moisture come into contact with any electrical parts.
•Never touch the device or its plug with wet hands.
•Before use, inspect the device and its parts. Do not use in case of damage.
•Before use, consult your specialist in case of a medical condition or treatment.
•Only use on soft tissues of the body and directly on clean, dry skin. Do not use through clothes.
•If you experience any discomfort or pain, stop use immediately and contact your specialist.
•Not suitable for diabetics or people with poor blood circulation.
•Do not use on sensitive, irritated or damaged skin.
•Some bruising may occur regardless of control setting or applied pressure. Frequently check the treated area during use.
•Always unplug before installing, removing or adjusting any of the accessories and for long-term storage.

5. Operation
Charging
1.Before first use, fully charge the battery. This will take about 2-3 hours.

6. Maintenance
If necessary, clean the device with a slightly damp towel and wipe it dry with a soft cloth. Store it in the supplied storage case and keep it a dry, cool place, away from direct sunlight.

7. Disposal
Dispose of the article and the packaging materials in accordance with current local regulations.

8. Environmental Protection
This product is subject to the provisions of the European Directive 2012/19/EU. The crossed out wheeled bin symbol means that this product shall not be disposed of with normal household waste. Electronic and Electrical Equipment not included in the selective sorting process are potentially dangerous for the environment and human health due to the presence of hazardous substances. Please dispose responsibly at an approved waste or recycling facility.

9. CE Marking
The CE marking is a certification mark that indicates conformity with health, safety and environmental protection standards for products sold within the European Union. The CE marking is the manufacturer's declaration that the product meets the requirements of the applicable EC directives.

10. Recycling
Recycling helps to reduce the consumption of resources and reduce the environmental impact. If you wish to discard the packaging, comply to the environmental legislation in your country.

NL
2.Plug the provided USB charging cable into the massage gun and connect it to an adapter (not provided).
3.The LED display flashes during charging and shows 100% when the massage gun is fully charged.

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