
Rechargeable Sonic Power Toothbrush
Ultrasoon Elektrische Tandborstel
Operating Instructions
Handleiding
Model: DT20

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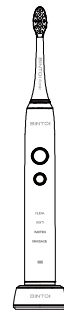ENG

IMPORTANT SAFEGUARDS

The Power toothbrush is safe and effective when used according to directions. Use all electrical appliances, however, basic safety precautions should be followed, especially when children are present.

DANGERS-WATER

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not place or drop the charger base where it can be pulled into or spilled with water.
- Do not use the product if the charger base is wet or after liquid spillage.
- Do not use the product if the charger base is wet or after liquid spillage.

DANGERS-VOLTAGE

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug the charger base into an outlet with a voltage other than specified on the label of the charger.
- The product is designed to operate within a range of 100-240 volts.
- Always disconnect the cord from the power source.
- Do not use the plug into an outlet.
- The plug does not fit easily into the outlet.

DANGERS-WARNINGS

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not use the product if the charger base is wet or after liquid spillage.
- Do not use the product if the charger base is wet or after liquid spillage.

SAVE THESE INSTRUCTIONS

CONFIGURATION

Professional brush head

Color ring

Power button

Mode indicator

Charging indicator

Charger base

Power cord

Power plug

HOW TO USE

Slide the brush head onto the handle. Push brush head down until it locks into place on handle. There is a small gap between color ring and the handle (Fig. 1).

Use the toothbrush and apply a small amount of toothpaste.

Brush the teeth and gums (Fig. 2).

Push the mode button to select mode.

If the toothbrush has not been used for a long time, it may need to be recharged.

CHARGING

1. **SOFT:** Gentle cleaning with polishing sensitive teeth.

2. **WHTEN:** Teeth whitening with polishing sensitive teeth.

3. **MASSAGE:** Pulse with gentle gum stimulation.

GETTING STARTED

Charging the toothbrush

Be sure to turn off the toothbrush before charging. Plug the power plug into the outlet. Then insert the charger. The charger indicator light will be on, indicating the toothbrush is charging. After the battery is fully charged, the charger indicator light will turn off, indicating the toothbrush is fully charged.

Notes:

- It takes at least 24 hours to charge the battery.
- The toothbrush must be charged for at least 24 hours the first time using the toothbrush.
- The toothbrush has not been used for a long time, it may need to be recharged.
- After the battery is fully charged, the charger indicator light will turn off, indicating the toothbrush is fully charged.

HOW TO USE

Slide the brush head onto the handle. Push brush head down until it locks into place on handle. There is a small gap between color ring and the handle (Fig. 1).

Use the toothbrush and apply a small amount of toothpaste.

Brush the teeth and gums (Fig. 2).

Push the mode button to select mode.

If the toothbrush has not been used for a long time, it may need to be recharged.

CHARGING

1. **SOFT:** Gentle cleaning with polishing sensitive teeth.

2. **WHTEN:** Teeth whitening with polishing sensitive teeth.

3. **MASSAGE:** Pulse with gentle gum stimulation.

TROUBLESHOOTING

Problem	Possible cause	Solution
The toothbrush does not work.	You have not purchased the correct model for your country.	Change the main unit for the correct model.
The toothbrush does not work.	The toothbrush is not fully charged.	Charge the main unit for at least 24 hours.
The toothbrush does not work.	The toothbrush is not fully charged.	Charge the main unit for at least 24 hours.

WARRANTY

Your toothbrush is guaranteed for 1 year against manufacturer's defects. The warranty is not valid if the toothbrush is damaged due to misuse, accident, fire, or other causes. The warranty is not valid if the toothbrush is damaged due to misuse, accident, fire, or other causes.

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MAINTENANCE

Be sure to remove the plug from the electrical outlet before cleaning the product.

Toothbrush

Remove the brush head. Wipe all surface areas with a damp cloth. Do not submerge the handle into water. There the brush head and handle into water. Do not use the brush head and handle into water.

Charger

Using a soft cloth, wipe off any dirt on the charger and the settings. Clean the plug periodically to ensure it remains in good working order. Do not use the charger if the plug is damaged. Do not use the charger if the plug is damaged.

TROUBLESHOOTING

Problem	Possible cause	Solution
The toothbrush does not work.	You have not purchased the correct model for your country.	Change the main unit for the correct model.
The toothbrush does not work.	The toothbrush is not fully charged.	Charge the main unit for at least 24 hours.
The toothbrush does not work.	The toothbrush is not fully charged.	Charge the main unit for at least 24 hours.

WARRANTY

Your toothbrush is guaranteed for 1 year against manufacturer's defects. The warranty is not valid if the toothbrush is damaged due to misuse, accident, fire, or other causes. The warranty is not valid if the toothbrush is damaged due to misuse, accident, fire, or other causes.

DISPOSAL

The entire surface of this product is designed to be recycled. Do not dispose of this product in the trash. Remove the battery before you dispose of the product. If you have recycling programs, contact your local waste management office.

Removing and recycling the battery

Remove the battery from the product. Do not dispose of the battery in the trash. Remove the battery from the product. Do not dispose of the battery in the trash.

BINTOI

Professional brush head

Color ring

Power button

Mode indicator

Charging indicator

Charger base

Power cord

Power plug

HOW TO USE

Slide the brush head onto the handle. Push brush head down until it locks into place on handle. There is a small gap between color ring and the handle (Fig. 1).

Use the toothbrush and apply a small amount of toothpaste.

Brush the teeth and gums (Fig. 2).

Push the mode button to select mode.

If the toothbrush has not been used for a long time, it may need to be recharged.

CHARGING

1. **SOFT:** Gentle cleaning with polishing sensitive teeth.

2. **WHTEN:** Teeth whitening with polishing sensitive teeth.

3. **MASSAGE:** Pulse with gentle gum stimulation.