# M A I C U R A **PRO**

HANDLEIDING | USER GUIDE

WELKOM

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
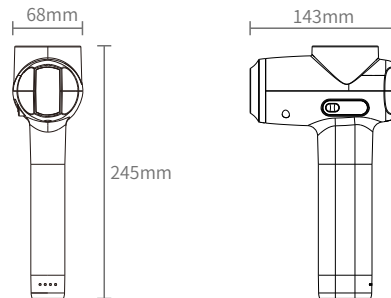
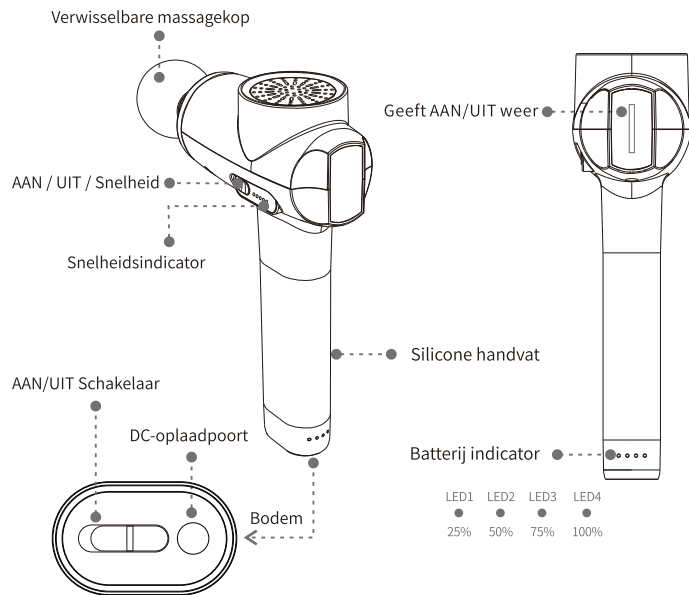
WELCOME

# Congratulations on your new Maicura

You are one step closer to moving better, feeling better, having less pain, and simply being a better version of yourself.

The Maicura massage gun was created to keep you doing what you love. Whether you're an athlete, a hard worker, want to relax, or simply trying to get your body to feel its best: the Maicura massage gun will help you get there.
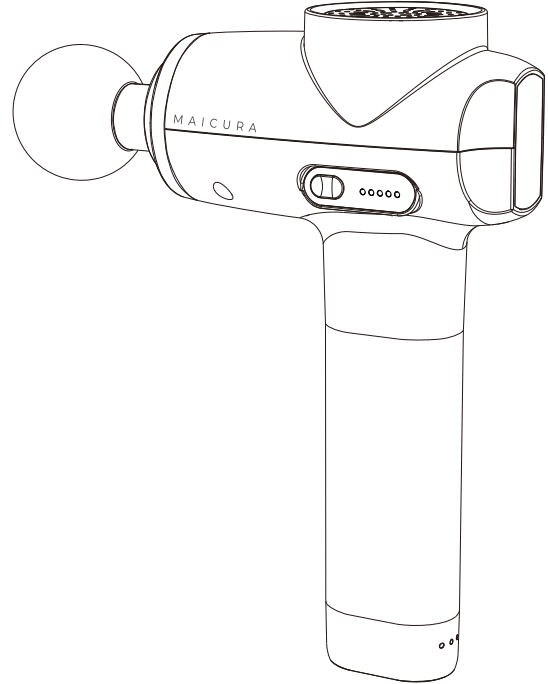
## Lets get started

**1. Power up**

Switch your massage gun to the ON position at the base of the handle. Once your device is powered ON multiple blue lights will illuminate.

**2. The preferred speed**

Choose the speed you want (1-5) by pressing the speed setting button on the side of the massage gun once per speed. The blue lights on the side of the massage gun will indicate the current speed.

**3. Changing attachments**

Choose the massage head you want and insert it into the shaft/opening while pushing firmly. Are you using the Fork attachment? Align the protrusion on the fork with the groove on the shaft/opening and push. Are you using the infrared head? It turns on automatically when the head is in the massage gun and you turn on the first speed.

**4. Recharging**

Only use the included Maicura charger to maintain the warranty. When the battery of the massage gun is lower than 20%, the left LED light begins to flash. Full charge is indicated when the 4 LED lights on the handle stop flashing and turn persistently blue. When the battery of the infrared head is lower than 20%, one of the LED lights starts flashing. when the center LED light turns persistently red.

**5. Taking care of your massage gun**

Use a damp, clean, cloth and gently wipe your massage gun. Make sure the battery charger is not attached and that the massage gun is turned off. Store in a clean, cool place when not in use.

## HEAD ATTACHMENTS

### Choosing the best massage head

#### **Ball**

Suitable for all muscle groups.

#### **Cushion**

Suitable for sensitive areas such as wrists, ankles, knees, and around the neck.

#### **Flat**

Suitable for large surface areas such as hamstrings, quads, back, and chest.

#### **Bullet**

Suitable for deep tissue and trigger points.

#### **Fork**

Suitable for longer muscles such as calves, forearms, biceps, and ankle/Achilles area.

#### **Infrared**

Suitable for all muscle groups. Helps with improving the circulation of blood, the relaxation of stiff and sore muscles, and relief from pain. Note: It takes a little while for the massage head to warm up.

## INSTRUCTIONS

WARNING: TO REDUCE THE RISK OF ELECTRIC SHOCK, FIRE, PERSONAL INJURY, OR PROPERTY DAMAGE, THIS PRODUCT MUST BE USED BY THE FOLLOWING WARNINGS, CAUTIONS, AND SAFETY INSTRUCTIONS:

### **IMPORTANT SAFETY INSTRUCTIONS**

The massage gun should be used by applying light pressure and moving it over the body for approximately (60 - 120) seconds per muscle group. The massage gun should only be used on soft tissue areas of the body and should not be used on any hard areas, such as bones, including the head. Please immediately stop using the massage gun upon experiencing any pain or discomfort, apart from light muscle soreness. When using the Maicura massage gun, regardless of the pressure exerted, there is a risk of bruising. If this is the case, the massage gun should not be used on the affected area until it has completely healed. This also applies to other bruises, contusions, rashes, and, other irritated or injured areas of the skin until they are completely healed.

The massage gun must be kept away from liquids and heat sources at all times. It must therefore not be immersed in water and be stored in a dry place. The ventilation holes on the top of the massage gun should be kept free of dust, sand, dirt, food, and/or other substances that could enter and clog the holes. In addition, the massage gun should never be used for more than 1 hour in a row and the massage gun should rest for at least 30 minutes after each use before switching on again.

THE MASSAGE GUN SHOULD ONLY BE CHARGED WITH THE CHARGER THAT WAS SUPPLIED WITH THE MASSAGE GUN.



After charging the massage gun, the charger should be unplugged from the massage gun before use. Before disposing of the massage gun the battery must be removed. The battery must be complied with safely at an appropriate recycling facility.

The massage gun should never be used for any purpose other than its intended use, should not be dropped, and should never be tampered with or altered in any way.

PLEASE DO NOT USE THE MASSAGE GUN, OR ANY PERCUSSION OR VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM THE PHYSICIAN IF ANY OF THE FOLLOWING APPLY:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery or injury, epilepsy or migraines, herniated disks, recent joint replacements or IUDs, metal pins or plates, or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any percussion or vibration device. These contraindications do not mean that you are not able to use a percussion or vibration device but we advise that before doing so you consult your physician or doctor.

## SPECIFICATIONS

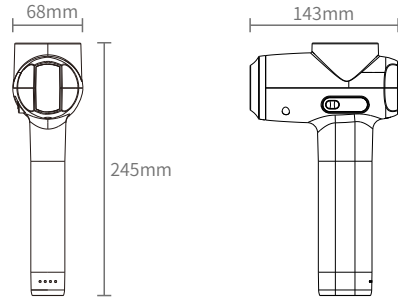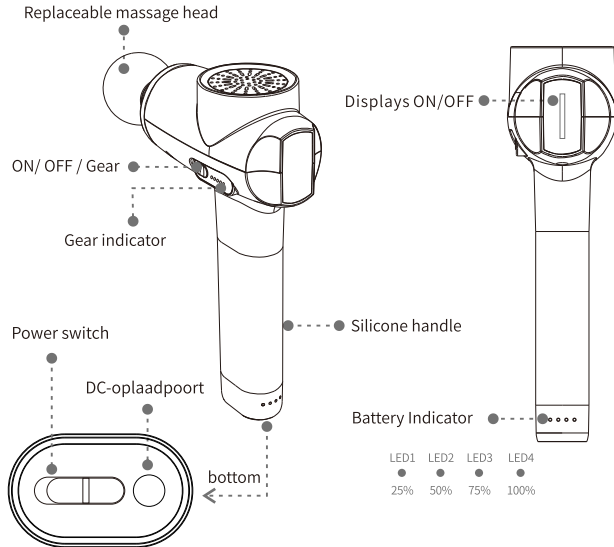Charger: Input: 100-240V 50/60Hz 0.4A Output: 21.9V 0.8A  
 Frequency: 1800-3200 R/MIN  
 Battery: 2600 mAh  
 Weight: 900g  
 Brand: Maicura  
 Address: Tourslaan 13, 5627KW, Eindhoven (NL)