
**Alecto BC-19**

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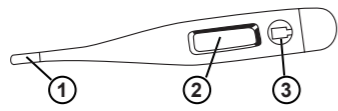
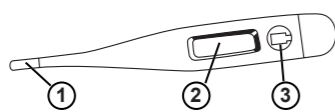
- 1: temperature sensor
- 2: display
- 3: on/start/stop/memory/off button

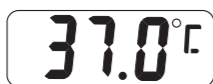temperature display

temperature display memory

battery empty

**INTRODUCTION**

This Alecto BC-19 thermometer is a digital thermometer developed to take your temperature as simply and accurately as possible. In order to get the most reliable results, these instructions must be followed accurately. For the most accurate temperature reading, we recommend using the thermometer rectally. Although the generally accepted 'normal' temperature is 37.0°C, body temperature can vary from 36.1°C to 37.2°C. Fluctuations in temperature can be due to activities such as sports, smoking, eating and drinking, and also depend

on the part of the body where the temperature is taken. When measured orally, the body temperature is generally 0.5°C lower than when measured rectally, while a temperature measured under the armpit is 1.0°C lower. Your temperature is lower in the morning than in the evening.

**SPECIFICATIONS**


Temperature range: 32.0°C - 42.9°C  
 Low temperature display: temperature ≤ 32°C, display: Lo°C  
 High temperature display: temperature ≥ 42.9°C, display: HI°C  
 Display resolution: 0.1°C  
 Accuracy: 35.5°C - 42.9°C: ± 0.1°C  
 ≤ 35.5°C ≥ 42°C: ± 0.2°C  
 Normal temperature ≤ 37.8°C: 10 beeps  
 Fever alarm ≥ 37.8°C: 10x3 beeps  
 Auto-OFF: 10 minutes ± 2 minutes  
 Battery: LR41, 192 or SR41  
 Power consumption: 0,15 mW  
 Battery life: 2 years  
 Dimensions: 136 x 22 x 14,5mm  
 Weight: ±14 gr (batt incl.)  
 Conditions for use:  
 Temperature: 10°C - 35°C  
 Humidity: 25% - 80% RH  
 Storage conditions:  
 Temperature: -25°C - 55°C  
 Humidity: 25% - 80% RH

**MEASURING YOUR TEMPERATURE**

1. Press the ON/OFF key once to activate the thermometer. You will hear a beep.
2. The display shows in succession: '188,8°C' the last measured temperature '37,0' (for internal calibration) 'LO' and flashing '°C' (\*)

3. Place the thermometer in the chosen body part to measure your temperature.
  4. When you hear beeps, your temperature has been recorded.
    - if your body temperature is normal, you will hear 10 beeps
    - if your body temperature is above normal (≥37.8°C), you will hear 10x3 beeps.
  5. Press the ON/OFF key to switch off the thermometer. (The thermometer will switch off automatically after approx. 10 minutes).
- \*: If the room temperature is higher than 32.0°C, the room temperature will be shown instead of 'Lo'.

**Measurement under the tongue:**  
 Place the mometer under the tongue, as far back as possible to the left or right.  
 A temperature taken under the tongue is approx. 0.5°C lower than when taken rectally.

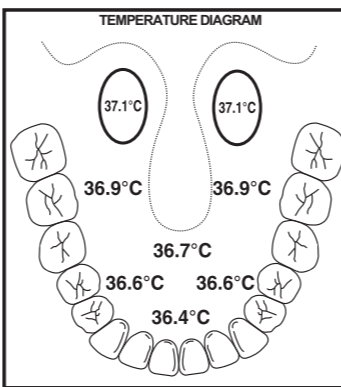**Measurement under the arm**  
 Make sure the armpit is dry and place the thermometer under the arm. The sensor must touch the skin. Keep your arm crossed in front of your chest. A temperature taken under the tongue is approx. 1.0°C lower than when taken rectally.

**Rectal measurement**  
 Insert the thermometer sensor gently (maximum 1.25 cm).

**BATTERY LOW**

The battery is low when the symbol appears in the display. The battery cannot be replaced. Please replace the complete thermometer.

**DISPOSAL OF THE DEVICE**

At the end of the product lifecycle, you should not throw this product into the normal household garbage but bring the product to a collection point for the recycling of electrical and electronic equipments.

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**ADVICE AND WARNINGS**

- Keep the thermometer and battery out of the reach of children.
- Do not talk or walk around while taking your temperature.
- Store the thermometer in the special case.
- Do not expose the thermometer to direct sunlight, dust or moisture.
- Do not dismantle the thermometer.

**STORAGE AND CLEANING**

Clean the thermometer before and after use, using a soft cloth and alcohol, or cold soapy water. The thermometer may be briefly immersed in water or alcohol (do not soak or sterilise). Store the thermometer in the special case.

**CE**

The declaration of conformity is available at:  
<http://DOC.hesdo.com/BC-19-DOC.pdf>

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Hangzhou Hua'an Medical & Health Instruments Co.,Ltd,  
 Building 2, 1# Fuzhu Nan RD,  
 Wuchang Town, Yuhang District,  
 Hangzhou, Zhejiang 310023  
 Model: DT-K101A

Shanghai International Holding Corp. GmbH,  
 Eiffestrasse 80,  
 20537  
 Hamburg Germany

Imported by: Hesdo bv.  
 Azielaan 12, 5232BA  
 's-Hertogenbosch  
 The Netherlands

**Alecto**  
 Service ↓ Help  
  
[WWW.ALECTO.NL](http://WWW.ALECTO.NL)  
[SERVICE@ALECTO.NL](mailto:SERVICE@ALECTO.NL)  
 Hesdo, Australiëlaan 1  
 5232 BB, 's-Hertogenbosch  
 The Netherlands

Hangzhou Hua'an Medical & Health Instruments Co.,Ltd,  
 Building 2, 1# Fuzhu Nan RD,  
 Wuchang Town, Yuhang District,  
 Hangzhou, Zhejiang 310023  
 Model: DT-K101A

Shanghai International Holding Corp. GmbH,  
 Eiffestrasse 80,  
 20537  
 Hamburg Germany

Imported by: Hesdo bv.  
 Azielaan 12, 5232BA  
 's-Hertogenbosch  
 The Netherlands