
USER MANUAL

HANDLEIDING ZAKLAMP

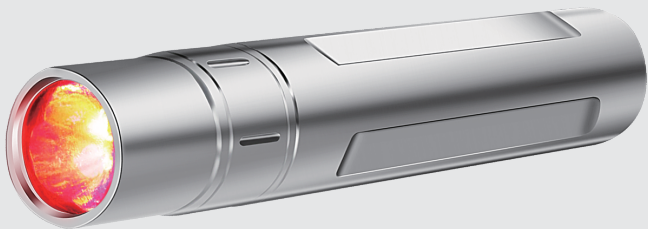
INDEX

1- 6.....English

8- 13.....Nederlands

INSTRUCTION MANUAL

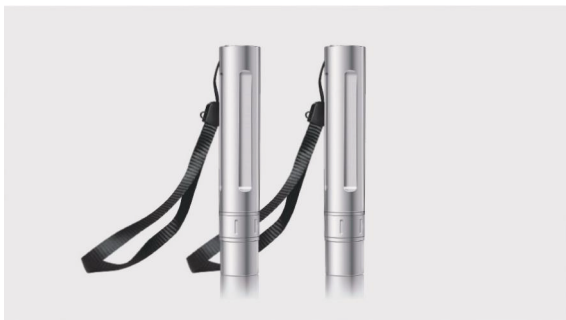Congratulations on taking control of your pain ! Now you can use the LED Red Light Therapy to provide direct, non-addictive, drug free pain relief.

Efficient Cells Can Enhance Health

Please read the instructional manual to get the most out of your device.

Uses a combo of red 630nm 660nm and Infrared 850nm LEDs which is chosen to give the ideal light intensity for biological effects anywhere on the body.

1. INCLUDED WITH THIS PACKAGE

- 1* Red Light Therapy Torch
- 1* USB Cable
- 1* User Manual

2. BENEFITS

Temporary relief of minor muscle and joints pain , arthritis and muscle spasm, relieving stiffness, promoting relaxation of muscle tissue.

3. DIRECTIONS

1. Commit to a 3-week Treatment Plan. While many people achieve pain relief from Red Light Therapy almost immediately, others see measurable pain relief incrementally over time. Depending upon injury type severity and different pain pathways involved, people vary in therapy response time.

This is why we recommend committing to a 3-week Treatment Plan. If you see results right away, great, continue treating as needed. However, if you do not feel relief as quickly, know that this is also normal, so try not to get frustrated. Simply treat your pain consistently with 3 quick treatments per day with the, For up to 3 minutes each time per treatment spot. After 3 weeks of consistent red light therapy, most people experience a diminishing of their pain.

2. Do not suddenly stop your other treatments upon beginning red light therapy. Although you may not be satisfied with current treatments, they may be having a partial effect on the level of your pain. While continuing with your current treatments, start your red light therapy, giving your red device time to provide optimal pain relief. After you have achieved maximum pain relief from your red device, you may then consider slowly weaning yourself off of other medications or treatments. If you plan to reduce or stop a prescription medication, first consult with your physician.

For many people, Works so well that they are able to come off pain pills completely. Others are able to reduce the amount of pain medications they are on, reducing the risk of drug side effects.

Hold the unit over the painful area at a distance anywhere from touching the skin to 1 inch away from the skin , Depending on the depth of treatment desired . The Red light penetrates up to 1 inch into your tissues or joints.

Use the Dosage Chart below to find the best Dosage to treat your injury. Dosage Chart.

Pain Level	Light Pain	Average pain	Deep pain
DOSAGE	1 DOSE Per treatment spot (1 minute treatment) Once a day	2 DOSE Per treatment spot (2 minutes treatment) Twice or 3x a day	3 DOSE Per treatment spot (3 minutes treatment) 3 times a day

Do not keep the device on for more than 5 minutes at a time. If you want to treat several parts of your body consecutively, turn off the device for at least 1 minute, after every 5 minutes on . This will keep the unit from overheating.

The treatment spot is an area covered by the red light beam. A small joint i.e. knuckle of a finger, should have only one treatment spot, because the light shines over the entire affected joint.

For wrist treatment

For finger treatment

For knee treatment

4. PRODUCT SPECIFICATION

Product Size	140 x25mm	Net Weight	0.2 kg
Material	Stainless Steel	Wavelength	630:660:850nm
Irradiance	Surface:217mw/cm2	Irradiance	3 inch:51mw/cm2
Illuminance	180 lux at 15cm	Battery	18650
Warranty	24-Month	Lifespan	50,000 Hours

5. Caution!

- Do not use the lamp near water
- Do not stare directly into the light use protective glasses while using the device
- Do not handle near the eyes

6. Maintenance of the equipment

Cleaning and disinfection:

Use only a soft dry cloth to clean the exterior. Never use an abrasive detergent or immerse it in water or any other liquid. This device cannot be disinfected by any liquid.

Storage: Store the device in a cool, dry place, away from direct light and extreme heat or humidity.

Warning: Do not allow moisture to come into contact with the housing or light of the device. Do not immerse the device in water.

Warning: If the device is not used for an extended period of time (1 month +), remove the batteries from the device to prevent corrosion.

To request permission, please write to the publisher at "Attention: Permissions Coordinator," at the address below.

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