

**MINI** *injection*  
DIAMOND DERMABRASION

# Diamond Exfoliation Spraying Device

Designed to Energize Your Skin

## User Manual

**MINI** *injection*  
DIAMOND DERMABRASION



# Introduction

---

The Diamond Exfoliation Spraying Device integrates two functions in one device including exfoliation and spraying. It is specially designed for your skin and gives you the best skin care experience.

Dear customer, thank you for purchasing this device!

Please read this product manual through before device operation considering your safety. Keep this user manual properly for reference next time.



# Contents

---

Product functions-----	1
Product display-----	3
Use guide-----	5
Istructions for-----	6
Caution-----	7
After sale service-----	8

# Product functions

## Why does skin get dehydrated?

Most skin problems caused by dehydration



Air pollution

Season change

Unsuitable skin care products

Improper skin care

In daily life, women are more likely to try various methods to prevent skin from dehydration by using a pile of skin care products like face mask, toner or essence lotion. However, many skin problems, like spot and acne, still emerge. Why? That is because skin pores are clogged by dirty things which fail skin to absorb skin care nutrition. Now with this device, you are going to clean your skin completely, addressing skin problems like acne, wrinkles.



### Exfoliation:


The exfoliation head mounted on the exfoliation pen comes with inlaid crystallite particles which is effective in removing dead skin and promoting skin metabolism. The exfoliation is a non-invasive way of skin care method which is safe and effective in addressing skin problems like acne, large pores and fine wrinkles.

### Spraying :

The diameter of skin pore is about 0.02--0.05mm which is easy for water molecule to enter in. With the help of the spraying function, the water molecules can be divided into smaller particles with diameters 0.3--0.5um(about 0.0003--0.0005mm) to help skin absorb in skin care nutrition and balance water intake for skin.

The spraying apparatus mounted on the device is designed to “process” the liquid through it. Specifically, this spraying apparatus can vaporize the liquid through it to make the particles of the liquid to be more effectively absorbed by skin.

# PRODUCT DISPLAY

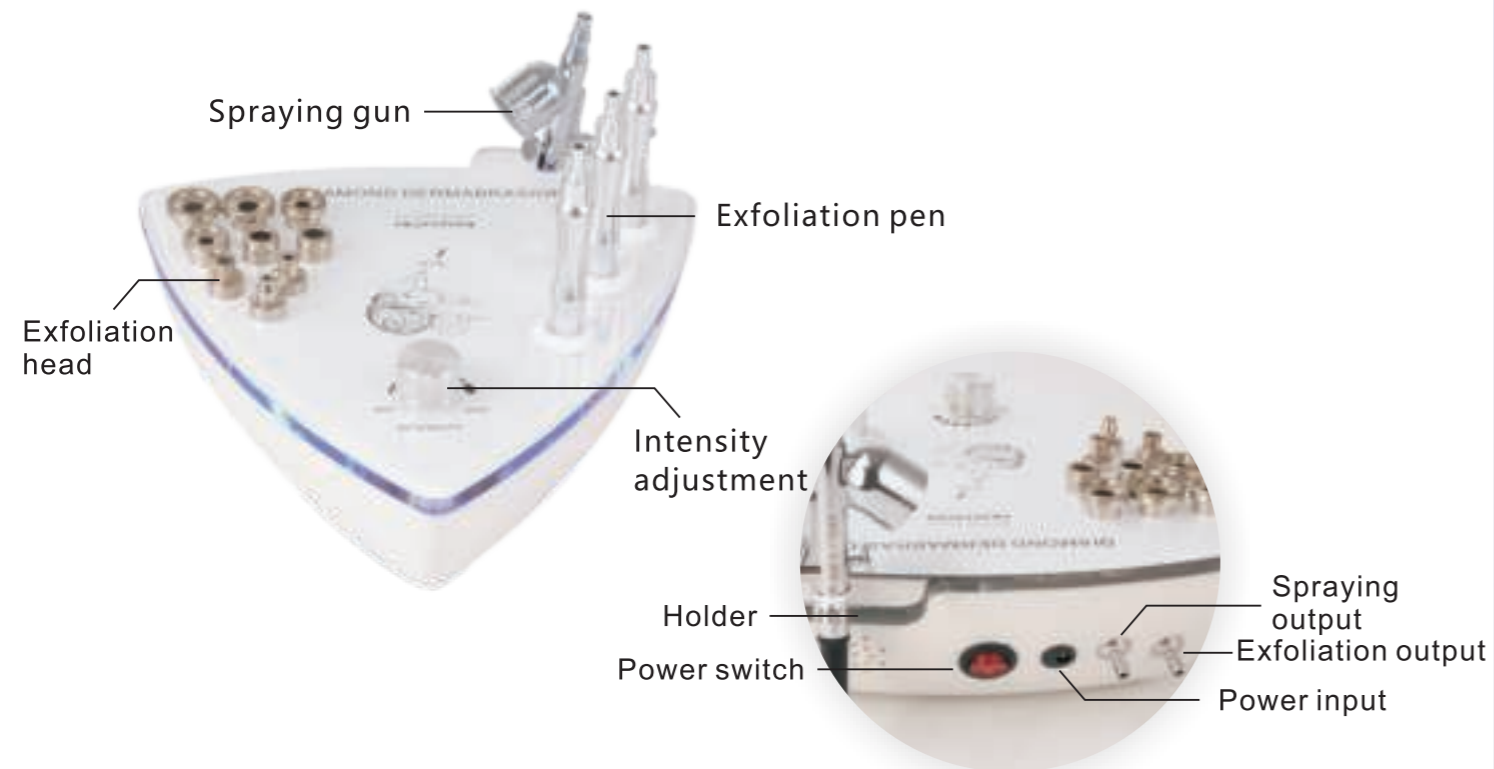

# PRODUCT PARAMETERS



Product name: Diamond Exfoliation Spraying Device

Product mode: LB197

Device size: 230\*245\*80mm

Packaging size: 335\*260\*105mm

Net weight: 2.7kg

Power input: 100-240V 0.8A 50/60Hz

Power output: 12V=4.0A

Device accessory:

1 Main device\*1

2 IPVC vacuum tube\*2

3 Exfoliation pen\*3

9 Exfoliation head\*9

1 Spraying gun\*1

1 Device cable\*1

1 Spraying gun packaging case\*1

1 Exfoliation head and sponge packaging case\*1

# USE GUIDE

1. please place this device stably in a dry working environment and meanwhile make sure the heat emission holes can ventilate heat without blocking.
2. Connect the spraying gun, exfoliation head and exfoliation pen rightly.
3. Connect the device to power and plug the device.
4. Press the red switch button on the back of the device.
5. After 30 minutes continuous working, pause the device for 20--30 minutes before second time using.
6. Adjust the device's working intensity to a proper degree according to client's condition.
7. After use, please switch off the device and unplug the device.



# INSTRUCTIONS FOR



1. Clean face first then smear some skin care gel on face and massage gently.



2. After face massaging, clean face with clean water.



3. Switch on the device and select the exfoliation function , adjust to proper intensity , to exfoliate face area about 10 minutes. After this, wear a facial mask to care and refresh skin.



4. Add essence lotion into the water chamber of the spraying gun, then switch on the spraying gun and close eyes.



## CAUTION

- 1.this device belong to sophisticated category which requires extra maintenance to extend the durability.
- 2.Move this device with care in case any damages to the device.
- 3.After device using, please clean this device carefully and keep it properly for the next time use.
- 4.People who have wounds on skin are advised not to use.
- 5.Keep face moisturized and stay away from strong sunlight
- 6.Do not wash face with hot water (above 30°C)in three days after treatment

### Use Prohibited

- 1.Device cable is overheat
- 2.The plug does not connect tightly
- 3.The plug becomes bent
- 4.Do not unplug the device during device in working
- 5.Power cable of other device can not be used as a alternative for this device
- 6.Put this device in places where unreachable for kids
- 7.Cleaning with chemical liquids prohibited
- 8.Unplug device by wet hands prohibited

## AFTER SALE SERVICE

This device(main body of this device) enjoys 6 months warranty since the day of purchase. Any device (main body of this device) problems detected in 3 months, charge-free repair service is offered, repair fee is required over 3 months. The device accessories enjoy 1 month charge-free repair for surfaced quality problems, repair fee is required over 1 month.

**Following device malfunctions are not included in the range of charge-free repair service.**

- 1.Dismantle or alter the device by users themselves.
- 2.Accidentally fall during device using.
- 3.Device malfunctions caused by poor maintenance.
- 4.Operate the device without following the product manual.

### Service card

Product name:\_\_\_\_\_ Date of purchase:\_\_\_\_\_

Buyer' sname:\_\_\_\_\_ TEL:\_\_\_\_\_

Address:\_\_\_\_\_

Date of repair:\_\_\_\_\_ Malfunction description:\_\_\_\_\_

Repair statement:\_\_\_\_\_

