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### IMPORTANT SAFEGUARDS

The Power toothbrush is safe and effective when used according to directions. Use all electrical appliances, however, before safety precautions should be followed, especially when children are present.

#### DANGERS-WATER

- Do not use the toothbrush in bathrooms, swimming pools, hot tubs or other places where water can fall on the toothbrush.
- Do not use the toothbrush in bathrooms, swimming pools, hot tubs or other places where water can fall on the toothbrush.

#### DANGERS-VOLTAGE

- Do not use the toothbrush in bathrooms, swimming pools, hot tubs or other places where water can fall on the toothbrush.
- Do not use the toothbrush in bathrooms, swimming pools, hot tubs or other places where water can fall on the toothbrush.

#### DANGERS-WARNINGS

- Do not use the toothbrush in bathrooms, swimming pools, hot tubs or other places where water can fall on the toothbrush.
- Do not use the toothbrush in bathrooms, swimming pools, hot tubs or other places where water can fall on the toothbrush.

### CONFIGURATION

Professional brush head


Color ring

Power button

Mode indicator

Charging indicator

Charger bases
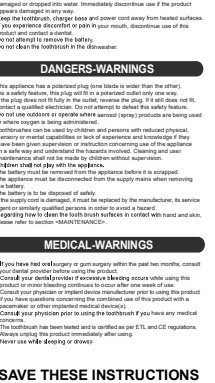
### SAVE THESE INSTRUCTIONS

### GETTING STARTED

#### Charging the toothbrush

Be sure to turn off the toothbrush before charging. Plug the power plug into the outlet. Then plug the toothbrush into the charger. The charger indicator light will be on, indicating the toothbrush is charging. After the battery is fully charged, the charger indicator light will turn off, indicating the toothbrush is fully charged.

#### HOW TO USE

1. Slide the brush head onto the handle. Push the brush head down until it locks into place on the handle.
2. Turn the handle and brush a small amount of toothpaste.
3. Turn the handle and brush a small amount of toothpaste.

### HOW TO USE

1. Slide the brush head onto the handle. Push the brush head down until it locks into place on the handle.
2. Turn the handle and brush a small amount of toothpaste.
3. Turn the handle and brush a small amount of toothpaste.

### HOW TO USE

1. Apply light pressure and gently move the brush head across the teeth in a semi-circle.
2. Apply light pressure and gently move the brush head across the teeth in a semi-circle.

### MAINTENANCE

Be sure to remove the plug from the electrical outlet before cleaning the product.

#### Toothbrush

1. Remove the brush head.
2. Wipe all surface areas with a damp cloth.
3. Wipe all surface areas with a damp cloth.

#### Charger

1. Using a soft cloth, wipe off any dirt on the charger and the settings.
2. Clean the plug periodically to ensure it remains in good working order.

### TROUBLESHOOTING

Problem	Possible cause	Solution
The toothbrush does not charge.	You have not purchased the correct power plug for your country.	Change the main unit and the power plug for your country.
The toothbrush is not fully charged.	The toothbrush is not fully charged.	Charge the main unit for at least 24 hours.

### WARRANTY

Your toothbrush is guaranteed for 1 year against manufacturer's defects.


### DISPOSAL

The entire surface of this product is designed to be recycled. At the end of the product's life, please remove the battery before you discard the product and recycle the remaining parts.

#### Removing and recycling the battery

1. Turn the handle and brush a small amount of toothpaste.
2. Turn the handle and brush a small amount of toothpaste.

### BINTOI

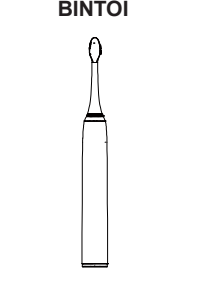Model: DT20

MADE IN CHINA