
C CARE

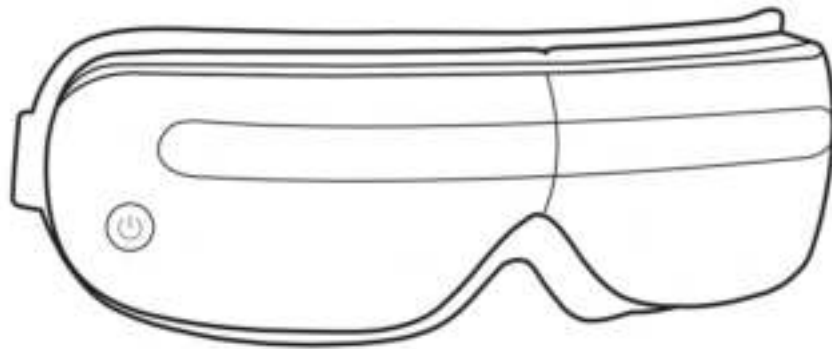SMC 110

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Service adres:

Cresta International B.V
Rolbrugweg 4 | NL-1332 AS ALMERE
info@cresta-international.com
www.cresta-international.com

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INSTRUCTIONS

We congratulate you on purchasing the C Care Eye Massager SMCT10 and thank you for the confidence placed in our brand and products.

We are convinced that with this purchase you have acquired a high quality product that you will enjoy using for many years to come.

Please read these operating instructions carefully before using the product, even if you are familiar with the type of product. It gives you useful tips to get the most out of your product.

If, after reading this manual, you still have questions about this product and how to use it, please send an e-mail to the address below; we will answer your question immediately

CRESTA INTERNATIONAL BV

Rolbrugweg 4

NL-1332 AS ALMERE

The Netherlands

info@cresta-international.com

www.cresta-care.com

Tel: +31 36 3030063

Product description

The intelligent air pressure eye massager applies the principle of traditional Chinese meridian and beauty, combined with the microchip control technology. The product improves the eye's blood circulation through gentle heat, vibrations and simulation of acupressure on the eye. With the built-in speaker and soothing music you can relax even more.

Scope

This product has a health function, it is suitable for young people with myopia, for students in school, employees who work intense with their eyes or work computers for many hours, writers, accountants, visually impaired people, people who sit in front of a screen of a computer or TV for many hours, people who stay up late, overuse of the eyes and the elderly who feel uncomfortable around the eyes.

Product features

Intelligent air pressure vibration to acupuncture points of the eye and eyelid, relaxes the eye and main nervous system.

The intelligent air pressure eye massager adopts the principle of traditional Chinese meridian and beauty, combined with the microchip control technology, the product improves eye blood circulation by gentle heat, vibration and simulation of acupressure on the eye. With the built-in speaker and smooth music, you can relax even more.

Intelligent air pressure vibration to acupuncture point of eye and eyelid, relaxes the eye and main nervous system

Constant soft, warm temperature.

The use of a carbon fiber element creates a constant temperature for the eye massager, which increases the blood circulation, increases the power of the eye nutrition and effectively relieves dry eyes.

Enjoy 10 minutes of relaxing music

Duration of use

After switching on, the eye massager will automatically switch off after 10 minutes. If you want to use the product longer, you should switch it on again.

Switch ON / OFF

Press and hold the ON / OFF button for 3 seconds to turn on the eye massager. By default, it is set in the automatic mode. When the eye massager is switched on, you can choose from 3 programs by briefly pressing the ON / OFF button.

When you switch on the eye massager, it will start in the automatic mode. There is pressure, heat and music.

C CARE

If you press the ON / OFF button again, the eye massager will switch to the dynamic mode. Here the massager works with pressure and music.

By pressing the ON / OFF button again, the massager will go into the sleeping mode. Here you will feel the soothing warmth and hear the relaxing music.

If you - temporarily - do not want to hear music, you can pause it by double clicking the ON / OFF. Repeating this action will restart the music.

Press and hold the ON / OFF button for 3 seconds and the eye massager will turn off.

Tips for use

Check that the massager has sufficient battery power. With low power the massager cannot be turned on. Charge the massager first and try turning it on again. Make sure to unplug the adapter after charging to avoid prolonged charging

Charging

Connect the massager to the USB port with the USB cable

The red light means that the massager is now charging. When the battery is fully charged, the light will glow blue.

The product must be switched off during the charging process. The self-protection program starts and the operation stops if you charge when the massager is turned on.

The massager is equipped with a built-in rechargeable lithium battery with a charging time of 2-3 hours.

Do not use the massager when the voice alarm indicates "low battery". When the batteries are fully charged, put the massager back into service (blue light).

Charging methods

You can connect the massager with the USB cable to a computer, a mobile power bank or the standard adapter. We advise you to use the standard adapter as much as possible

Tips for use

Check that the massager has sufficient battery power. With low power the massager cannot be turned on. Charge the massager first and try turning it on again. Make sure to unplug the adapter after charging to avoid prolonged charging

Precautions

Before use, check the precautions to make sure you are using the eye massager correctly.

Do not use the eye massager while walking or driving.

Do not use the eye massager under high temperatures, in a flammable environment, in an environment with electromagnetic radiation or under high humidity conditions.

If you have an eye disease, do not use the eye massager without consulting your doctor beforehand. This also applies to patients who have had eye surgery, who suffer from glaucoma, cataracts, retinal detachment or other eye conditions.

Patients with a severe stroke or suffering from high blood pressure should use the product with caution. Consult your treating physician in advance. If you experience any discomfort during use, stop using the product immediately.

Persons who cannot express themselves or who are insensitive to heat should only use this product under supervision.

Do not drink alcoholic beverages during use.

The eye massager has a warm compress function. If you are temperature sensitive, handle this with care.

First place the eye massager on the eyes and then turn it on to prevent damage to the air cushions.

First remove any contact lenses, hair clips, hair accessories, etc. that you wore before use.

Clean the eye protectors before each use.

This massage device is a personal health product and cannot replace medicines and medical equipment.

If you feel pain or severe discomfort while using, stop using immediately and consult a physician.

This eye massager glasses is designed as an electronic home care product and are not intended for use in diagnostic, therapeutic and other commercial applications.

Maintenance and cleaning

If cleaning is necessary, wipe it with a cotton cloth dampened with medical alcohol. Do not wash the product with water.

Do not disassemble this product yourself.

C CARE

Use the power adapter configured by our company and pay attention to electrical safety.

Please use this product with your eyes closed. After use, the eyes may be blurry. After closing the eyes for 1-3 minutes, the eyes will slowly recover.

Do not immerse the eye massager water or other liquid.

Do not bend the eye massager with force to avoid damage to the frame.

Clean the eye massager with a dry, soft cloth if there is dirt on the surface.

After use, keep the product in a safe and dry place.

Take measures to prevent moisture and pollution

Consult a professional technician if a repair of the unit is necessary or return it to our company. You will find the contact details at the front of this manual. Do not repair the product yourself, this will void your warranty claims.

Operating and transport conditions

Atmospheric pressure : 86 kpa - 106 kpa
Ambient temperature : -10 ° C - 65 ° C
Relative Humidity : ≤80% RH,

Use only in a well-ventilated area

Transportation: This product is equipped with simple anti-shock features, suitable for air, rail, road and marine transportation.

Prevent water on the product, shocks and impacts.

If you store the massager for a period longer than 6 months, check its effectiveness regularly.

Store the massager in a well-ventilated area and avoid direct sunlight.

TROUBLESHOOTING

PROBLEM	CAUSE	SOLUTION
PRODUCT DOES NOT SWITCH ON	LOW BATTERY POWER	CHARGE THE BATTERY
	PRODUCT SWITCHED OFF BY SELF-PROTECTION PROGRAM	CHARGE 5 SECONDS TO UNLOCK
	PRODUCT CANNOT BE CHARGED	ADAPTER DEFECTIVE
LOW MESSAGE POWER	LOW BATTERY POWER	CHARGE THE BATTERY

Product specifications

Product name : Eye massage glasses SMC110

Power: : 4.5 W

Voltage : DC 3.7V

Charger : 5V / 1A

WARRANTY

CRESTA INTERNATIONAL B.V. warrants this product to be free from manufacturing defects and material errors for a period of 36 months from the date of purchase. In the event of a defect within this period, CRESTA will at its discretion repair or replace the device free of charge for the same or an equivalent device. This warranty does not apply in the event of normal wear and tear or if, in the opinion of CRESTA INTERNATIONAL B.V. the defect was caused by improper use, falling or impact, water damage, leaking batteries, if applicable, and the like. The warranty is only applicable if the warranty claim can be demonstrated with the purchase receipt.

Do not try to open the product and repair it yourself. Your warranty claims will lapse

Requirements:

If a repair is necessary, you can return the product to the store where you bought it. You can also send the product directly to us. Send a message in advance with your complaint description and do not forget to state your full address details - not a PO box - when sending so that we can return the product to you after repair or replacement. If you send the product to us, you must frank it sufficiently. After repair or exchange you will receive your product from us delivered free.

Service address:

Cresta International B.V

Rolbrugweg 4 | NL-1332 AS ALMERE

info@cresta-international.com

www.cresta-international.com

ENVIRONMENT FRIENDLY DISPOSAL

The European Directive 2002/96/EC on Waste Electrical and Electronic Equipment - WEEE - requires that old household electrical appliances should not be disposed of in the normal unsorted municipal waste stream. Old appliances should be collected separately in order to optimize the recovery and recycling of the materials they contain and reduce the impact on human health and the environment.

 The crossed out "wheeled bin" symbol on the product reminds you of your obligation that when you dispose of the appliance, it should be separately collected. Contact your local authority for information concerning the correct disposal of your old appliance.

PRINTING ERRORS AND MODIFICATIONS RESERVED

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