
C CARE

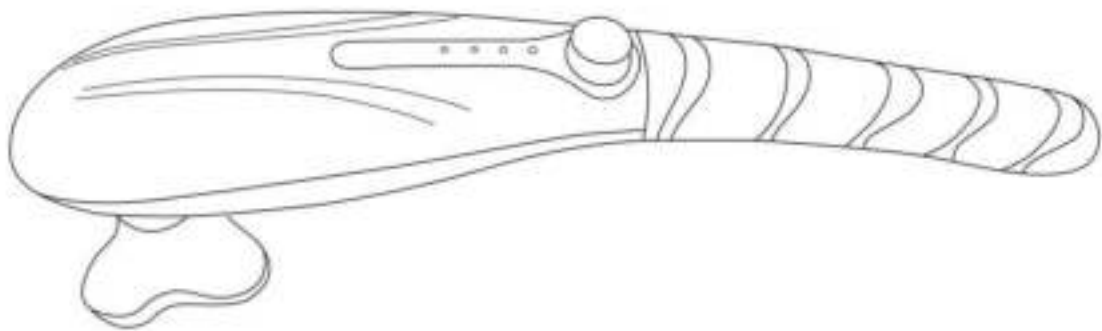SMC 100

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Service adres:
Cresta International B.V
Rolbrugweg 4 | NL-1332 AS ALMERE
info@cresta-international.com
www.cresta-international.com

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INSTRUCTIONS

We congratulate you on purchasing the C Care Handheld Massager SMC100 and thank you for the confidence placed in our brand and products.

We are convinced that with this purchase you have acquired a high quality product that you will enjoy using for many years to come.

Please read these operating instructions carefully before using the product, even if you are familiar with the type of product. It gives you useful tips to get the most out of your product.

If, after reading this manual, you still have questions about this product and how to use it, please send an e-mail to the address below; we will answer your question immediately

CRESTA INTERNATIONAL BV

Rolbrugweg 4

NL-1332 AS ALMERE

The Netherlands

info@cresta-international.com

www.cresta-care.com

Tel: +31 36 3030063

Product description

This electronic handheld massager will provide you a penetrating and relaxing tapping massage with the intention to relax and loosen your muscles. The use of the product will relieve muscle aches and pain, will reduce stress, relieve tension and will improve your blood circulation.

Scope

This product has a health function, it is suitable to be used by all people over the age of 8 years.

Safety Instructions

To ensure a safe and proper use of this product read this instruction manual carefully first

Usage environment

- Do not use the appliance in a wet or dusty environment as this may result in malfunction or an electric shock.
- Do not use the appliance in the bathroom or in bed
- Do not place or store the appliance where it can fall or be pulled into a tub or sink. Store it in a safe location.
- Do not sit on the product
- Do not expose the appliance to a heater or direct sunlight.
- Do not drop or insert any object into any opening
- Handle the heated surfaces with great care . The unattended use of the heated surfaces by children or disabled may be dangerous
- Do not lift or carry this appliance by the power cord
- This appliance is designed for personal use only, not for professional purposes
- Do not cover the appliance during the operation
- Do not crush the appliance and avoid any sharp folds
- The appliance is not intended for persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge. They can use the product under supervision of a person responsible for their safety.
- Extensive use can lead to excessive heating and physical damage. Should this occur discontinue the use and allow the unit to cool down before operating it again.
- For people under medication or medical treatment please consult your doctor before using the product. Do not use the product if you do not feel well.
- If the product is stored in a cold environment and is being taken into a warm environment allow the product to adopt an hour to the warmer environment before use to get the optimal massaging performance
- Do not open the product, do not repair by yourself.

Main unit & accessories

The handheld massager can be used with 4 different massage heads

- 4 point intensive massage head
- Comb massage head
- Shiatsu massage head
- Triangle intensive massage head

Product operation guide

Before using the product select the proper massage head and plug the appliance to the power socket.

Turn the control button clockwise to switch on the appliance, you will hear a soft "click" to confirm the activation. The massage head will start the tapping movement.

If you wish to increase the intensity of the vibration turn the control button gradually in the clockwise direction. To decrease the intensity gradually turn the control button anti clockwise.

To switch of the product turn the control button anti clockwise till you hear a "click" sound.

After the use unplug the appliance from the power socket.

Battery condition and charging

To charge the battery attach the charging adapter to the charging port on the handle of the product and plug the adapter into the power socket. In the control panel a light will indicate the battery condition.

If the light is solid red it indicates the battery has less than 15% power

If the light is blinking blue it indicates the battery is charging

If the light is solid blue it indicates the battery is charged 100%

When the battery is fully charged unplug the adapter from the power socket and remove the charging plug from the main unit.

IMPORTANT

When the appliance is charging do not touch the head

The unit should be switched off during charging.

Product maintenance

Storing the unit

- Place the massager preferably in the original giftbox and store it in a safe, dry and cool place. Avoid contact with sharp edges or pointed objects that might cut or puncture the surface.
- Wind up and store the power cord under the appliance bottom frame. Do not hang the unit by the power cord.

Cleaning

Unplug the unit and allow it to cool down before cleaning

Clean only with a soft, damp sponge

Do not allow water or any other liquid to go inside the unit

Avoid immersing the unit into any liquid during the cleaning

Do not use abrasive cleaners, brushes, gasoline, kerosene, glass/furniture polish or paint thinner to clean the product

Caution

- Avoid tearing on the appliance surface
- Only gentle force should be exerted against the unit in order to eliminate any risk of injury
- If you are experiencing any health complications or health conditions like suffering from a heart disease, pregnancy or any other health concern consult your doctor before using the product
- Do not exceed the use of more than 15 minutes at one time

C CARE

- The product should not be used by children or disabled without adult supervision.
- Do not use the product just before going to sleep
- Do not use the product if it is damaged
- Do not use the product if you do not feel well or experience severe pain during the massage

Troubleshooting

PROBLEM	SOLUTION
Not operating	<ol style="list-style-type: none">1. Power failure, use when the electricity resumes2. Check the connection with the power socket
During the operation the controller is working but no massage	The massager could be overheated after prolonged use. It will automatically cool down to resume its normal function. Please turn off the power and allow the massager to cool down for 30 to 50 minutes
No massage function	<ol style="list-style-type: none">1. Turn on the power switch2. Reconnect the power cord
Sudden stop	<ol style="list-style-type: none">1. Power failure2. Reconnect the adapter
Massager does not work	<ol style="list-style-type: none">1. Make sure the massager is fully charged, the control panel light should be solid blue2. Reset the massager by plugging it into the charger for at least 3 seconds then unplug and turn it on again3. The lithium-ion battery life may be expired
The massager does not start immediately	<ol style="list-style-type: none">1. It can take up to 3 seconds for the massager to start. Ensure you turned the control button enough2. Hear a click3. Check the power light is blue,

Product specifications

Product name	:	Rechargeable Handheld Massager
Reference	:	SMC100
Power Source	:	100-240V
Power Frequency	:	50-Hz
Power consumption	:	9 W
Adapter output	:	DC 9V, 1A
Product dimensions	:	412 x 97 x 70 mm

ENVIRONMENT FRIENDLY DISPOSAL

The European Directive 2002/96/EC on Waste Electrical and Electronic Equipment - WEEE - requires that old household electrical appliances should not be disposed of in the normal unsorted municipal waste stream. Old appliances should be collected separately in order to optimize the recovery and recycling of the materials they contain and reduce the impact on human health and the environment.

The crossed out "wheeled bin" symbol on the product reminds you of your obligation that when you dispose of the appliance, it should be separately collected. Contact your local authority for information concerning the correct disposal of your old appliance.

WARRANTY

CRESTA INTERNATIONAL B.V. warrants this product to be free from manufacturing defects and material errors for a period of 36 months from the date of purchase. In the event of a defect within this period, CRESTA will at its discretion repair or replace the device free of charge for the same or an equivalent device. This warranty does not apply in the event of normal wear and tear or if, in the opinion of CRESTA INTERNATIONAL B.V. the defect was caused by improper use, falling or impact, water damage, leaking batteries, if applicable, and the like. The warranty is only applicable if the warranty claim can be demonstrated with the purchase receipt.

Do not try to open the product and repair it yourself. Your warranty claims will lapse

Requirements:

If a repair is necessary, you can return the product to the store where you bought it. You can also send the product directly to us. Send a message in advance with your complaint description and do not forget to state your full address details - not a PO box - when sending so that we can return the product to you after repair or replacement. If you send the product to us, you must frank it sufficiently. After repair or exchange you will receive your product from us delivered free.

Service address:

Cresta International B.V

Rolbrugweg 4 | NL-1332 AS ALMERE

info@cresta-international.com

www.cresta-international.com

PRINTING ERRORS AND MODIFICATIONS RESERVED

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