




## U-shaped electric toothbrush

Well-designed · easy to operate · safe and convenient



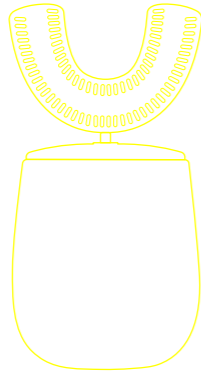

(This manual is a general manual, and the following single models are the main ones)

## Instructions

- 


1 Evenly on both sides of the brush head Squeeze on mousse toothpaste

2 Put the brush head in the mouth, up and down two Put the teeth in the tooth socket at the same time, hold gently
- 


3 Don't bite hard, use left and right hands Shake the body to touch the sides Teeth, the best frequency is 63 beats per minute

4 After brushing your teeth, rinse the toothbrush Just wipe dry

## Product Features


- USB contact charging:  
The charging time is no more than 120 minutes, and it can last for more than a week on a single charge, which is energy-saving and environmentally friendly.
- Convenient touch:  
One-button operation, six gear options: A blue light whitening, B standard cleaning, C excellent cleaning, D white cleaning, E comfortable massage, F pulse massage.
- 360 degree cleaning:  
After the switch is activated, it only takes 45 seconds (8 series)/59 seconds (k series) to complete a 360° brushing action, which is quick and convenient.
- Automatic memory function:  
Smart memory user habits, efficient and convenient.

## Features


- The toothbrush adopts an all-round waterproof and sealed design, IPX7, which is convenient, durable and clean.
- U-shaped brush head design, using food-grade silicone antibacterial material, 360° contact with gums and teeth, every time Brushing teeth can bring a quick and comfortable experience to the oral cavity.
- Using this product can prevent gum atrophy, promote periodontal blood circulation, and reduce the occurrence of periodontal diseases.

## Charging Instructions


- Insert the DC end of the charging cable into the back of the toothbrush vertically In the round hole at the bottom, the USB end is connected to the charging plug Port for charging.
- When charging, the LED light on the body flashes with breathing, The LED light goes out when it is full.

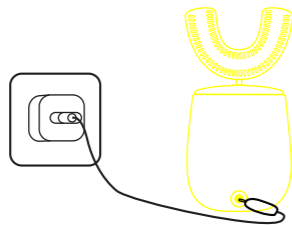


## Product Structure

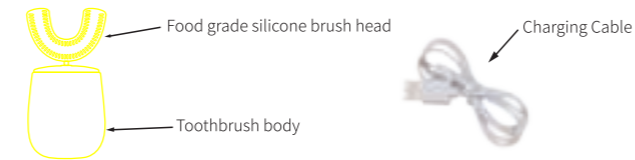

## Brush head replacement method

- 

1 Use scissors to put in the middle of the bottom of the brush head (shaft) Evenly exert force on both sides of the body, slowly tilting upwards
- 

2 Hold the brush holder with your thumb and index finger Central position, pull out slowly

## Precautions

- The normal working temperature of the product is 0°C-50°C.
- When charging the product, due to the particularity of the product, please charge it at room temperature and keep the body dry. The power supply is less than 5V1A.
- When the product is not used for a long time or stored for a long time, it is recommended to charge and discharge once within three months.
- Store the product in a dry and cool place as much as possible, and avoid sunlight.
- After brushing your teeth, rinse the brush head with clean water, gently shake off the surface water droplets, try to avoid soaking in water.
- The output power of the product is large, and it is forbidden to work for a long time, generally less than 10 minutes.

## Product parameter

Rated voltage	5V	Rated power	<5W
Working current	<1.6A	Charging time	<120 Minute

## Certificate of conformity

Grade: First grade Test: 08

Executive standard: GB4706.1-2005/GB4706.59-2008