panacea

USER MANUAL

FOOT MASSAGER
WITH HEAT

Model: CM-003F

Index

| | |
|------------|-------|
| Nederlands | 4-7 |
| English | 8-11 |
| French | 12-15 |
| Deutsch | 16-19 |

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Thank you for purchasing the Panacea foot massager with heat. Please read all parts of this manual carefully before use and keep them for future reference. If you need further assistance, please contact our customer service: info@panacea.nu or via WhatsApp on +31630972659

Table of contents

1. Features
2. Safety measures
3. Warnings
4. During use
5. When the machine is not in use
6. When not to use the device
7. Explanation of the device
8. How to use?
9. Maintenance and attention points
10. Package contents
11. Troubleshooting
12. Specifications

1. FEATURES

1. Breathable fabric & toe-touch controller.
2. Six massage heads with 24 rotating massage heads for deep penetration of the muscles, for a pain relieving and relaxing experience.
3. Kneading, rolling & shiatsu massage from toes to heels to treat all acupuncture points.
4. There is a heating function, which can be turned on/off manually; two massage speeds; two working modes available (clockwise & anti-clockwise rotation)
5. Easy to carry with the built-in handles. Enjoy massage at home, in the office or when travelling.
6. Powered by DC12V Adapter, safe and reliable.

2. SAFETY MEASURES

The product must be used correctly according to the instructions in this manual in order to avoid physical and material damage.

3. WARNINGS

Persons with one or more of the following conditions or undergoing medical treatment should consult their physician before using the Panacea foot massager:

- Pregnant or in childbed;
- Have a fever;
- Not recommended by your doctor, or not feeling well;
- Recovering on doctor's prescription;
- You are suffering from a chronic illness or cancer;
- Have an injury to the brain (clots, cuts, traumatism, etc.);
- You are recovering from an operation on the upper body, buttocks, thighs or calves;
- You suffer from a heart disease, heart problems or you use an electronic medical device;
- You have a pacemaker;
- You are injured or have a skin disease;
- You have broken or malformed bones;
- This device is not suitable for children, disabled persons or invalids.

4. DURING USE

- Do not place the device in a place which is humid or very dusty as this may lead to mechanical failure;
- Do not place the appliance in a heat source or in direct sunlight;
- Place the massager on a flat surface;
- Do not use the device together with other therapeutic devices or electrical medical equipment or similar equipment;
- Avoid splashing liquid/water on the device when in use;
- Place the device on a flat, hard surface before use and avoid tilting, vibrating or shaking the device;
- Do not use the device with wet feet or other wet body parts; as this may lead to mechanical failure or even an electric shock.
- Your body should be relaxed during the treatment;
- Do not use the device immediately after a meal, wait at least ½ hour;
- Limit the daily use to a maximum of 40 minutes.

Attention! Make sure that the voltage of your outlet is suitable for the device and never unplug the device while it is on.


5. WHEN THE MACHINE IS NOT IN USE

- Make sure the machine is switched off after use;
- Unplug it after use and when you want to clean it;
- Store the appliance in a clean, dry and even place and do not place heavy objects on it.

6. WHEN NOT TO USE THE DEVICE

- If the device falls into water: The device should not come into contact with water or other liquids;
- If you feel strange or very uncomfortable while using the device, turn it off immediately and consult your doctor;
- If the power plug or cord is damaged or does not work properly.

7. EXPLANATION OF THE DEVICE

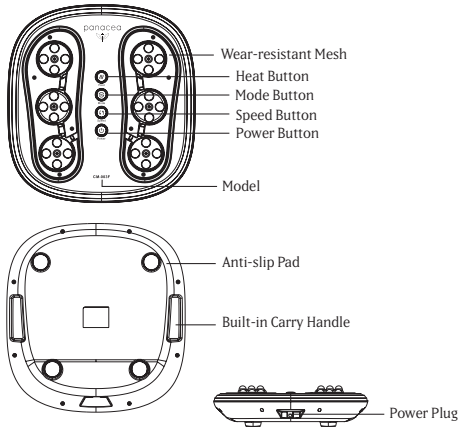8. HOW TO USE?

1. Press the button to start the massage (all button lights on). The massager will start in slow speed mode 1 by default and will last for 20 minutes

Note: The heating function is on by default, you can press the heat button to turn it off.

2. Press the button to select the massage mode: Mode 2 Switch the kneading direction: clockwise or anti-clockwise

3. Press to select the speed of the massage: Slow or Fast.

4. Press the button to turn off the heating function.

5. Press the button to end and switch off the massage.
Note: The unit switches off automatically after 20 minutes, you can restart it at any time if you wish.

9. MAINTENANCE AND ATTENTION POINTS

- Please unplug the device after use;
- Make sure the device is unplugged while cleaning;
- Use a soft cloth for cleaning. Never use a damp cloth, brushes or hard materials to clean the device;
- Do not use the device for more than 40 minutes per day.

10. PACKAGE CONTENTS

- Panacea foot massager
- Power adapter
- User Manual
- Washable protective cover

11. TROUBLESHOOTING

| Problems | Possible Causes | Solution |
|--------------------------------|---|--|
| The device does not work | 1. Power loss 2. Poor connection of power connector | 1. Turn on the device 2. Insert the power plug firmly into the socket |
| Stopped while using the device | If you use the device for too long, the device will overheat and stop automatically to protect the device | Turn the device off and disconnect the power supply, and then turn the device back on after an hour. |
| | The 20 minute timer stops | Press the power button to restart the device if you want a longer massage. |
| Noise during use | It is normal for the massager to make noise during use. This does not affect the use of the device. | You do not need to do anything. If you have any questions, please contact info@panacea.nu |

12. SPECIFICATIONS

- Product Name: Panacea foot massage device
- Model No. CM-003F
- Voltage: DC12V/3.0A
- Power Output: 36W

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