# Oral-B®

## VITALITY™








powered by **BRAUN**

Charger Type 3757  
Handle Type 3710

91571849

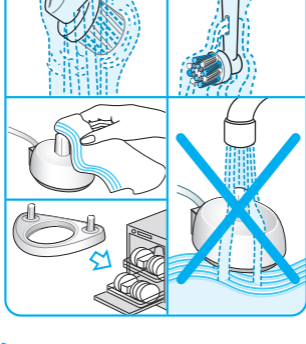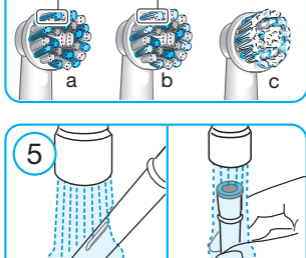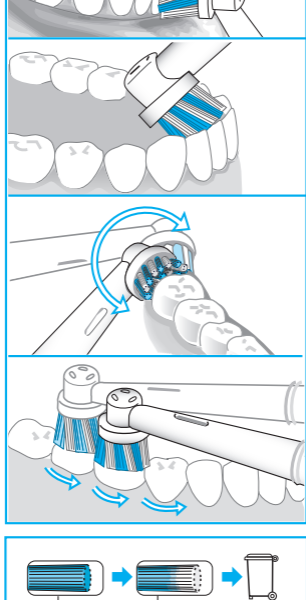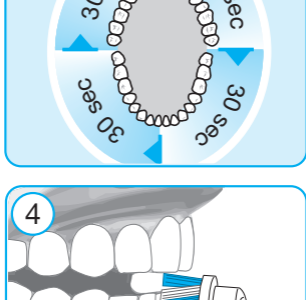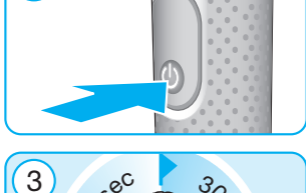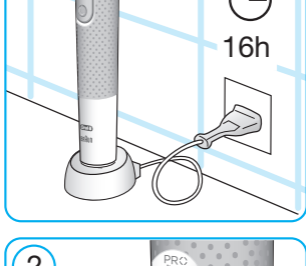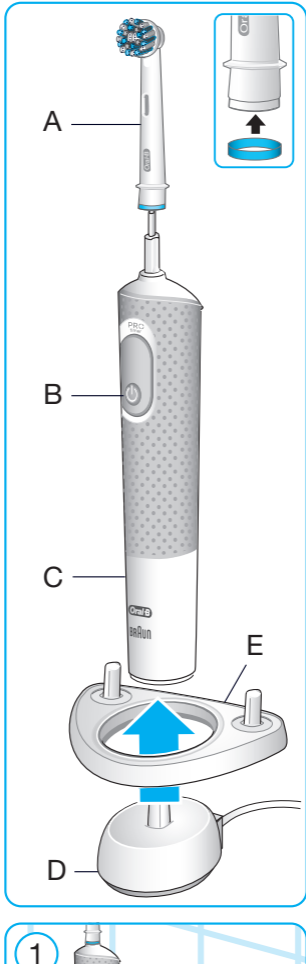- (DE) (AT) 00 800 27 28 64 63**
- (CH) 08 44 - 88 40 10**
- (GB) 0 800 731 1792**
- (IE) 1 800 509 448**
- (FR) 0 800 944 803**  
(service et appel gratuits)
- (BE) 0 800 14 592**
- (LU) 800 21172**
- (ES) 900 814 208**
- (PT) 808 20 00 33**
- (IT) 800 124 600**
- (NL) 0800-4 45 53 88**
- (DK) 70 15 00 13**
- (NO) 22 63 00 93**
- (SE) 020 - 21 33 21**
- (FI) 020 377 877**

**Internet**  
[www.oralb.com](http://www.oralb.com)  
[www.braun.com](http://www.braun.com)  
[www.service.braun.com](http://www.service.braun.com)

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Welcome to Oral-B!

Before operating this toothbrush,

please read these instructions and

save this manual for future reference.

## IMPORTANT

- Periodically check the entire product/cord for damage. If the product/cord is damaged, take it to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
  - Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
  - Cleaning and maintenance shall not be performed by children.
  - Children shall not play with the appliance.
  - Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- WARNING**
- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
  - This appliance contains batteries that are non-replaceable. Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
  - When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
  - If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
  - This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

To avoid brush head breakage which may generate a small parts choking hazard or damage teeth:

- Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fit properly. Never use without a brush head.

- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible

- Replace brush head every 3 months or sooner if brush head becomes worn.

- Clean brush head properly after each use (see "Cleaning recommendations" paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.

## Description

A Brush head

B On/off button

C Handle

D Charging unit

## Accessories

(depending on model):

E Brush head holder (depending on model)

**Note:** Content may vary based on the model purchased.

## Specifications

For voltage specifications please refer to the bottom of the charging unit.

Noise level: ≤ 65 dB (A)

## Charging and Operating

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet and place the toothbrush handle (C) on the charging unit. A full charge takes about 16 hours and enables up to five days of regular brushing (twice a day, 2 minutes) (picture 1).

- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible.

- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

## Using your toothbrush

### Brushing technique

- Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on the appliance (picture 2).

- When brushing your teeth with one of the **Oral-B oscillating-rotating brush heads** guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface.

- When using the **Oral-B «TriZone» brush head** place the toothbrush bristles against the teeth at a slight angle towards the gumline. Apply light pressure and start brushing in back and forth motions, just like you would do with a manual toothbrush.

- With any brush head start brushing the outsides, then the insides and

finally the chewing surfaces. Brush all four quadrants of your mouth equally. Do not press too hard, simply let the brush do all the work (picture 4). You may also consult your dentist or dental hygienist about the right technique for you.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

#### **Timer** (depending on model)

An incorporated timer memorises the elapsed brushing time (depending on model) (picture 3).

**«2-Minute» timer:** A short stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

**«Professional» timer:** A short stuttering sound at 30-second intervals reminds you to brush equally all four quadrants of your mouth. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

### **Brush heads**

#### **Oral-B CrossAction brush head**

Our most advanced brush head. Angled bristles for a precise clean. Lifts and powers away plaque.

#### **Oral-B FlossAction brush head**

features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.

#### **Oral-B 3D White brush head**

features a special polishing cup to naturally whiten your teeth. Please note that children under 12 years old should not use the Oral-B «3D White» brush head.

#### **Oral-B Sensi Ultrathin brush head**

features a combination of regular bristles to clean teeth surfaces and ultrathin bristles for a gentle brushing experience on the gumline.

Most Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums (picture 4).

We do not recommend using the Oral-B «FlossAction» or Oral-B «3D White» brush head with braces. You may use the Oral-B «Ortho» brush head, specifically designed to clean around braces and wires.

### **Cleaning recommendations**


After brushing, rinse your brush head for several seconds under running water with the handle switched on. Switch off the handle and remove the brush head. Rinse both parts separately under running water then wipe them dry before reassembling the toothbrush. From time to time also clean the charging unit with a damp cloth; never place the charging unit in water.

Brush head holder (E) is dishwasher-safe (picture 5).

Subject to change without notice.

#### **Environmental Notice**

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

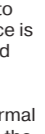Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

#### **Guarantee**

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

#### **UK only:**

This guarantee in no way affects your rights under statutory law.

#### **Warranty replacement brush heads**

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads. Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

#### **30 Day Risk Free Trial**

Satisfaction guaranteed or your money back!

Full details at [www.oralb.com](http://www.oralb.com) or call:  
**GB 0800 731 1792**  
**IE 1800 509 448**

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