A i r F y x

Red Light Therapy

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Instruction manual

Manuel d' Instructions

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Congratulations with the purchase of your AirFyx Red Light Therapy. The first step on taking control of your pain!

Airfyx uses LED Red Light Therapy to provide direct, non-addictive, drug free pain relief. Simplify your life with AirFyx Red Light therapy.

Please read the instructional manual to get the most out of your AirFyx Red Light Therapy. Airfyx uses a combo of red (630nm, 660nm, and near infrared 850nm) LEDs which is chosen to give the ideal light intensity for biological effects anywhere on the body.

● What is Red Light Therapy?

Red light therapy (RLT) is a treatment that may help skin, muscle tissue, and other parts of your body heal. It exposes you to low levels of red or near-infrared light. Infrared light is a type of energy your eyes can't see, but your body can feel as heat. Red light is similar to infrared, but you can see it.

● How Does Red Light Therapy Work?

With red light therapy, you expose your skin to a lamp, device, or laser with a red light. A part of your cells called mitochondria, sometimes called the “power generators” of your cells, soak it up and make more energy. Some experts think this helps cells repair themselves and become healthier. This spurs healing in skin and muscle tissue. Red light therapy uses very low levels of heat and doesn't hurt or burn the skin. It's not the same type of light used in tanning booths, and it doesn't expose your skin to damaging UV rays.

● What are the benefits?

This red-light device provides powerful pain relief without causing harmful side effects. Unlike traditional pain relief methods, or ways that are expensive and can cause a myriad of side effects, our device is an affordable way to reduce aches and pains in your body.

Infrared light is best for superficial skin, reduce wrinkles, smooth skin. Proven to reach the pain point directly and relieve pain more safely and effectively.

Temporary relief of:

- Minor muscle and joints pain
- Arthritis,
- Muscle spasm
- Relieving stiffness
- Promoting relaxation of muscle tissue

● Included with this package

- 1* Red Light Therapy Torch
- 1* Carry Bag
- 2* Rechargeable Battery 18650
- 1* Battery Charger
- 1* USB Cable
- 1* Protective glasses
- 1* User Manual
- 1* Microfiber cloth

● Important to know

- Commit to a 3-Week Treatment Plan. While many people achieve pain relief from red-light therapy almost immediately, others see measurable pain relief incrementally over time. Depending upon injury type, severity and different pain pathways involved, people vary in therapy response time. This is why for we recommend committing to a 3-Week Treatment Plan. If you see results right away, great, continue treating as needed. However, if you do not feel relief as quickly, know that this is also normal, so try not to get frustrated. Simply treat your pain consistently with 3 quick treatments per day with AirFyx Red Light Therapy, for up to 3-5 minutes each time per treatment spot. After 3 weeks of consistent red-light therapy, most people experience a diminishing of their pain.
- Do not suddenly stop your other treatments upon beginning red-light therapy. Although you may not be satisfied with current treatments, they may be having a partial effect on the level of your pain. If you stop suddenly, you may experience a large increase in pain. This can be a sort of “rebound” pain effect or pain “flare up” which may be very difficult to get back under control. Some medications may also cause withdrawal like side effects if stopped too quickly, especially prescription medicines. While continuing with your current treatments, start with your AirFyx Red Light Therapy, giving it time to provide optimal pain relief.

After you have achieved maximum pain relief from your AirFyx Red Light Therapy, you may then consider for slowly weaning yourself off of other medications or treatments. If you plan to reduce or stop a prescription medication, first consult with your physician. For many people, red-light therapy works so well, that they are able to come off pain pills completely. Others are able to reduce the amount of pain medications they are on, reducing the risk of drug side effects.

● Directions for use

Hold the unit over the painful area at a distance anywhere from touching the skin to 3 cm away from the skin, depending on the depth of treatment desired. The red-light penetrates up to 2.5 cm into your tissues or joints. Hold the unit on the treating area for 5 minutes, guides as below:


Pain Level	Light Pain	Average Pain	Deep Pain
Dosage	Per Treatment Spot	Per Treatment Spot	Per Treatment Spot
	Once a day	Twice or 3 a day	Three times a day

Important

Do not keep the device **ON** for more than 5 minutes at a time. If you want to treat several parts of your body consecutively, **TURN OFF** the device for at least 1 **MINUTE**, best is to turn the device after 5 minutes **ON**. This will keep the unit from overheating.

“The treatment spot is an area covered by the red-light beam.” A small joint ie. knuckle of a finger, should have only one treatment spot, because the light shines over the entire affected joint.

↪ **For wrist treatment**

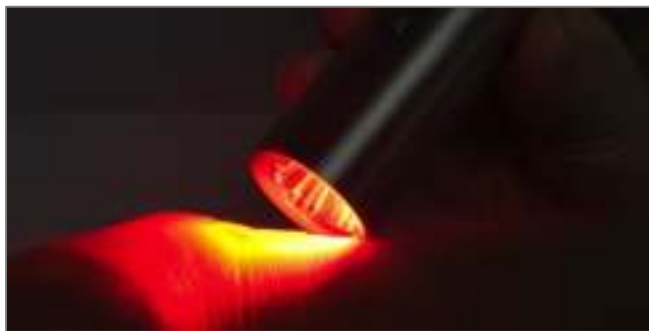↪ **For knee treatment.**

↪ **For finger treatment.**


● How to install the battery?

Insert the Battery provided on the battery compartment

- Positive Pole against the TOP.
- Negative Pole against the BOTTOM

Please only use 18650 batteries

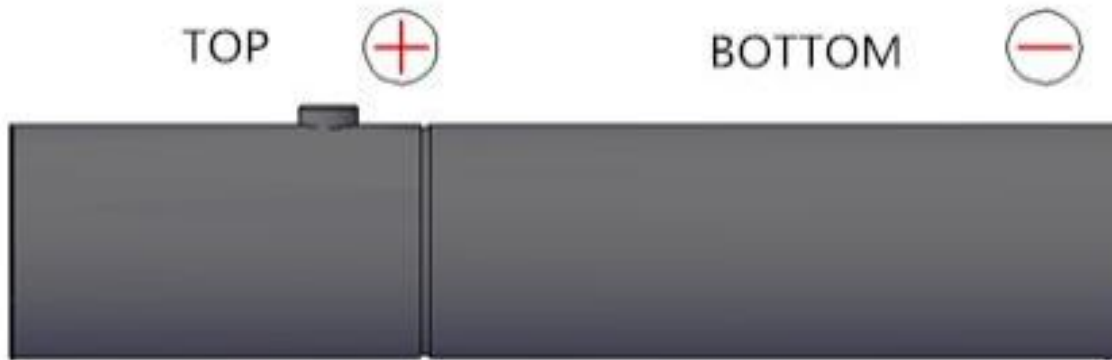● Product specification

Product Size	D26*120.5mm	Net Weight	0.21 KG
Material	Stainless steel	Wavelength	630:660:850nm
Irradiance	2.54cm: 141mW / cm²	Irradiance	7.62cm: 30.1mW / cm²
Illuminance	2.54cm: 141mW / cm²	Battery	18650
Warranty	12-Month	Lifespan	50,000 Hours

● Caution

- Don't use the light around water.
- Don't stare directly at the light. Use protective glasses while using the device
- Don't treat near the eyes.
- The light should NOT be used by children under the age of 13 and pregnant woman.
- Don't use it more than 5 minutes at a time for 3 times a day.

● Equipment maintenance

- **CLEANING AND DISINFECTING:** Use only a soft dry cloth to clean the exterior. Never use an abrasive cleanser or submerge in water or any other liquid. This device cannot be sterilized by any liquid or autoclave method.
- **STORAGE:** Store the unit in a cool, dry place, away from direct light and extreme heat.
- **WARNING:** Do not allow moisture to come in contact with the case or light of the unit. Do NOT immerse the unit in water.
- **WARNING:** If unit will not be used for an extended period of time (1 month +), remove the batteries from the unit to avoid their Corrosion.

“The symbol indicates that according to national laws and regulations, your products or batteries must be treated separately from household garbage. Please follow these disposal instructions when discarding products or batteries”

● Warranty card

Thank you for purchasing our products. To safeguard your right, please keep this warranty service card. Please feel free to contact our dealers in case of quality issues.

Product name:	AirFyx Red Light Therapy
Model:	Focus
Name:	
Purchase date:	
Order number:	
Phone:	

Notice

- Users receive a 12-month warranty from the date of purchase.
- Show this warranty and proof of purchase or invoice when requesting a repair.
- Extra parts and spare parts are not covered by the warranty.
- The sellers are not liable for any medical or other type of damage when using this product. Do not allow water or moisture to enter the device. You cannot claim under the warranty in the event of negligence and incorrect use of the product.
- The warranty conditions may change over time. As a result, only the most recent warranty conditions are valid for claiming your warranty. The most recent terms and conditions can be found on the official website of AirFyx.com

Non Warranty

- For improper or inappropriate use and self-caused damage.
- When this warranty certificate is lost or damaged.
- Damage due to natural disasters or other causes of force majeure.

Thank you for using our products!

We are committed to providing much better products, more attentive service, more favorable price and top-class service to consumers. We have a wide product range to meet your demand.

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