# sanitas

## SMG 11

<b>DE</b>	<b>Mini-Massager</b> <i>Gebrauchsanweisung</i> .....	2	<b>NL</b>	<b>Mini-massageapparaat</b> <i>Gebruiksaanwijzing</i> .....	48
<b>EN</b>	<b>Mini massager</b> <i>Instruction for Use</i> .....	11	<b>PL</b>	<b>Miniaparát do masažu</b> <i>Instrukcja obsługi</i> .....	58
<b>FR</b>	<b>Mini masseur</b> <i>Mode d'emploi</i> .....	20			
<b>IT</b>	<b>Mini massager</b> <i>Istruzioni per l'uso</i> .....	29			
<b>RU</b>	<b>Мини-массажер</b> <i>Инструкция по применению</i> ..	38			

CE

## DEUTSCH

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## ENGLISH

**Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.**

### **WARNING**

- The device is intended only for domestic/private use, not for commercial use.
- This device can be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and are aware of the consequent risks of use.

- Children must not play with the device.
- Never open or repair the device. Failure to observe this will invalidate the warranty.
- Cleaning and user maintenance must not be performed by children unless supervised.

## Contents

1. Included in delivery.....	12	6. Initial use .....	17
2. Signs and symbols .....	13	7. Cleaning and maintenance.....	18
3. Warnings and safety notes.....	15	8. Disposal .....	18
4. Proper use .....	15	9. Technical data .....	19
5. Description de l'appareil.....	16	10. Warranty / Service .....	19

## 1. Included in delivery

Check that the device packaging has not been tampered with and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

- 1 x mini massager
- 3 x 1.5 V battery (AAA)

## 2. Signs and symbols

The following symbols appear in the instructions for use and on the device:

	<b>WARNING</b>	Warning notice indicating a risk of injury or damage to health.
	<b>IMPORTANT</b>	Safety note indicating possible damage to the device/accessory.
	<b>Note</b>	Note on important information.
	Marking to identify the packaging material. A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard	
	Separate the packaging elements and dispose of them in accordance with local regulations.	
	Importer symbol	
	Separate the product and packaging elements and dispose of them in accordance with local regulations.	
	United Kingdom Conformity Assessed Mark	

	<p>Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).</p>
	<p>Manufacturer</p>
	<p>CE labelling This product satisfies the requirements of the applicable European and national directives.</p>
	<p>Read the instructions for use</p>
	<p>Direct current The device is suitable for use with direct current only</p>
	<p>The products demonstrably meet the requirements of the Technical Regulations of the EAEU.</p>
	<p>Do not dispose of batteries containing hazardous substances with household waste</p>

### 3. Warnings and safety notes

#### **WARNING**

- Do not use any attachments that are not recommended or offered as accessories by the manufacturer.

### 4. Proper use

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper use.

This device is intended only for the massaging of individual parts of the human body. Do not use the massage device if any of the following warnings apply to you.

#### **WARNING**

Use of the massage device is prohibited in the following cases:

- On children.
- Not with changes caused by disease or injury in the back region (e.g. burns, open wounds, scars, skin rashes or other injuries).
- Not in case of cardiovascular diseases and headaches.
- Not on the head, bones (e.g. joints, spine) or other sensitive parts of the body.
- Not after operations.
- The massage device must not be used during pregnancy.
- Never use the device after taking medication or consuming alcohol (reduced levels of awareness).
- Do not use the device while operating machinery.

- For external use.

Use of the device does not replace medical consultation and treatment.


You should consult your doctor before using the massage device if:

- you are fitted with a cardiac pacemaker, implant or other aid,
- in case of pains of unclarified origin,
- if you **are not** certain whether the device is suitable for you.

### **Note**

The massage must be pleasant and relaxing at all times. Stop the massage or change your position or the contact pressure if you find it is painful or unpleasant.

## **5. Description de l'appareil**

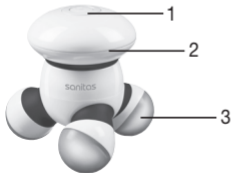1. ON/OFF button (1 x on, 2 x off)
2. Battery compartment
3. Massage head with LED

## 6. Initial use

Inserting or replacing batteries:

- Twist off the battery compartment lid counter-clockwise.
- Insert the batteries into the battery compartment correctly.
- Place the battery compartment lid on the battery compartment.

### Notes on handling batteries

- If your skin or eyes come into contact with battery fluid, flush out the affected areas with water and seek medical assistance.
- ⚠ **Choking hazard!** Small children may swallow and choke on batteries. Store the batteries out of the reach of small children.
- Observe the plus (+) and minus (-) polarity signs.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Protect the batteries from excessive heat.
- ⚠ **Risk of explosion!** Never throw batteries into a fire.
- Do not charge or short-circuit batteries.
- If the device is not to be used for a long period, take the batteries out of the battery compartment.
- Use identical or equivalent battery types only.
- Always replace all batteries at the same time.
- Do not use rechargeable batteries.
- Do not disassemble, split or crush the batteries.

## 7. Cleaning and maintenance

Your massager is maintenance-free. Only the outer housing needs to be cleaned

regularly. Remove the batteries from the battery compartment for cleaning. Clean your Mini Massager carefully using only a slightly moistened cloth.

- Do not use any detergents or solvents.
- Never hold the instrument under water as otherwise liquid can penetrate and damage the device.

## 8. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its service life.

Dispose of the device at a suitable local collection or recycling point in your country. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries. The codes below are printed on batteries containing harmful substances:

Pb = Battery contains lead,

Cd = Battery contains cadmium,

Hg = Battery contains mercury.

## 9. Technical data

Weight	Approx. 162 g
Device dimensions	9.8 x 9.8 x 9.5 cm

## 10. Warranty / Service

Further information on the guarantee and guarantee conditions can be found in the guarantee leaflet supplied.

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Beurer UK Ltd., Suite 16, Stonecross Place, Yew Tree Way
Hans Dinslage GmbH, Riedlinger Str. 28,