beurer

stress releaZer®

DE	stress releaZer Gebrauchsanweisung.....	2
EN	stress releaZer Instructions for use.....	15
FR	stress releaZer Mode d'emploi	27
ES	stress releaZer Instrucciones de uso	40
IT	stress releaZer Istruzioni per l'uso	53
TR	stress releaZer Kullanım kılavuzu.....	66
RU	stress releaZer Инструкция по применению	79
PL	stress releaZer Instrukcja obsługi	92

CE

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89231 Neu-Ulm

Germany

Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Dear customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, gentle therapy, blood pressure/diagnosis, weight, massage, beauty, air and baby-care.

With kind regards,
Your Beurer team

Table of contents

1. Getting to know your device	17	6.1 Positioning the stress releaZer ...	23
1.1 How does the stress releaZer work?	17	6.2 Switching on the stress releaZer/ selecting the breathing cycle	23
1.2 "beurer CalmDown" app.....	18	6.3 Heat function	24
2. Intended use	18	6.4 Light function	24
3. Notes.....	19	7. Cleaning and maintenance	24
4. Device description	21	8. Disposal	25
5. Initial use.....	22	9. Technical data	25
6. Use	23	10. Warranty / service	26

Included in delivery

- 1x stress releaZer
- 1x micro USB charging cable
- 1x mains adapter

WARNING

- The device is intended only for domestic/private use, not for commercial use.
- This device can be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and are aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.

Signs and symbols

The following symbols appear in these instructions for use:

	WARNING Warning notice indicating a risk of injury or damage to health
	IMPORTANT Safety note indicating possible damage to the device/accessory
	Note Note on important information

	Observe the instructions for use
	Device not to be used on individuals with a pacemaker or other implants, such as an insulin pump or an implanted defibrillator
	The device is double protected and therefore corresponds to protection class 2.
	Only use indoors.
	Manufacturer
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Dispose of packaging in an environmentally friendly manner

1. Getting to know your device

1.1 How does the stress releaZer work?

The beurer stress releaZer was specially developed for the prevention and treatment of stress. By covering the diaphragm area with the stress releaZer, primary nerve pathways in this region are relaxed and your natural breathing rhythm is restored.

The built-in vibration motor transmits pulses to the user, which allows them to breathe normally. The stress releaZer also has an illuminated ring, which functions as a mood light. While the device is in use, the illuminated ring pulsates in time with the user's breathing rhythm.

The stress releaZer also has an optional heat function on the underside for a more comfortable application.

1.2 “beurer CalmDown” app

You can also connect the stress releaZer to your smartphone using *Bluetooth*[®]. To do this, all you need to do is download the free “beurer CalmDown” app from the Apple App Store (iOS) or from Google Play (Android[™]).

The free “beurer CalmDown” app offers the following advantages:

- Detailed instructions for use with videos
- Audiovisual breathing exercises with an HRV (heart rate variability) training programme
- Relaxing music to accompany use
- A calendar function for recording periods of use
- An option to add personal notes
- System requirements: 
 - iOS ≥ 10.0
 - Android[™] ≥ 5.0
 - Bluetooth*[®] ≥ 4.0
- Data transfer via *Bluetooth*[®] wireless technology.

The stress releaZer[®] uses *Bluetooth*[®] low energy technology, frequency band 2.400 - 2.480 GHz, the maximum transmission power emitted in the frequency band < 20 dBm, compatible with *Bluetooth*[®] 4.0 smartphones/tablets.

Technical information is subject to change without notification to allow for updates.

2. Intended use

WARNING

The stress releaZer is intended for private use only. The stress releaZer is not a medical device, rather it is a relaxation device. The stress releaZer may only be used for the purpose for which it is designed and in the manner specified in the instructions for use. Any form of improper use can be dangerous. The manufacturer is not liable for damage resulting from improper or careless use.

3. Notes

Safety note regarding the use of the stress releaZer

- Keep the packaging material away from children. choking hazard!
- Before use, ensure that the device and accessories have no visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.
- Keep the mains adapter away from water.
- Do NOT use the device if you have a pacemaker or other implants, such as an insulin pump or an implanted defibrillator.
- Do NOT use the device on the head.
- Do NOT use the device on the face.
- Do NOT use the device on the genitals.
- Do NOT use the device if you suffer from a medical abnormality or injury in the diaphragm area (e.g. open wound).
- Do NOT use the device during pregnancy.
- Do NOT use the device while sleeping.
- Do NOT use the device in a vehicle.
- Do NOT use the device on animals.
- Do NOT use the device whilst undertaking any activity where an unexpected reaction could be dangerous.
- Do NOT use the device after the consumption of substances that lead to impaired cognitive function (e.g. pain medication, alcohol).
- The device features an optional heat function. Persons not sensitive to heat must be careful when using the heat function.
- Use of the device must feel pleasant at all times.
- Only use the device with the supplied accessories.
- Repairs must only be carried out by Customer Services or authorised retailers. Under no circumstances should you open or repair the device yourself, as faultless functionality could no longer be guaranteed thereafter. Failure to comply with this instruction will void the warranty. Never attempt to repair the device yourself!

If you are unsure whether the device is suitable for you, please consult your doctor.


Please consult your doctor before using the stress releaZer:

- If you suffer from a severe illness or have recently had surgery on your upper body.
- In the case of feverish colds, thrombosis, diabetes, neurological diseases or severe inflammations.
- If you are suffering pain, the cause of which is unknown.

Notes on handling batteries

- If your skin or eyes come into contact with fluid from a rechargeable battery cell, flush out the affected areas with water and seek medical assistance.
- Choking hazard! Small children may swallow and choke on rechargeable batteries. Store rechargeable batteries out of the reach of small children.
- Protect batteries from excessive heat.
- Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, split or crush the rechargeable batteries.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use.
- In order to achieve as long a battery service life as possible, fully charge the battery at least twice per year.

4. Device description

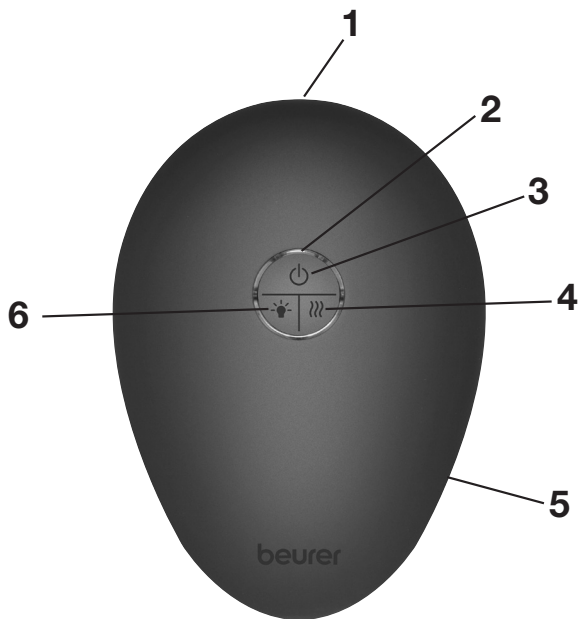1 Micro USB connection

4 Heat button

2 Function light

5 Illuminated ring

3 ON/OFF button ⏻

6 Illuminated ring button

5. Initial use


Note

Before using the stress releaZer for the first time, let it charge for a minimum of 4 hours. Proceed as follows:

Charging the stress releaZer

- Connect the micro USB charging cable to the mains adapter and the stress releaZer (Figure 1). Then insert the mains adapter into a suitable socket. Alternatively, you can also charge the stress releaZer with the charging cable via the USB port on your computer or laptop (Figure 2).

When the device is charging, the function light will light up alternately in a circle. When the stress releaZer is fully charged, the function light will be permanently illuminated. The battery will last for approx. 3 hours when fully charged (with heat and light function activated) or approx. 11 hours without heat and light function. Charge the stress releaZer again as soon as the function light flashes blue during operation.

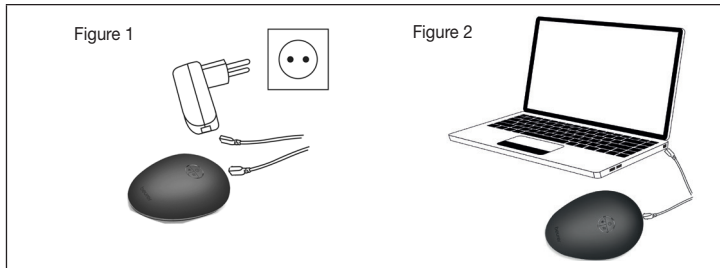Further information and FAQs available at:
<http://www.beurer.com>


6. Use

6.1 Positioning the stress releaZer

i Note

You can use the stress releaZer either on bare skin or over clothing. We recommend that you use the stress releaZer for 30 minutes every day.

1. Lie flat on your back. Alternatively, you can also sit on a reclined chair.
2. Place the stress releaZer on the diaphragm area (resting on the bottom rib bones/chest). If you find that use directly on the diaphragm is uncomfortable, position the stress releaZer above the diaphragm.

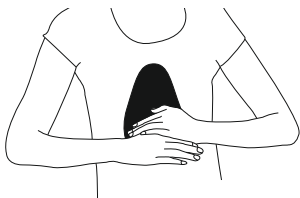i Note

Do not use the device near your solar plexus (area between your chest and your stomach). This may cause dizziness and nausea.

6.2 Switching on the stress releaZer/selecting the breathing cycle

i Note

Select the breathing cycle of your choice. If the stress releaZer is used over an extended period of time, it is possible to increase your breath duration.

Press  once	=	Press the ON/OFF button  to switch on the stress releaZer. <i>Bluetooth</i> [®] is activated when three sections of the function light light up blue. You can now connect the stress releaZer to your smartphone using the “beurer CalmDown” app if you wish.
Press  twice	=	The first breathing cycle is activated. In this breathing cycle, breathe in for 5 seconds, and then breathe out for 5 seconds.
Press  three times	=	The second breathing cycle is activated. In this breathing cycle, breathe in for 6 seconds, and then breathe out for 6 seconds.

Press  four times	=	The third breathing cycle is activated. In this breathing cycle, breathe in for 7 seconds, and then breathe out for 7 seconds.
Press  five times	=	The stress releaZer is switched off.


Note

You can also switch off the stress releaZer, regardless of the setting, by pressing and holding down the ON/OFF button.

Note

The stress releaZer features an automatic switch-off function, which is set to the maximum operating duration of 30 minutes. After 30 minutes, the stress releaZer switches off automatically.

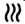6.3 Heat function

You can switch on a soothing heat during use to increase relaxation. To do this, press the heat button  on the stress releaZer. The function light starts to glow blue. Press the heat button  again to switch off the heat.

6.4 Light function

You can also switch on a blue light during use to increase relaxation. To do this, press the light button  on the stress releaZer. The illuminated ring around the stress releaZer will now start to pulse in time with the user's breathing rhythm. Press the light button  again to switch off the blue light. When the light function and the heat function are both switched on, the illuminated ring will glow red, not blue.

7. Cleaning and maintenance

The service life of the stress releaZer depends on careful use.

IMPORTANT

- Clean the stress releaZer after every use using a soft, slightly damp cloth (water or a mild cleaning solution). Dry the stress releaZer using a lint-free cloth. Do not use any abrasive cleaning products.

- Protect the stress releaZer from knocks, dust, chemicals, drastic changes in temperature, electromagnetic fields and nearby sources of heat (ovens, heaters).
- Do not submerge the stress releaZer in water or any other fluids and make sure that no fluids can get into the stress releaZer.

8. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life.

Dispose of the device at a suitable local collection or recycling point in your country. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

Rechargeable battery disposal

The empty, completely flat rechargeable batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.

To remove the battery from the device, pry open the cover above the illuminated ring using a screwdriver. Remove the battery from inside the housing.

The codes below are printed on rechargeable batteries containing harmful substances:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

9. Technical data

Supply	100-240V~; 50/60Hz; 0.2A
Dimensions	18 x 12.8 x 5 cm
Weight	Approx. 380g
Battery: Capacity	2200mAh
Nominal voltage	3.7V
Type designation	Lithium-ion

Data transfer	<i>Bluetooth</i> [®] low energy technology, 2.400 - 2.480 GHz frequency band transmission power < 20 dBm 
---------------	---

Technical information is subject to change without notification to allow for updates.

We hereby guarantee that this product complies with the European RED Directive 2014/53/EU. Please contact kd@beurer.de to obtain detailed information, such as the CE Declaration of Conformity.

10. Warranty / service

In case of a claim under the warranty please contact your local dealer or the local representation which is mentioned in the list “service international”. In case of returning the unit please add a copy of your receipt and a short report of the defect.

The following warranty terms shall apply:

1. The warranty period for BEURER products is either 3 years or- if longer- the country specific warranty period from date of Purchase.
In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Repairs (complete unit or parts of the unit) do not extend the warranty period.
3. The warranty shall not be valid for damages because of
 - a. improper treatment, e.g. nonobservance of the user instructions.
 - b. repairs or tampering by the customer or unauthorised third parties.
 - c. transport from the manufacturer to the consumer or during transport to the service centre.
 - d. The warranty shall not be valid for accessories which are subject to normal wear and tear (cuff, batteries etc.).
4. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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The *Bluetooth*[®] word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Beurer GmbH is under license. Other trademarks and trade names are those of their respective owners.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC.

Android is a trademark of Google LLC.

Beurer GmbH • Söflinger Str. 218 • 89077 Ulm, Germany
www.beurer.com • www.beurer-healthguide.com
www.beurer-gesundheitsratgeber.com