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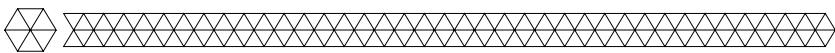
## Warnings and Safety Instructions

Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!

- The Acupressure mat is NOT a substitute for treatments administered by doctors or physiotherapists.
- Do NOT use the Acupressure Mat, or consult your physician if you are pregnant, have a skin disease or hemophilia or if you take anticoagulants.
- The Acupressure Mat should NOT be used by children under the age of 5 or persons with physical, sensory or mental disabilities, unless they are supervised or have been given instruction concerning the use by a person responsible for their safety.

## Contents of Package

- 1x LifeGoods Acupressure Mat and Pillow; 1x User Manual; 1x Instruction manual; 2x massage ball; 1x Carry-on Bag

## Product Specifications

Size	66,00 x 40,00 x 2,00 CM	Pads	Nr. of spikes: 27
Weight	0.685 KG	Nr. of pads	Mat: 230   Pillow: 66
Material	100% cotton, foam, ABS		

## Maintenance - Cleaning - Storage

- The Acupressure Mat is hand wash only. Use lukewarm water and air dry. Store the mat in a clean and dry place. Do not store the device in extreme high or low temperatures.

## Disposal - Recycle


- Contribute to a cleaner living environment! Do not dispose of the mat with the household waste. For proper recycling, please contact your local authority or your household waste disposal service for further details on your nearest designated collection point.

## Service and Warranty

- If you need service or information regarding your product, please contact our Customer Support at [service@lifegoods.nl](mailto:service@lifegoods.nl)
- LifeGoods provides a 2 year warranty on its products. To obtain service during the warranty period, the product needs to be returned with a proof of purchase. Product defects have to be reported within 2 years from the date of purchase. This warranty does not cover: damages caused by misuse or incorrect repair; parts subject to wear; defects that the customer was aware of at the time of purchase; damages or defects caused by customer neglect; damages or defects caused by third parties.

## Disclaimer

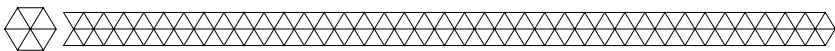- Changes reserved; specifications are subject to change without stating grounds for doing so.

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© LifeGoods B.V.