TH930

NL. CONTACTLOZE VOORHOOFD /
OORTHERMOMETER
FR. THERMOMÈTRE SANS CONTACT
FRONTAL ET AURICULAIRE
DE. BERÜHRUNGSLOSES STIRN /
OHR THERMOMETER
EN. NON-CONTACT FORHEAD /
EARTHERMOMETER

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MODEL: PG-IRT1603

INSTRUCTION MANUAL

We thank you for buying our C-CARE TH950 digital Ear/Forehead thermometer. It shows your confidence in our brand and our products. Please be assured you have bought a high quality product which you will enjoy using for many years to come.

Please read this instruction manual thoroughly before first use, even if you are familiar with this type of product. It gives you tips for optimal use.

If you require technical support and/or have questions about the functioning of the product please send an e-mail to our after sales service: the address is listed here under. We will give you a fast reply.

CRESTA INTERNATIONAL B.V.
ROLBURGWEG 4
NL-1332 AS ALMERE
The Netherlands
info@cresta-international.com
www.cresta-care.com
Tel: + 31 20 6932932

CONTENTS

Safety warnings

Product Introduction

Safety instructions

About body temperature

Product layout

Display layout

Battery installation

Basic instructions

Introduction of measurement methods

Cleaning

Frequently asked questions and solutions

Troubleshooting

Product specifications

SAFETY WARNINGS

The warning signs and illustrations in this manual are intended to enable you to use the product safely and correctly and to prevent any harm to you and others. Please find the several signs, illustrations and their meaning listed here under:

⚠ ATTENTION: Here under the signs and their corresponding meaning.

Legends	
	This mark means a warning, instructions you have to follow
	This mark means a general warning.
	This mark means prohibited things. Things that are not allowed to do.
	This mark means "prohibited to disassemble"
	Type BF Applied person
	This mark means Water Resistance, Moisture Resistance
	This mark means "Read the manual carefully before use"
	This mark means that the packaging of this unit can comply with the requirements of a green environmental protection
	This mark means that the material of the product or the product itself is made of renewable material so that we can cycle this product, which benefits the environment and our earth.
	This mark means "Do not remove with the normal municipal waste"
	IP CLASSIFICATION: IP22

WARNING

- Please ask professional doctors to explain the measured value of your body temperature.
- The product is only to be used for the measurement of the human body temperature.
- Do not use this product for purposes other than body temperature measurement.
- Do not use mobile phones around this product.
- Please do not use any equipment that generates electromagnetic fields near the product.

- Do not disassemble or repair this device. Do not change the product. Your guarantee will be void.
- Do not bend or stretch the product.
- Please do not touch or drop the product.

SAFETY INSTRUCTIONS

- This unit is only to be used for the measurement of the human body temperature without the use of a disease diagnosis; It cannot be used for emergency and continuous measurement during surgery.
- Please keep children away from this product. For children ranging from zero to twelve years old they may only use this product under supervision of their parents.
- Patients cannot self diagnose the disease and get treatment by themselves on the base of the measurement results. Contact your doctor and follow his/her instructions.
- Children under 12 years of age and people who are unable to express themselves are prohibited to use this product without supervision.
- If you use or store this product beyond the range of the specified temperatures and humidities the product may not reach the original performance specifications.

Environment temperature	:	from -10 °C - +40 °C,
Humidity	:	from 15% RH - 93% RH
Storage temperature	:	from -25 °C - +55 °C,
Humidity	:	from 0% RH - 93% RH

PRODUCT INTRODUCTION

Intended use: Infrared Ear/Forehead Thermometer intended to measure the human body temperature by measuring in the ear canal or on the forehead.

Scope of application: It is suitable to display the measured body temperature by measuring the heat radiation in the ear canal or on the forehead.

Features:

- Non-touching type of infrared measurement.
- Multi-color and backlight display: White, Green, Orange and Red.
- Memory for 9 measurements.
- Measuring unit selectable in degree Celsius °C and degree Fahrenheit °F. The original factory setting is degree Celsius °C.
- Instant measurement result within 1 second.
- The design is convenient and economical without the use of any accessory, which saves subsequent using costs.
- Sound support can be switched on or off.
- Automatic shut down after 30 seconds of no operation.

ATTENTION:

The measurement results of this unit are at any time for reference only. The result can not replace the medical diagnosis of a professional doctor. If you have any questions about your personal measurement result please contact your doctor for a professional diagnosis and treatment.

SAFETY MEASURES

WARNING

- It is dangerous for patients to judge the measurement result by themselves and adjust a treatment. Always contact your doctor to determine any treatment and follow his / her instructions.
- Self-judgment may lead to worsening of the patients condition.
- Do not touch the infrared sensor with your hands or blow the sensor "clean" with your mouth.

- If the infrared sensor is damaged or dirty this may lead to abnormal measurement results.
- If there is a temperature difference between the storage site and the measuring site, please place the thermometer in the measuring site to adjust to the temperature for about 30 minutes before you measure.
- Large temperature differences may lead to incorrect results.
- Please keep this product out of reach of children.
- Pay attention to small parts. If children by accident swallow such a part, like a battery contact, your doctor immediately.
- Do not measure your temperature near an airconditioning. The airconditioning influences the measuring results.
- Use a cotton swab dipped in 75% alcohol each time you use the thermometer to clean the infrared sensor surface so the sensor glass is dirt free. When you see stains, a small glass lens surface or water on the sensor clean the sensor glass carefully.
- Do not wipe the infrared sensor with toilet paper, kitchen paper or facial tissues as this can cause scratches on the sensor lens and influence the measurement results.
- If the product is damaged the measurement results may be incorrect.
- If the product accidentally touches water or is immersed into water please fully dry it before you re-use the product. Specially clean the water on the sensor lens with a cotton swab. This prevents that water that was left behind influences the measurement results.
- We aim to avoid causing safety accidents and affecting the measurement accuracy.

ATTENTION:

- Please do not use this product for people who suffer from otitis externa, tympanitis and other ear diseases. It is possible to worsen the affected area.
- Do not use the product immediately after swimming, bathing or with wet ears. The measured temperature may be lower than in reality.
- Do not place a waste battery in a danger zone, the battery may break. Put it in a safe place before taking it to the chemicals depot of your municipality. Batteries may leak.
- When measuring the temperature in the ear canal of the human body the thermometer should be operating in the ear temperature mode, hence the readings may be incorrect.

Suggestions

- When you inform your doctor about the measured body temperature explain you have measured with an ear thermometer.
- Do not force the thermometer to impact, fall, trample and shake the product.
- Do not disassemble, repair or modify the product, your guarantee will be void.

- Please do not allow liquids such as alcohol, water drops, hot water and so on to enter the product body as this is not water resistant.
- Store the product in a clean, dry place.
- If you have any problem with the product please contact your store or online company where you bought the product. Do not repair the product yourself, it will make your guarantee void.
- Please do not use the product in an environment with electromagnetic interference such as near a microwave, PC, wireless telephone etc. Please deal with the waste and residue of this product at the end of its life cycle according to the local laws and regulations. Contact your municipality for your specific local regulations.

ABOUT BODY TEMPERATURE

The comparison of different measurement methods.

The measured values are different if we use different measurement methods. The World Health Organization (WHO) provides normal human body temperature reference values, see the table here under, for the specific temperature differences.

MEASUREMENT METHOD	NORMAL BODY TEMPERATURE
Anal temperature	36,6°C -38,0°C
Oral temperature	35,5°C -37,5°C
Axillary temperature	34,7°C -37,3°C
Cochlear temperature	35,8°C -38,0°C
Cochlear temperature Measured with TH730	35,5°C -37,5°C

The changes in the human body temperature.

Humans belong to the creatures with a constant temperature. The body temperature is basically constant but this does not mean it's totally changeless. The human body temperature is constantly changing during the day. The details are as follows.

At night:

Lowest. The body temperature is lowest during your sleep and decreased activities. (Under 37 °C)

In the morning

Higher. From your warm bed to the lower temperature room will make your whole bodies muscles get contractions and produce heat.


At noon

Highest. After lunch the human body reaches its highest temperature and the body will adjust naturally.


Three or four o'clock in the afternoon.
Lower. Due to physical exertion and low blood sugar.

In the evening
Lowest. Due to the sunset the environmental room temperature goes down.

PRODUCT LAYOUT

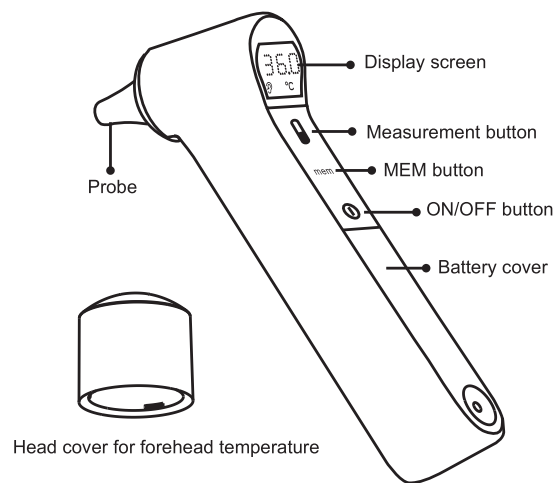DISPLAY LAYOUT

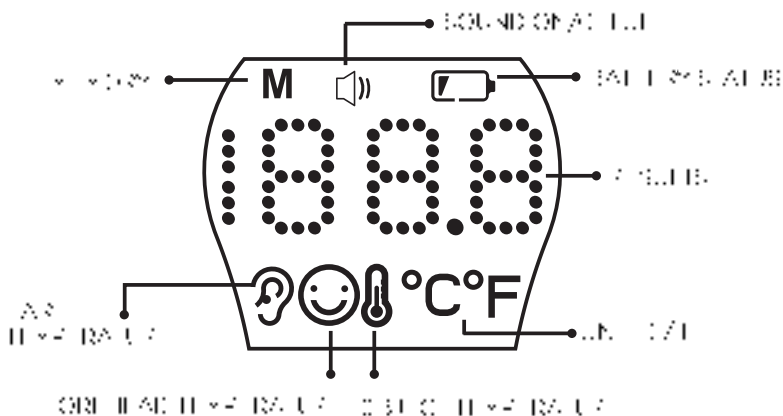
BATTERY INSTALLATION


1. Press the battery cover, it will bounce open automatically
2. Put 2 pieces 1,5 V AAA batteries carefully watching the correct + and - polarity

We recommend the use of alkaline batteries.

Press the battery cover to open, it will bounce automatically

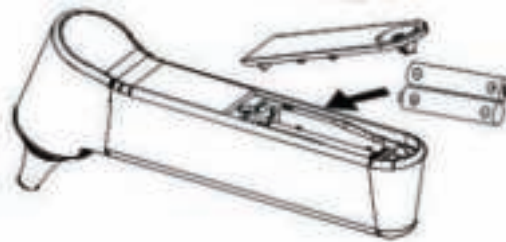Press with your fingers the poles of the batteries to remove them

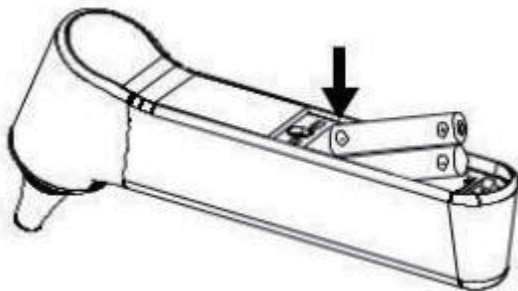LOW BATTERY POWER

When the battery power is too low "LO" will appear in the display and the battery mark will constantly be on. This is the sign that the batteries need to be replaced.

ATTENTION

If you intend not to use this product for a longer period we advise to remove the batteries to assure the longevity and prevent battery leakage. The leakage liquid may cause damage to the product and may pollute the environment.

For information how to dispose of empty batteries please contact your municipality for the local regulations.

BASIC INSTRUCTIONS

Sound function on / off

Press in the power "ON" status the "MEM" key to set the sound "ON" or "OFF".


Press the "MEM" - key, the Led-display shows "🔊". This mark means the sound function is "ON". You will hear a short beep sound. Press the "MEM" key once more and "🔊" will disappear. This means the sound function is "OFF".

° Setting the measurement unit: Celsius or Fahrenheit.

In the "OFF" status long press – around 6 seconds – the "MEM" - key. You can now switch between degree Fahrenheit (° F) and degree Celsius (° C). Wait for 8 seconds to turn on the product automatically or press "🔘" to turn off the product directly.

Memory Storage Function

In the "OFF" status press the "MEM" - key. The product can save and read 9 measurement values in order – see picture below. The Thermometer will turn off automatically after 30 seconds without operating or press "🔘" to turn off the product manually.

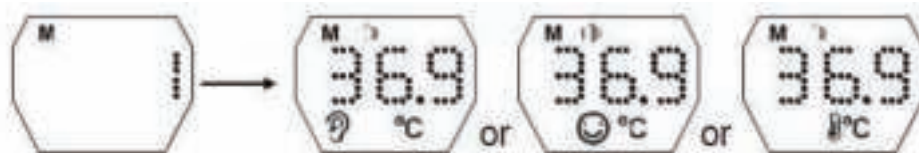BACK LIGHT STATUS INSTRUCTIONS

When the measured temperature is lower than 34 ° C, the display shows "LO" with red backlight.

When the measured temperature is between 34 ° C and 37,1 ° C, the display will show the measured temperature with green backlight.

When the measured temperature is between 37,2 ° C and 38,1 ° C, the display will show the measured temperature with orange backlight.


When the measured temperature is between 38,2 ° C and 43,0 ° C, the display will show the measured temperature with red backlight.

When the measured temperature is higher than 43°C, the display shows "HI" with red backlight.

ATTENTION: This function is for reference only.


C CARE

INTRODUCTION OF THE MEASUREMENT METHODS

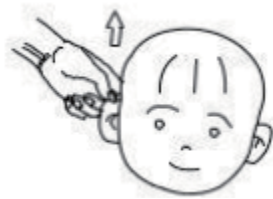POSITION WITH CHILDREN WITH UNDER 1 YEAR OLD

EAR TEMPERATURE

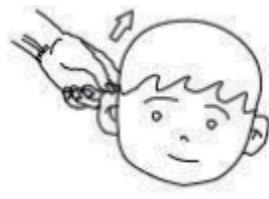POSITION WITH CHILDREN OVER 1 YEAR AND ADULTS


EAR TEMPERATURE

FOREHEAD

FOREHEAD TEMPERATURE

EAR TEMPERATURE MEASUREMENT

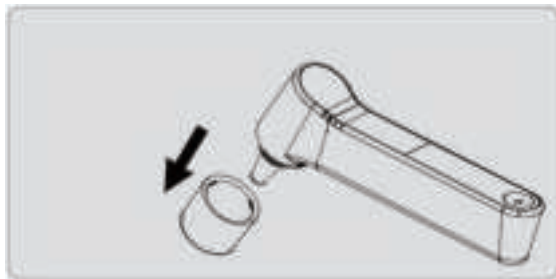
REMOVE THE HEAD COVER OF THE FOREHEAD THERMOMETER TO MEASURE THE EAR TEMPERATURE

After removing the head cover press the "0" - button to turn the thermometer on and to enter the ear temperature measurement mode. In the LCD display the last measurement result will appear. Then put the thermometer in the ear and press the probe thoroughly against the eardrum. The temperature can be measured immediately by pressing the "M" button. After 1 second you can see the measurement result.

ATTENTION

IF YOU DO NOT PRESS THE SENSOR DEEP ENOUGH INTO THE EAR CANAL AGAINST THE EARDRUM YOU WILL GET AN INCORRECT - OFTEN TOO LOW - MEASUREMENT RESULT.

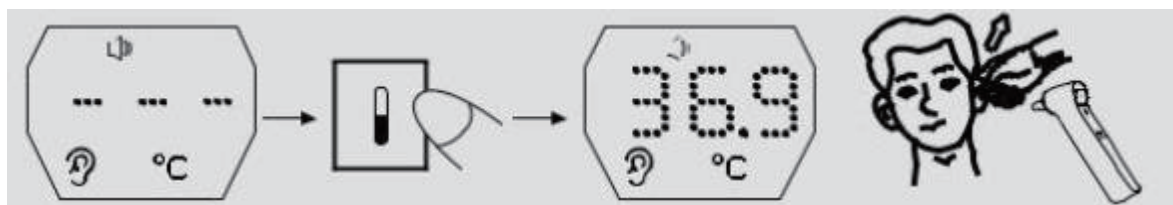
When the LCD display is flashing we are waiting for the measurement to start.

Measurement result

Attention: If you did not hear the beep sound, which represents the temperature measurement, the measuring is not completed. Please do not move the thermometer from the ear until you hear the beep sound. If you have switched off the sound indication you will not hear the beep sound.

FOREHEAD TEMPERATURE MEASUREMENT

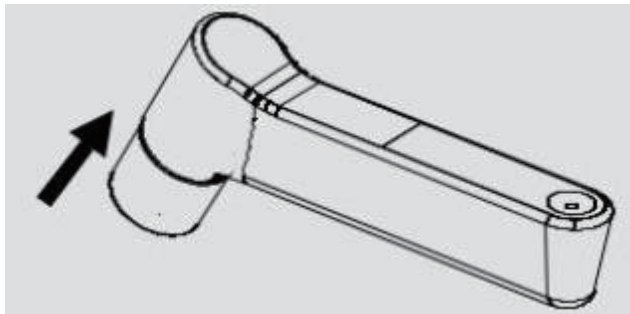Install the head cover of the forehead thermometer so you can measure the forehead temperature. After installing the head cover press the "Ⓚ" button to turn on the thermometer and to enter the forehead temperature measurement mode. The Lcd display shows the last reading. Then aim the forehead thermometer at the forehead with a measurement distance of 0 - 5 mm. The forehead temperature can be measured directly by pressing the "Ⓚ" button.

After 1 second you can see the measurement result.

When the lcd display is flashing we are waiting for the measurement to start.

Measurement result


Attention: If you did not hear the beep sound, which represents the temperature measurement, the measuring is not completed. Please do not move the thermometer from the forehead until you hear the beep sound. If you have switched off the sound indication you will not hear the beep sound.

OBJECT PATTERN MEASUREMENT

The installation of the head cover of the forehead thermometer does not influence the object pattern measurement.

Long press the "M" button to enter the object measurement mode for 5 seconds, and then aim the thermometer at the object. The temperature of the target object can be measured directly by pressing the "M" button.

After 1 second you can see the measurement result.

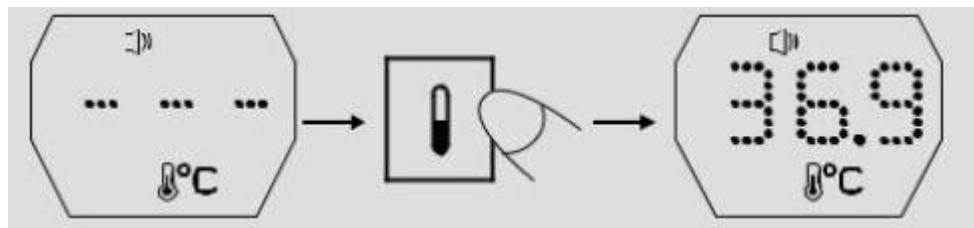When the LCD display is flashing '...' and '°C', we are waiting for the measurement to start.


Measurement result

Attention: If you did not hear the beep sound, which represents the temperature measurement, the measuring is not completed. Please do not move the thermometer from the forehead until you hear the beep sound. If you have switched off the sound indication, you will not hear the beep sound.

CLEANING

In order to ensure the precise measurements, we recommend to clean the device after each use.

Please use a cotton swab and dip into 75% alcohol to wipe the thermometer probe and remove any residue and dirt. We recommend to clean the probe at least 5 seconds. Then wipe the other components with a soft, dry cloth for at least 10 seconds. Please make sure that after the cleaning, no visible dirt or residue is left on your thermometer. Wait 10 minutes before measuring.

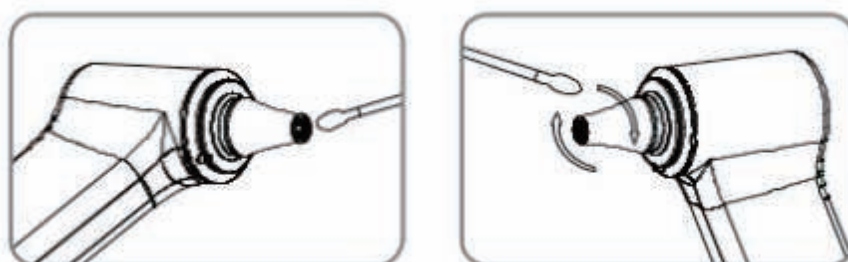ATTENTION:

Why do you need to clean the device after each use?

As the infrared temperature adopts a highly sensitive technique to detect the temperature of the target object any dust or skin will affect the accuracy of the measurement. We therefore recommend to clean the probe after each use as shown in the picture.

FREQUENTLY ASKED QUESTIONS AND SOLUTIONS

DISPLAY	REASON	SOLUTION
HI	When the measured temperature is higher than the measurement range "HI" is shown in the display. The range for the ear / forehead mode is $< 3^{\circ}\text{C}$, and $93,2^{\circ}\text{C}$ for the object temperature mode.	
LO	When the measured temperature is lower than the measurement range "LO" is shown in the display. The range for the ear / forehead mode is 34°C , and 0°C for the object temperature mode.	When the measurement probe is not properly into the ear canal / ear drum or the measurement distance is too far during the measurement process the measured result may be low. When the measurement probe is dirty the measured result may be low. Please clean the probe properly with a cotton swab dipped in 75% alcohol.
Er.H	The upper limit of the operating environment temperature for this product is 40°C . When the environment temperature exceeds this limit "Er.H" is shown in the display.	Make sure you measure in an environment with a temperature lower than 40°C .
Er.L	The lower limit of the operating environment temperature for this product is 10°C . When the environment temperature exceeds this limit Er.L is shown in the display.	Make sure you measure in an environment with a temperature higher than 10°C .
Err	When the environment temperature changes rapidly about 5 degrees the display will show "Err" during the measurement in object measurement mode and will then shutdown automatically.	When "Err" is displayed please place the product and keep it steady in the environment for 30 minutes so the product can have a constant temperature before the next measurement.

When the operation is wrong and the sound is turned on, the value will be read and you will be reminded about the error by 3 short beeps sounds.

TROUBLESHOOTING

PROBLEM	REASON	SOLUTION
No display when the power is switched on	<ul style="list-style-type: none"> The batteries are empty The battery polarity is incorrect 	<ul style="list-style-type: none"> Replace the new batteries Put the batteries with the correct polarity as shown in the battery case.
The measured temperature is LOW	<ul style="list-style-type: none"> The measurement position is incorrect. The probe is dirty 	<ul style="list-style-type: none"> Measure again according to the instructions Clean the probe as described earlier and measure again
Big temperature fluctuations with continuous measurement	The measurement interval is too short	The interval between 2 measurements should be at least 10 seconds

PRODUCT SPECIFICATIONS

Manufactured for	
Name	: CRESTA INTERNATIONAL BV Rohbrugweg 4 - 1352 AB ALMERE The Netherlands
Product name	: Infrared Ear/Forehead thermometer
Model number	: TH930
Product dimensions	: 31 x 175 x 72 mm
Product weight	: approx. 77 gram - except. batteries
Measuring range	: 34.0°C - 43°C / 93.2°F - 109.4°F
Object temperature	: 0°C - 93.2°C / 32°F - 199.7°F
Resolution/res.:	: 0.1°C / °F
Measurement location	: Laboratory
Accuracy	: 35.0°C - 42.0°C \pm 0.2°C 95.0°F - 107.9°F \pm 0.4°F Other temperatures \pm 0.3°C
Operation temperature	: 10.0°C - 40.0°C / 50.0°F - 104.0°F
Relative humidity	: Maximum 15% RH - 93% RH
Atmospheric pressure	: 70 kPa - 106 kPa
Transportation / storage temperature	: -25°C - 55°C / -13°F - 131°F
Relative humidity	: Maximum 50 kPa - 106 kPa
Display screen	: Lcd display, 4 bits numbers and special icons
Sound	: Short beep with turning on and ready to measure, long beep when measurement ready, 3 short beeps with system error, 10 short beeps with fever alarm
Memory	: Memory for 9 measurements
Auto shut off	: After 30 seconds of no operation
Batteries	: 2 x 1.5 V AAA batteries. The use of alkaline batteries recommended
Period of use	: Minimum 5 years
ISO 9001 Forehead mode	: Clinical bias Dcb: 0.078 Limits of Agreement LA: 0.243 Clinical repeatability: 00.069 Reference body site: forehead Measuring site: forehead
Packing parts list	: Thermometer Product manual

Appendix 1 Guidance and Manufacturer Declaration Tables

Guidance and manufacturer's declaration–electromagnetic emissions		
The Model PG-IRT1603 Infrared Thermometer is intended for use in the electromagnetic environment specified below. The customer or the user of the Model PG-IRT1603 Infrared Thermometer should assure that it is used in such an environment.		
Emissions	Compliance	Electromagnetic environment-guidance
RF emissions CISPR 11	Group 1	The Model PG-IRT1603 Infrared Thermometer uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The Model PG-IRT1603 Infrared Thermometer is used in home and it's powered by DC 3V
Harmonic emissions IEC 61000-3-2	N A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	N. A.	

Guidance and manufacturer's declaration – electromagnetic immunity

The Model PG-IRT1603 Infrared Thermometer is intended for use in the electromagnetic environment specified below. The customer or the user of the Model PG-IRT1603 Infrared Thermometer should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15KV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 KV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m, 50/60Hz	30 A/m, 50/60Hz	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE U_T is the a.c. mains voltage prior to application of the test level

Guidance and manufacturer's declaration – electromagnetic immunity			
The Model PG-IRT1603 Infrared Thermometer is intended for use in the electromagnetic environment specified below. The customer or the user of the Model PG-IRT1603 Infrared Thermometer should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz 6 Vrms 150 kHz to 80 MHz outside ISM bands	N/A	Portable and mobile RF communications equipment should be used no closer to any part of the Model PG-IRT1603 Infrared Thermometer, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = \frac{e \cdot \dot{u}}{r \cdot \sqrt{P}}$ $d = \frac{e \cdot \dot{u}}{r \cdot \sqrt{P}} \quad 80\text{MHz to } 800\text{MHz}$ $d = \frac{e \cdot \dot{u}}{r \cdot \sqrt{P}} \quad 800\text{MHz to } 2.7\text{GHz}$
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz	10 V/m	where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres(m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range ^b Interference may occur in the vicinity of equipment marked with the following symbol: 

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations.

Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a The ISM (industrial, scientific and medical) bands between 0,15 MHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz. The amateur radio bands between 0,15 MHz and 80 MHz are 1,8 MHz to 2,0 MHz, 3,5 MHz to 4,0 MHz, 5,3 MHz to 5,4 MHz, 7 MHz to 7,3 MHz, 10,1 MHz to 10,15 MHz, 14 MHz to 14,2 MHz, 18,07 MHz to 18,17 MHz, 21,0 MHz to 21,4 MHz, 24,89 MHz to 24,99 MHz, 28,0 MHz to 29,7 MHz and 50,0 MHz to 54,0 MHz.

b The compliance levels in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2,7 GHz are intended to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas. For this reason, an additional factor of 10/3 has been incorporated into the formulae used in calculating the recommended separation distance for transmitters in these frequency ranges.

c Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Model PG-IRT1603 Infrared Thermometer is used exceeds the applicable RF compliance level above, the Model PG-IRT1603 Infrared Thermometer should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Model PG-IRT1603 Infrared Thermometer.

d Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Recommended separation distances between portable and mobile RF communications equipment and the Model PG-IRT1603 Infrared Thermometer

The Model PG-IRT1603 Infrared Thermometer is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Model PG-IRT1603 Infrared Thermometer can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Model PG-IRT1603

Infrared Thermometer as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = \frac{3.5}{V_1} \sqrt{P}$	80 MHz to 800 MHz $d = \frac{3.5}{E_1} \sqrt{P}$	800 MHz to 2.7 GHz $d = \frac{7}{E_1} \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1. 2	2.3
10	3.8	3. 8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.
NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

ENVIRONMENT FRIENDLY DISPOSAL

The European Directive 2002/96/EC on Waste Electrical and Electronic Equipment - WEEE - requires that old household electrical appliances should not be disposed of in the normal unsorted municipal waste stream. Old appliances should be collected separately in order to optimize the recovery and recycling of the materials they contain and reduce the impact on human health and the environment.

The crossed out "wheeled bin" symbol on the product reminds you of your obligation that when you dispose of the appliance, it should be separately collected. Contact your local authority for information concerning the correct disposal of your old appliance.

PRINTING ERRORS AND MODIFICATIONS RESERVED

Europe - EU declaration of conformity

According to Medical Device Directive. Classification: (MDD. Annex IX): II a, Rule 10

UMDNS CODE: 17887 Thermometer Infrared, Ear
UMDNS CODE: 17888 Thermometer, infrared, Skin

Manufactured for:

Name: CRESTA INTERNATIONAL B.V.
Address Rolbrugweg 4 - 1332 AS Almere
The Netherlands

Equipment:
EAN number 8711974340024
Report reference PG-IRT1603
Marketing name: C CARE TH930
Class Class IIa, rule 10

Hereby, Cresta International BV declares that the above-mentioned products meet the provisions of the following EC council Directives and Standards and conformity with the Swedish legislation. All supporting documentations are retained under premises of the manufacturer and the notify body.

Directives:

Medical Device Directive: Council Directive 93/42/EEC Annex II, excluding Section 4

RoHS Directive (2011/65/EU)
REACH Directive (1907/2006/EC)

Standard:

All applied harmonized Standards (published in the Official Journal of the European Communities)

Certificate:

Registration : HD 60146953 0001

Manufacturer:

Shenzhen Pango Electronics Co., LTD
[Tel:8675533826988](tel:8675533826988) Fax8675533825989 Main Site No 25., 1ste Industry zone, Fenghuang Rd., Xikeng Village.
Henggang Town. Longgang District. Shenzhen China. Additional site 1 2-4 Floor no5, Shanzhuang RD. Xikeng
Village Henggan Town, Longgan District, Shenzhen City, Guangdong Province, China.

Notified Body:

TUV Rheinland LGA Products GMBH
Tillystrasse 2
90431 Nurnberg - Germany - CE 019

EC REP

Authorised representative: Lotus NL B.V.
Koningin Julianaplein 10, 2595AA The Hague, Netherlands
Tel: + 31 645171879

Signed for and on behalf of: CRESTA INTERNATIONAL B.V.
Place: Almere
Date: August 4th, 2020
Name: A. Benay
Title: Director