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for back and neck

Gray

Read all the instructions carefully before use and keep them for future reference. We reserve the right to make technical and optical changes without prior notice. Pictures and accessories/-functions may vary depending on the model. Some accessories/functions are optional. The use of attachments or accessories which aren't supplied with the product isn't recommended. This could result in electrical shock or injury or permanent damage to the product. This product is for household use only..

EK-01-NL

## **IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!**

Read the instruction manual carefully before using this product, especially the safety instructions, and keep the instruction manual for future reference. Should you give this device to another person, it's vital that you also give them these instructions.

Don't puncture the product!

Don't use the product when it's folded!

The product shouldn't be used by very young children (0-3 years old)!

Only use the product indoors!

The product can be washed by hand!

Don't bleach the product!

Don't dry the product in a drier!


Don't iron the product!

Don't dry clean the product!

## PRODUCT OVERVIEW

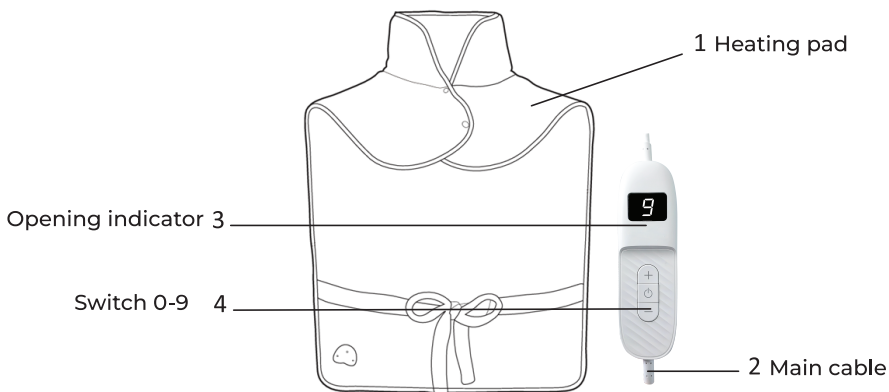automatic shutdown 90 min.

## FIRST TIME USE

Before using the product, inspect it for completeness and damages. Should the delivered contents be incomplete or damaged, contact your sales outlet immediately.

## OPERATION

- 1) Connect the power plug to a suitable power socket.
- 2) Switch on the product.

## TEMPERATURE SETTING

Easily set the temperature from low heat to maximum heat through 9 different heat settings.

## NOTE

Setting 1 is recommended for continuous use. The heating pad isn't designed for use in hospitals or for commercial purposes. It mustn't be used near babies, toddlers, heat sensitive people, disabled people or animals.

## AUTOMATIC SWITCH-OFF FUNCTION

The heating pad switches off automatically after approx. 90 minutes of continuous operation. Switch it off (setting 0) or unplug the power plug from the power socket if you're not going to use it after the automatic switch-off.

## CLEANING AND CARE

- Remove the cable and controller from the heating pad, and clean it according to the laundry care symbols on the label.
- Don't machine-wash, spin or machine dry, chemically clean, wrung out, press or iron the heating pad.
- Remove smaller stains with a cloth or moist sponge together (with a little liquid detergent for delicate fabrics, if necessary).
- Only use the heating pad again once it's completely dry.
- It's recommended to wash the heating pad only if it's dirty or dusty. If you wash the it too often, it could wear quickly.

## SPECIFICATION

Voltage: 220-2440V~50 Hz

Power: 100 W

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed, separately from your household waste, at the end of its lifespan. There are separate collection systems for recycling in the EU.

For more information, please contact us:

[klantenservice@lifeproductsstore.nl](mailto:klantenservice@lifeproductsstore.nl)

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Lifeproducts B.V.