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
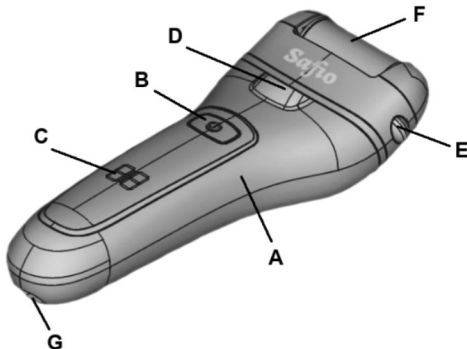
## SAFIO XV-7 CALLUS REMOVER INSTRUCTION

### CONTROL XV-7

- A. ABS body
- B. On-Off / 2-level speed switch
- C. Battery Capacity with LED Display
- D. LED lighting
- E. Roller head changing button
- F. Roller head
- G. Charging port
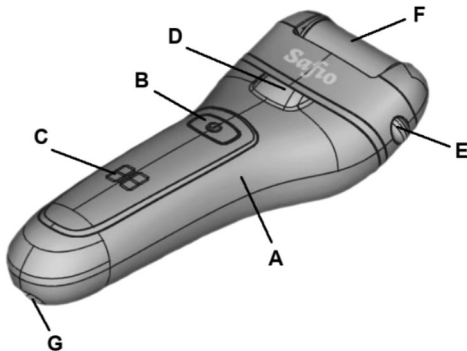
### CALLUS REMOVAL

This product is designed for dry and rough skin on hand, foot or heel to remove and should not be used on other parts of the body.

1. Ensure the product is enough charged.
2. Clean and soften the skin for about 10 minutes in warm water, then dry the skin to ensure the best results.
3. Gently turn the roller head every 2 or 3 seconds, to stay in the dry and rough skin area. Do not push too hard, let the roller head gently slide over the skin.
4. When the desired results are achieved, the skin can be rinsed, dried and rubbed with soothing cream.
5. Clean the device under running water and let it dry well. Store this product always in a dry cool place.

### 3 DIFFERENT ROLLER HEADS

**Rough:** Suitable for stubborn thick callus and dead skin (heel) foot

**Regular:** Suitable for the newly formed callus

**Fine** Suitable for fine callus like (hands) and normal skin care

## WARNING

Read the manual of the Safio XV-7 carefully before using the product and keep it for future reference.

1. This product cannot be used as toys for young children.
2. Children should use this product only under the guidance of an adult person.
3. This product cannot be used in spaces where explosive gases can be released.
4. This product cannot be disassembled, repaired or modified by yourself.
5. Do not store in a place exposed to direct sunlight, or in an environment that is too hot or humid.
6. Keep away from water when charging.
7. Do not submerge the product in water for a long time. Rinse the roller head under running water.
8. Only suitable for old callus, hard and thick dead skin to remove. Do not use on healthy skin.
9. For hygienic reasons, it is recommended to use one roller head per person.
10. Use this product after consulting your doctor if you have skin problems, or other medical problems.
11. Do not stay on one area for more than 3 seconds while using, make smooth movements over the calluses.
12. Do not use if the skin is red, irritated or infected, on ulcers or wounds on the skin.

### USE INSTRUCTIONS XV-7

#### CHARGING

This product needs to be fully charged before using it for the first time. Plug the cable in the device and a USB port of a laptop, smartphone charger or a USB connector. The LED display will light up while charging. After 2 hours charging the indicator will stop flashing when charging is completed. (charge time should not exceed 12 hours).

#### ROLLER HEAD DISSASSEMBLY

Replacing the roller head as follows. Choose the suitable roller head according to the roughness of the skin and the expected effect. To replace the roller head, ensure the device is turned off, the roller head will come loose by pressing on the roller head changing button on the right side of the device, Gently pull the roller head away from the device.

#### ROLLER HEAD ASSEMBLY

Place the new roller head by pressing on the roller head changing button on the right side of the device, you feel and hear the roller head click in the device. Before use check always if the roller head is well mounted.

### SPECIFICATIONS SAFIO XV-7

- Material: ABS body
- Charging voltage: 5V
- Battery: 1200 mAh
- Charge time: 2 hours
- Use time: 2,5 hours
- USB chargeable
- IPX7 waterproof

### ENVIRONMENT-FRIENDLY & RECYCLE

Broken electrical household devices can be returned for free with the purchase of a similar new product. Or hand in at various return points such as shops and hardware stores.

### AFTER SALE SERVICE

If there are any technical questions, you can contact our online customer service through the order page of your order.