# INFINITY GOODS MASSAGE GUN

PRODUCT HANDLEIDING

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---

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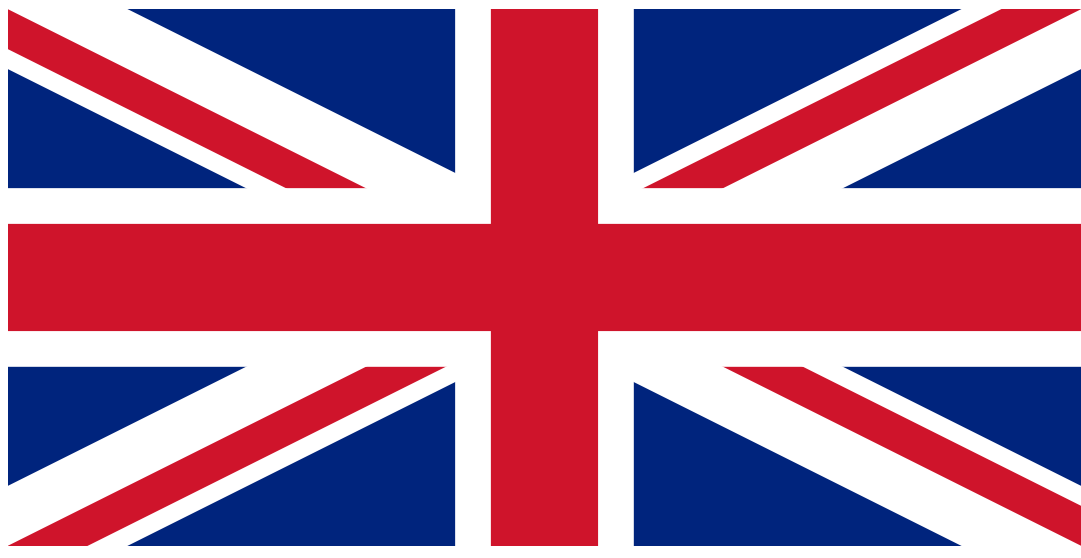
## **WARNINGS & CAUTION**

---

1. For adult use only, Do not use if injured. Consult your doctor before using this product.
2. Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
3. Do not remove screws or attempt to disassemble.
4. Do not operate continuously for more than one hour.
5. Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 second per region.
6. Use the massage gun only on the body's soft tissue as. desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body.
7. Keep fingers, hair or any other body part away from the shaft and backside of the applicator as pinching may occur
8. Do not place any objects into ventilation holes of Massage Gun

## **Disclaimer**

---

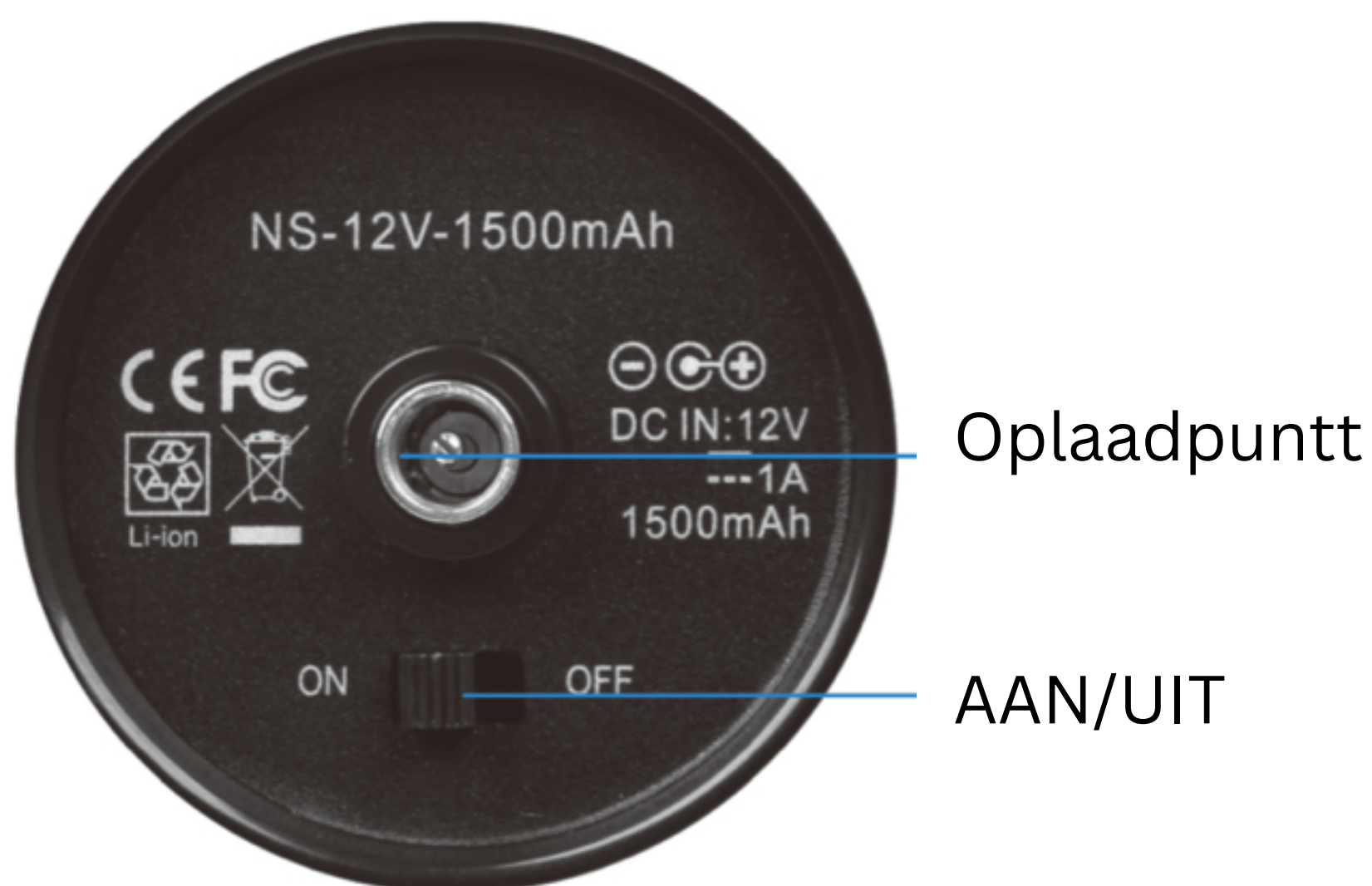
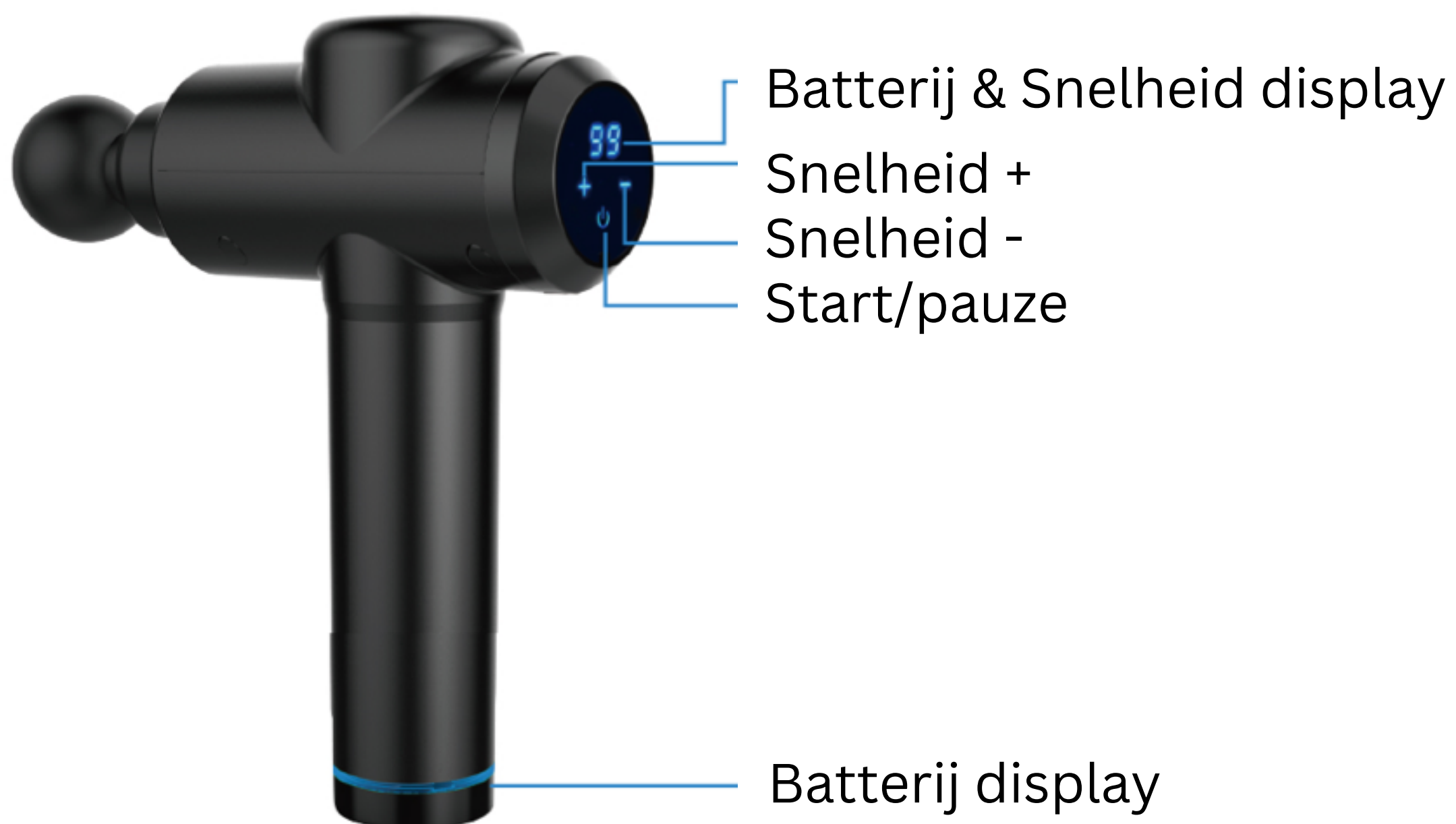
PLEASE DO NOT USE MASSAGE GUN WITHOUT APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY

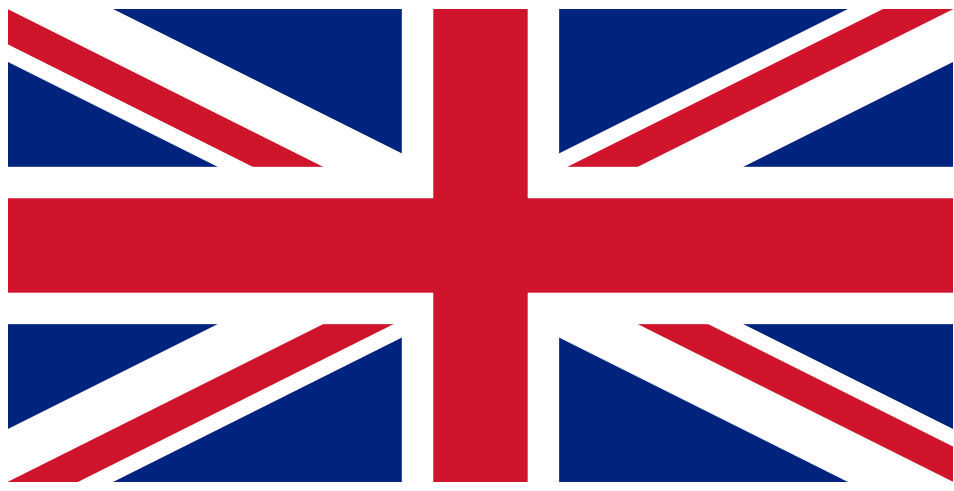
*Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylol-ysis or spondylosis, recent joint replacements or IUD'S, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we advice you to consult a doctor first*

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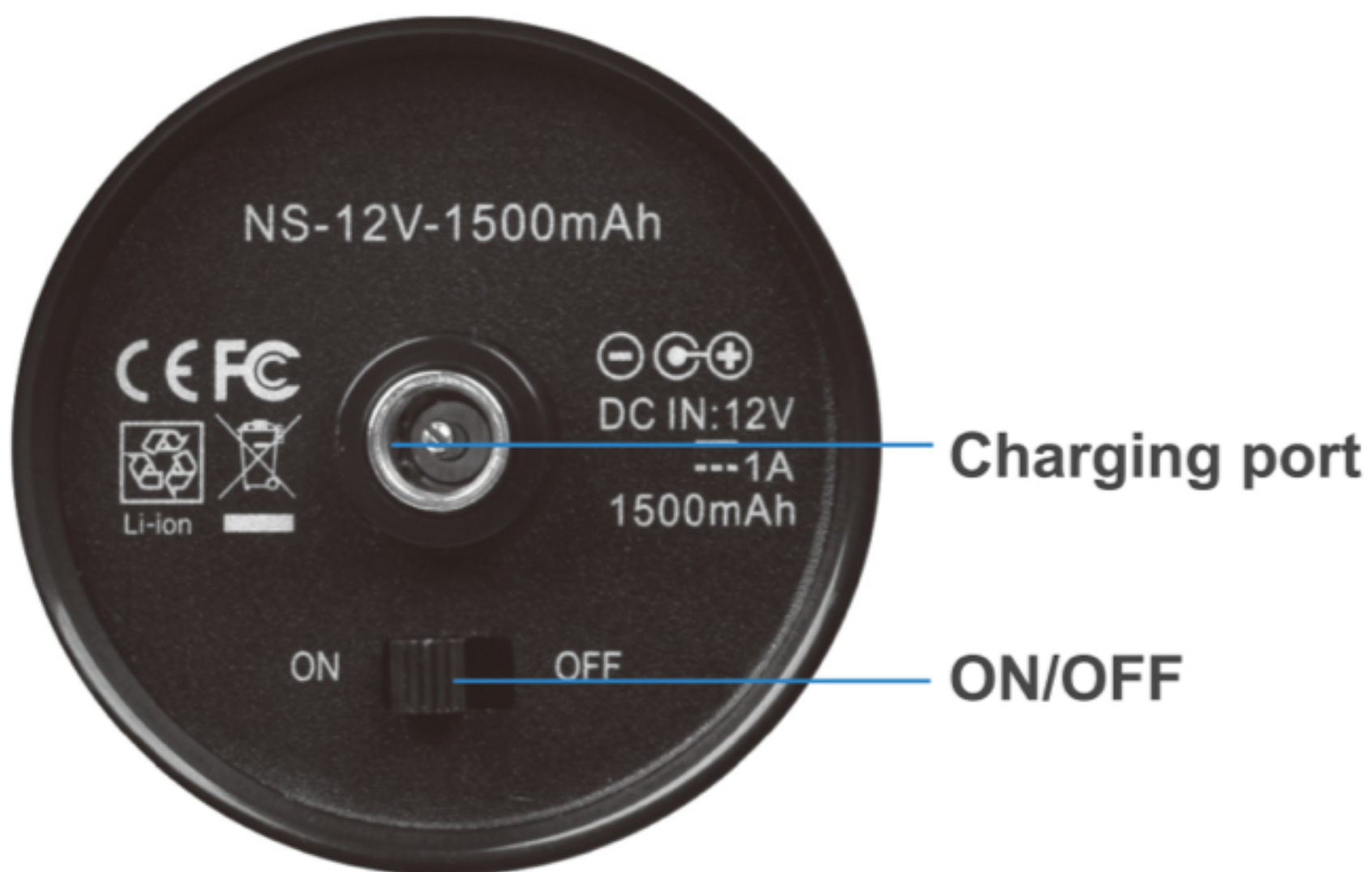
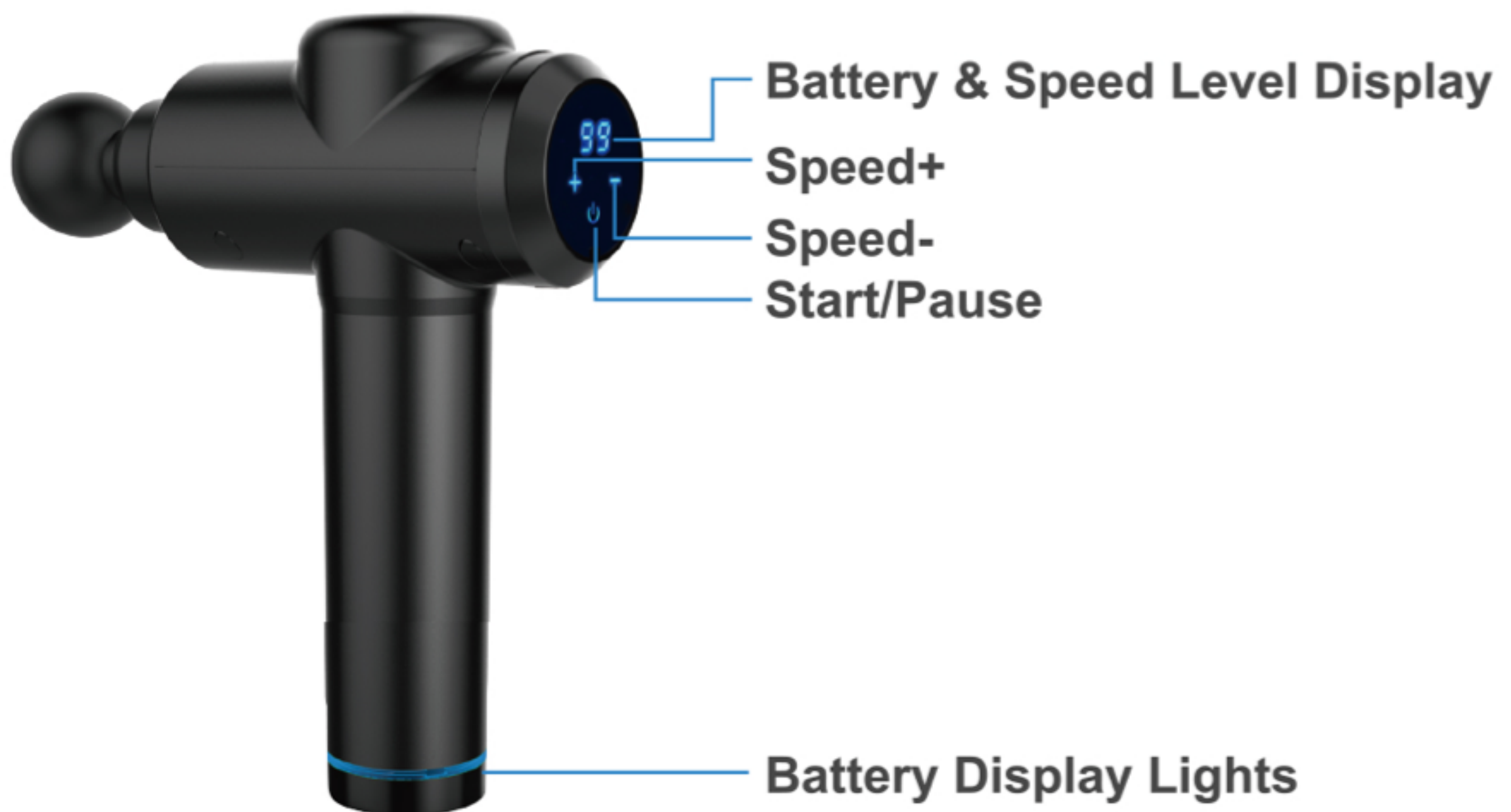
1. Switch to "ON" at the bottom to turn on.
2. Press "Start/Pause" to start to work.
3. Adjust to comfortable speed with "+" and "-".
4. Switch to "OFF" to turn massage gun off.

INFINITY

G O O D S