MAICURA **GO**

HANDLEIDING | USER GUIDE

WELKOM

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

SPECIFICATIONS

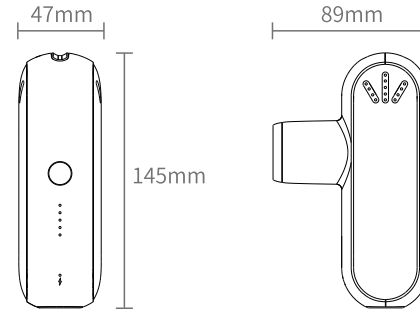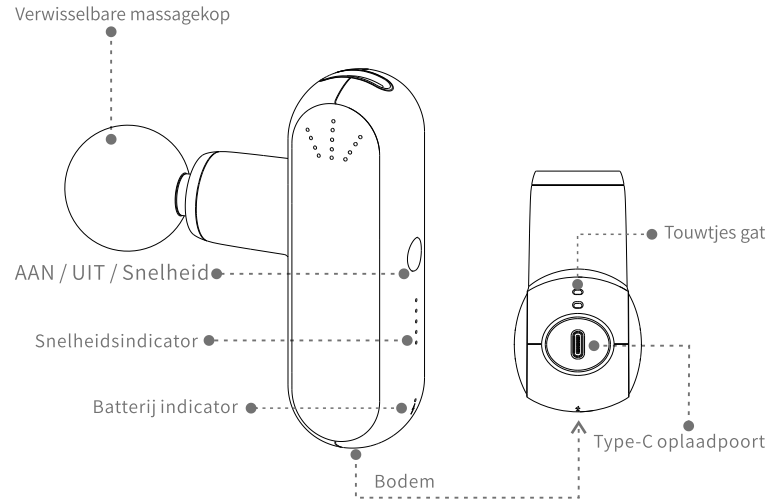Frequency: 1800-3200 R/MIN
Battery: 2600 mAh
Weight: 450g
Brand: Maicura
Address: Tourslaan 13, 5627KW, Eindhoven (NL)


WELCOME

Congratulations on your new Maicura

You are one step closer to moving better, feeling better, having less pain, and simply being a better version of yourself.

The Maicura massage gun was created to keep you doing what you love. Whether you're an athlete, a hard worker, want to relax, or simply trying to get your body to feel its best: the Maicura massage gun will help you get there.

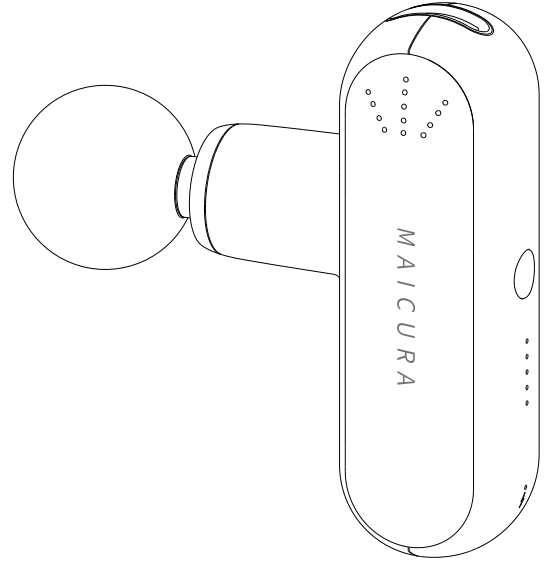QUICK START GUIDE

Lets get started

1. Turn ON/OFF

Turn on the massage gun by pressing and holding the middle button for 2 seconds, the massage gun will automatically start at speed (1). Then turn off the massage gun by pressing and holding the middle button for 2 seconds.

2. The preferred speed

Choose the desired speed (1-5) by pressing the middle button once. The blue lights below the middle button indicate the current speed.

2. Changing attachments

Choose the desired massage head and insert it into the shaft/opening while firmly pressing it. If you are using the fork head, align the protrusion on the fork head with the slot on the shaft/opening and push. If you are using the infrared head, it will automatically turn on when the head is inserted into the massage gun and you turn on the massage gun.

4. Recharging

Use the provided Maicura charging cable to charge the massage gun. Charge the massage gun to 100% when using it for the first time. The massage gun is charging when the bottom light on the massage gun turns red, and it turns off when it's fully charged. When the battery of the massage gun is below 30%, the bottom light turns red, when the battery is below 20%, the bottom light starts flashing red. When the battery of the infrared head is below 20%, the middle light on the infrared head starts flashing. A full charge is indicated when the middle light stays continuously red.

5. Taking care of your massage gun

Use a damp, clean, cloth and gently wipe your massage gun. Make sure the battery charger is not attached and that the massage gun is turned off. Store in a clean, cool place when not in use.

HEAD ATTACHMENTS

Choosing the best massage head

Ball

Suitable for all muscle groups.

Cushion

Suitable for sensitive areas such as wrists, ankles, knees, and around the neck.

Flat

Suitable for large surface areas such as hamstrings, quads, back, and chest.

Bullet

Suitable for deep tissue and trigger points.

Fork

Suitable for longer muscles such as calves, forearms, biceps, and ankle/Achilles area.

Infrared

Suitable for all muscle groups. Helps with improving the circulation of blood, the relaxation of stiff and sore muscles, and relief from pain. Note: It takes a little while for the massage head to warm up.

INSTRUCTIONS

WARNING: TO REDUCE THE RISK OF ELECTRIC SHOCK, FIRE, PERSONAL INJURY, OR PROPERTY DAMAGE, THIS PRODUCT MUST BE USED BY THE FOLLOWING WARNINGS, CAUTIONS, AND SAFETY INSTRUCTIONS:

IMPORTANT SAFETY INSTRUCTIONS

The massage gun should be used by applying light pressure and moving it over the body for approximately (60 - 120) seconds per muscle group. The massage gun should only be used on soft tissue areas of the body and should not be used on any hard areas, such as bones, including the head. Please immediately stop using the massage gun upon experiencing any pain or discomfort, apart from light muscle soreness. When using the Maicura massage gun, regardless of the pressure exerted, there is a risk of bruising. If this is the case, the massage gun should not be used on the affected area until it has completely healed. This also applies to other bruises, contusions, rashes, and, other irritated or injured areas of the skin until they are completely healed.



The massage gun must be kept away from liquids and heat sources at all times. It must therefore not be immersed in water and be stored in a dry place. The ventilation holes on the top of the massage gun should be kept free of dust, sand, dirt, food, and/or other substances that could enter and clog the holes. In addition, the massage gun should never be used for more than 1 hour in a row and the massage gun should rest for at least 30 minutes after each use before switching on again.

After charging the massage gun, the charger should be unplugged from the massage gun before use. Before disposing of the massage gun the battery must be removed. The battery must be complied with safely at an appropriate recycling facility. The massage gun should never be used for any purpose other than its intended use, should not be dropped, and should never be tampered with or altered in any way.

PLEASE DO NOT USE THE MESSAGE GUN, OR ANY PERCUSSION OR VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM THE PHYSICIAN IF ANY OF THE FOLLOWING APPLY:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery or injury, epilepsy or migraines, herniated disks, recent joint replacements or IUDs, metal pins or plates, or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any percussion or vibration device. These contraindications do not mean that you are not able to use a percussion or vibration device but we advise that before doing so you consult your physician or doctor.

SPECIFICATIES

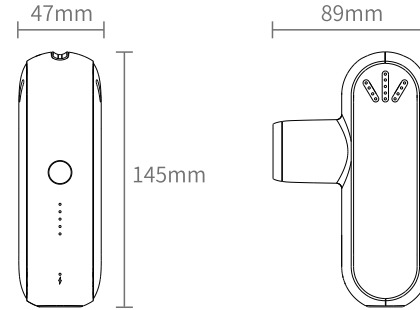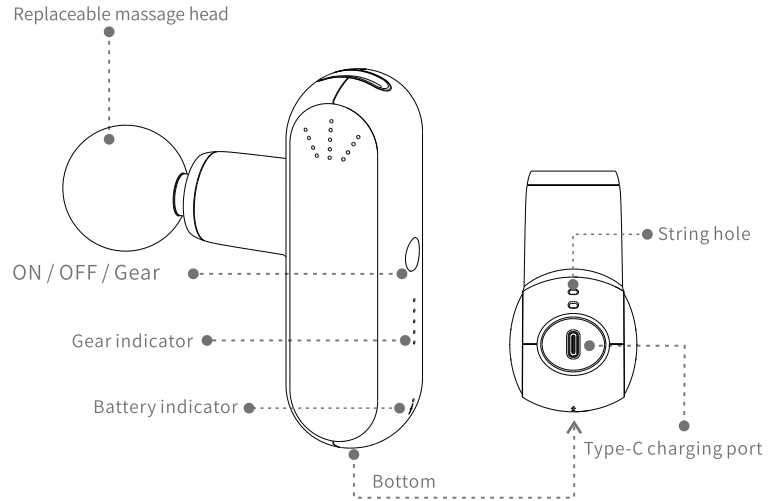Frequentie: 2000-3200 R/MIN
Batterij: 2600 mAh
Gewicht: 450g
Merk: Maicura
Adres: Tourslaan 13, 5627KW, Eindhoven (NL)