

Working principle of laser hair removal apparatus

Laser hair removal device is a personal care appliances, so for the sake of safety, should be used carefully. Please read this manual carefully before using the laser hair removal device for the first time.

This product is directly worked to the hair root by laser. Melanin in the hair roots absorbs light pulses from the product. Light pulse is a kind of soft and non-continuous light source with good selective photo thermal effect, which is mainly used in life and beauty, and it is a mild light source of beauty and skin care with high safety, effectiveness, non-invasion and non-trauma. The darker the hair, the more light it absorbs. This illuminant can stimulate hair to enter dormant period. During irradiation, hair naturally falls off and is inhibited from regenerating. It can take a week or two for the hair to start falling out. Good hair removal cannot be achieved in one step because only growing hair responds to light hair removal. During depilation, the hair at the depilation site is usually at different stages of the hair growth cycle. To be effective on all hair, repeat irradiation at least once a week for the first two months. The product is integrated with an optical filter to ensure that uv rays do not reach the skin and that the product is safe


The laser hair removal device will flash when running. Do not look at the light source directly!

Do not use laser hair remover if your skin is already tanned or will be in the sun recently. Using it in these situations can cause severe burns

or serious damage to your skin. Avoid exposure to the sun 4 weeks before using

a laser hair removal device. Do not use a laser hair remover near your eyes! Laser depilates are not effective on natural white, gray, gold, and red hair. If your own hair is of these colors, do not use.

Do not use laser hair remover for natural dark skin! Before using the laser hair removal device, please read carefully and follow all safety warnings and instructions of the laser hair removal device

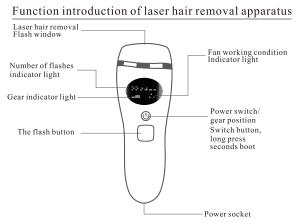

Product parameters	
Product model: IPL-DI9	rated power: 24W
Input voltage: 12V	current: 3A

Carton content When you open the carton, you can find that it contains the

following parts 1.Laser hair removal equipmen

2.A C adapter

4.Protective glasses 5.Shaving cutter

Preparation of laser hair removal apparatus

Before use, thoroughly clean the skin and shave off the hair (1mm or less) in the hair removal area, the skin is hairless, completely dry, and free from oily residues (perfume, lotion, body fluids, Residues such as sunscreen) to ensure that light energy is directed to the hair roots in the most efficient way possible.

Note: Some poor condition can happen when laser hair removal is performed on unshaven skin.

- 01. The dust in the hair may absorb the irradiation energy, which may cause sever skin pain or redness and swelling. Forked hair or dust may pollute the light window and block the normal light output, which may reduce the service life of the produc
- 02. Hair on the skin can get burnt and can severely scald the top layer of the skin. Check the machine body and light outlet for cleanliness, hair, dust and fiber etc Contaminant, if have, please wipe carefully clean.
- When you use the product for a period of time, please pay attention to check whether your skin color has changed significantly -- only pointing to the bad direction change, such as white to black, purple, etc., you need to go to the hospita

If you use this product for the first time, your skin color is within the available range, and you can use it in the first gear. If there is no obvious pain, you can continue to add the gear until there is slight pain in the irradiated skin. This is the

Use of laser hair removal apparatus

Laser hair remover is used to remove excess hair. It can also permanently reduce hair regurgitation. It can effectively remove excess body hair through four cycle hair removal steps 1-3 times a week

Post-treatment care

Skin care

- 1. It is recommended to avoid water contact within 6 hours after light operation. After that, moisturizer can be used to moisturize the irradiated part 2.Do not use perfume, facial mask and cosmetics containing frosted beads after
- 3. Avoid exposure to uv rays within 24 hours after light treatment. Apply sunscreen
- 4.Consult your surgeon if you feel unwell.
- 5. Sunbathing is likely to cause skin blackening or sunburn. Irradiation at this time will cause redness, inflammation and other conditions. Therefore, before using this product, pay attention not to be insolated.
 6. It is strictly prohibited to insolate the irradiated area within 48 hours after the laser hair removal operation. If insolating is required after 48 hours, please always pay attention to whether there are adverse conditions on the skin. Please
- se sunscreen to prevent sunburn after 72 hours.
- 7. Artificial tanning is not allowed during light treatment 8. If need plastic surgery after light (exposure area only) please consult a doctor.

• Possible side effects

After shaving \(\) exposing to light and ice compress treatment for a similar period of time, the skin becomes dry and itchy. Rest assured, this reaction is harmless and won't take long to wear off. You can apply a mild, non-irritating moisturizer.

Safety of laser hair removal apparatus

Safety is the first criterion for laser hair removal apparatus

IPL pulse technology provides excellent safety and low energy levels

Intense pulsed light technology can achieve long-term hair removal, so a small amount of this energy can be used in other optical hair removal equipment. The low energy of laser hair removal equipment can reduce the possibility of injury or complications, and contribute to overall safety.

The laser hair removal device has a built-in eye protection function, and the light outlet of the product must touch the skin to work, which can effectively protect both eyes and skin. The laser hair removal device is turned on but not used. It will automatically shut down after 10 minutes.

Important safety information - please read before use!

consult your physician or dermatologist)

Laser hair removal is not for everyone. Please read and consider the information in the following sections before using. 1.If you are unsure whether the laser hair removal device is suitable for you

2.If you are naturally dark, do not use a laser hair removal.)

The laser depilator removes excess hair by selectively removing pigment from ne hair, and different amounts of pigment can also exit in the surrounding tissue: of the skin. The amount of pigment in a person's skin is determined by the color of their skin and the degree of risk that using a laser depilator on a dark complexion can lead to burns, bubbles and changes in skin tone (pigmentation or hypo

3. Do not use a laser hair remover on tanned skin or after recent exposure. Tanned skin, especially after exposure to the sun, contains a large amount of melanin, the presence of a large number of melanin skin using laser hair remova will be at a higher risk. Side effects such as burns, blisters, and changes in skin color can result from the use of a laser hair remover on recently exposed skin.

What situation does not use laser depilator

Do not use the device if you have any of the following conditions

- 1. If you are pregnant or nursing If you have been exposed to the sun or had an artificial skin graft in the past 4
- 3. If you have a tattoo or permanent epilepsy in the area to be treated.
- 4.If you have dark spots, such as freckles, birthmarks, moles in the area to be
- 5. If you have eczema, psoriasis, lesions, exposed wounds or active infections
- such as lip blisters in the treatment area, wait until the wound heals before using
- 6. If you have A knot in one's history of formation, A known light (light), or drugs make the skin more sensitive, including the blade body anti-inflammatory agents (such as aspirin, ibuprofen and acetaminophen), tetracycline, command
- thiamethoxam become, thiamethoxam become, diuretics, urine, sulfa drugs, three by the mi makes amine, fluorouracil, VLB, griseofulvin, fruit acid, all the vitamin A acid, accutane or topical retinoids.
- 7.If you have a skin condition caused by diabetes, such as another systemic or metabolic disease.
 8.If you are currently taking alpha hydroxy acid, betadate, retin-a, acne wonder
- rug (isotretinoin), topical retinoic acid or azelaic acid. .if you have been taking A acne specific (isotretinoin) for the past 6 months 10. If you have had a hmong compound regimen in the past 3 months.

 11. If you have had an outbreak of herpes in the area you are treating, at least
- onsult your physician or receive preventive treatment before using the device. 12.If you have epilepsy
- 13.If you have an active implant, such as a pacemaker, an incontinence device, ar slet pump, etc
- 14.If you have related photopenic diseases, such as bugonorrhea, polymorphic
- photopolysis, solar urticaria, lupus erythematosus, etc.

 15.If you have a history of skin cancer or potential skin malignancies, etc.

 16.If you've had radiation therapy or chemotherapy in the last 3 months. 17. If you have a condition where your doctor says you're not fit for a depilation
- 18. If you are unsure whether the laser hair removal device is safe for you to use, nsult your physician or dermatologist.

How to use laser hair removal device safely

- .Choose the intensity level carefully
- The intensity level refers to the highest level released to your skin using the flash. The brightness level of the indicator can be set on the instrument. As the intensity level increases, so does the risk of side effects.
- 2. Always use a laser depilator from a low intensity level
 If you feel no or very little discomfort during the laser hair removal level, you can
 go up one level for the next treatment, and so on for the next treatment. 3. Avoid complications after use

1. Cover birthmarks and tattoos before use 2. Cover dark brown or black spots, such as freckles, birthmarks, and moles, before

- 2.Do not look directly at the device when light is emitted from the window. 3. Please do not use this device on the nipple and reproductive organs (including
 - 5. Please do not use this device where you want to grow hair in the future.
 - 6.Please do not use this device for any purpose other than hair removal.
 7.Please do not use the flash in the open area. During treatment, please ensure tha
 - the skin contact surface is in full contact with the skin
 - 3.If the skin contact is too hot, remove the skin. 9.Do not use alcohol or other flammable liquids (including perfume, disinfectant
 - or other substances containing alcohol) or adit to clean the skin when using.

 10. The use of this device may cause temporary pigment changes
 - 11. Keep this device away from children 12.lease do not use the device on or for children.

Laser hair removal device treatment plan

- Freatment plans and outcomes
- A normal full hair growth cycle may take 18 to 24 months, during which time multiple laser hair removal treatments may be required to achieve reduced or permanent hair loss. The effects of hair removal vary from person to person, depending on the body, hair color and how the laser hair removal device is used

Dehydrogenation purifies the skin for 8 weeks According to the growth of hair use many times, the general 8 weeks can achieve



4-5 week 2 times

delicate muscle!

air decreases in growth, maintain light rebound, maintain the Depending on hai

After 8 weeks

TO


say goodbye to

6-7 week Once a week

body smooth and delicate!





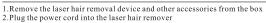

Application method of laser hair removal apparatus

- 3.Plug the adapter into the socket
- 4. Use a shaving knife to clean the surface of the skin 5. Aim at the skin, keep the flash mouth perpendicular to the skin for 90 minutes, press the flash button to remove hair

Procedures of laser hair removal apparatus





