

MASSAGE PILLOW

OPERATION MANUAL



**Please keep this instruction manual properly,
And read it carefully before using the product.**

The company reserves the right to improve the product and instruction manual.
Any change is subject to notice. The appearance of the product is subject to the
actual product. 本公司保留对产品进行改进、修改及更新的权利，产品外观以实物为准。

Thank you for purchasing our massage pillow products. This massage pillow adopts ergonomic design to massage back, waist, buttocks, legs, arms and other parts. It is a massage product with warm feeling and hot compress function. Before use, please read this manual carefully and pay special attention to safety precautions to ensure the correct operation and use of the product. Take good care of the manual for future reference.




The pictures are for reference only.
The products are subject to the actual situation

Product parameters

Product Name	Massage pillow
Rated Voltage	12V ---
Rated Power	24W
Working Hours	10-30min

Schematic diagram of manual controller

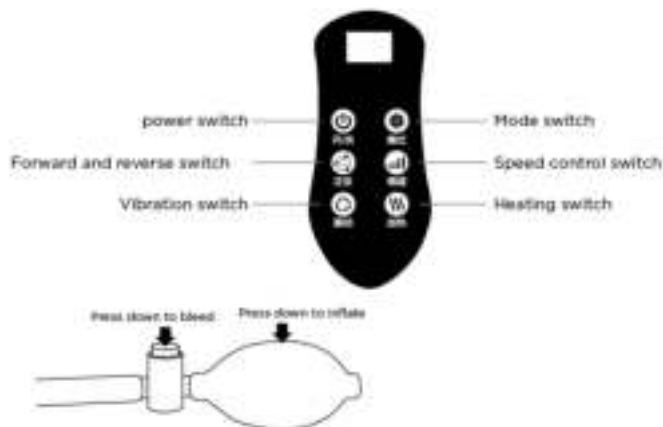

Instructions

1. Fix the massager on the part to be massaged.
2. Connect the power adapter with the power DC socket of the massager, and the massager will enter the standby mode.
3. Press the power on key to enter the automatic massage (forward and reverse automatic switch).
4. After start-up, warm start, first gear warm down, second gear warm off.
5. Press the forward and reverse switch to change the massage direction of the massage head. Press for forward massage, and then press for reverse massage.
6. Press the speed switch to change the massage speed of the massage head.
7. The vibration mode switch has three modes, and the fourth gear turns off the vibration function.
8. When the massager works, press the power on key to turn off the massage function, and the massager enters the standby state.

Packing list

Massage pillow × 1

Instruction / Certificate / warranty card × 1

Repair and maintenance

- ① Please use neutral detergent and soft cloth for daily cleaning. Do not use corrosive detergent.
- ② If the power cord is damaged, it must be replaced by the full-time personnel of the maintenance shop designated by the manufacturer in order to avoid danger.
- ③ Please unplug the power supply when not using or cleaning to avoid danger.
- ④ Do not use any sharp objects to cut the surface of the product.
- ⑤ Do not put the product in a humid environment, but in a dry and ventilated place.
- ⑥ Do not press or apply pressure on the product.

Safety precautions and warnings

- People with long hair need to roll up their hair to avoid being pinched by massage head and scratching scalp.
- Massage head cloth cover or external leather cover, mesh cloth damage, strictly prohibited to use.
- Please do not use the hand-held power cord or manual control cable to let the body hang or swing, which will cause the wire core to be disconnected and unable to work normally.
- After use, please hold the plug and then pull it out. Do not pull the power cord directly to avoid accidents.
- Prevent water or other corrosive liquid from entering into the machine, causing malfunction and damaging the machine.

Simple fault and treatment method

- When the machine fails or fails to start, please check the following items. If the product is still in trouble, please send the product to the designated maintenance unit or authorized technical personnel for inspection and repair.
- Do not try to remove or repair this product by yourself, otherwise it may cause personal injury, and we will not be responsible for it.

Common faults	Possible causes	Processing method
The whole machine does not work	Is the adapter plugged into a power outlet? Does the manual controller turn on the switch? Is the adapter DC plug and the massage DC socket paired?	Insert the adapter into the power socket correctly and firmly. Turn on the switch on the manual controller first. The adapter DC plug and massage DC head are firmly connected.
Massage and rolling are slow and weak	Power failure If the temperature of the movement is too high, it will stop automatically.	Use after call. When the temperature of the movement is reduced, the force of massage shall be carried out.
Load noise	Running in period of new machine	After 10 hours of use, the noise will naturally decrease.

If the above methods can not be solved, please contact our local distributors.





All services for this product must be provided by our professional staff.
If the power line is damaged, in order to avoid danger, it must be replaced by the manufacturer or its maintenance department or similar full-time personnel.

People who are not suitable to use this product

- (1) Malignant tumor patients;
- (2) Patients with acute diseases such as heart disease;
- (3) Pregnant women or in menstruation;
- (4) Children or those who have undergone surgical operation;
- (5) Weak and sickly;
- (6) The surface of the product is hot, and it should be used with caution by the heat sensitive person.

Please read all instructions carefully before use

- If you have health concerns, please consult your doctor before using this product;
- Consult your doctor before using;
- Do not let children use this product alone;
- Do not use the same part for more than 15 minutes;
- If the product is used for a long time, please stop using it until the product cools down;
- Do not use this product directly on swollen, inflamed or rotten skin areas;
- This product is a non professional instrument and can not be used to replace drug therapy;
- Please do not use it before going to bed, because this product has stimulating effect and can delay sleep;
- Do not use this product when sleeping;
- Do not use this product if you have difficulty in movement or paralysis of lower limbs due to physical defects;
- When children or people with disabilities use this product, other adults must be present;
- If the surface of the appliance is heated, attention should be paid to those who are not sensitive to heat.

Certificate

Product Name: Message Pillow

Inspection Date: _____

Inspector No: _____

