



M1825 YUNMAI-X Mini 2 Smart Body Fat Scale User Manual

[Home](#) » [YUNMAI](#) » M1825 YUNMAI-X Mini 2 Smart Body Fat Scale User Manual

M1825 YUNMAI-X Mini 2 Smart Body Fat
Scale User Manual



YUNMAI-X

User Manual

Dear customer,

Thank you for purchasing our product. Please read the following instructions carefully before first use and keep this user manual for future reference. Pay particular attention to the safety instructions. If you have any questions or comments about the device, please contact the customer line.



www.alza.co.uk/kontakt



+44 (0)203 514 4411

Importer Alza.cz a.s., Jankovcova 1522/53, Holešovice, 170 00 Praha 7, www.alza.cz

Contents [[hide](#)]

- [1 Warning](#)
- [2 Important](#)
- [3 Notice before use](#)
- [4 Note](#)
- [5 Weighing Notice](#)
- [6 Maintenance](#)
- [7 Product Introduction and Installation](#)
- [8 Install YUNMAI App](#)
- [9 Connect to the scale via Bluetooth](#)
- [10 Start Using](#)
- [11 YUNMAI App](#)
- [12 Product Specifications](#)
- [13 Warranty Conditions](#)
- [14 EU Declaration of Conformity](#)
- [15 WEEE](#)
- [16 Documents / Resources](#)
- [17 Related Posts](#)

Warning

Do not disassemble, squeeze, puncture the product or put it into water, fire or expose it to temperature higher than 60°C (140°F)

Important

YUNMAI Smart Scale should not be used to diagnose or treat any medical situation.

Notice before use

- Please ignore the first weighing reading at your first time use or after moving the scale.
- After the scale power on or reset, please wait for the LED display flash “ ” for several seconds, and then start your weigh-in.
- After each move of the scale, please wait for 5 seconds for the sensor to self recalibrate with the floor, and then weigh for accurate reading.

Note

Following groups may not get accurate results:

- Juveniles
- Professional athletes
- Pregnant women
- Users with Cardiac Pacemaker device

Weighing Notice

- Do not weigh during charging.
- The scale surface may become slippery in humid places, please watch your steps.
- Before meal or before baths is the most recommended time for weighing, because eating and bathing may affect the body moisture and temperature, which might affect the sensor accuracy.
- Recommend to weigh yourself at the same time and same condition in each day for more accurate data record.


Maintenance

- Use a damp cloth to wipe the surface of the scale.
- Do not use any detergent to clean the scale, because detergent might damage the metal electrode sheet which may affect the accuracy of detecting body fat data.

Product Introduction and Installation


Thank you for choosing YUNMAI X Smart Scale.

Concealed display screen

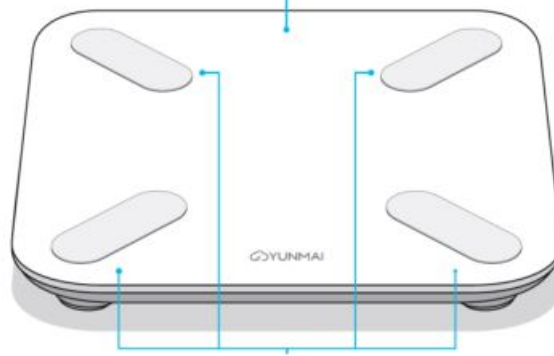

Metal electrode sheet

Scale foot

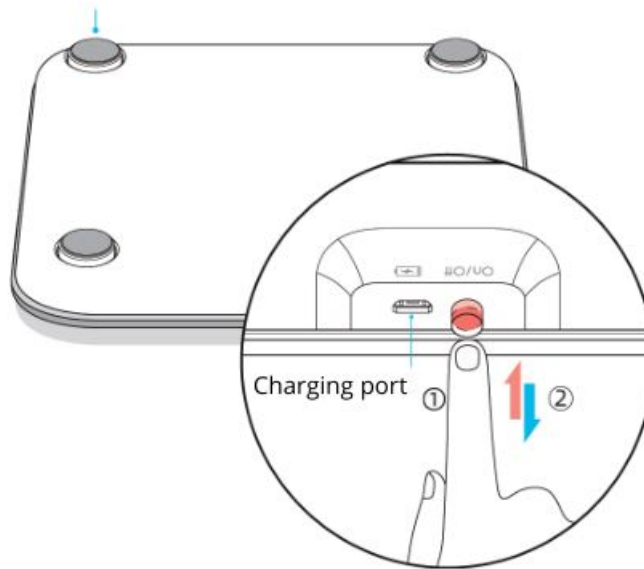

1. Note: Charging port. When in low power, plug the charger cable into the charging port to recharge the battery, red light indicates the battery is charging, while green light indicates the battery is fully charged.
2. Note: Power Switch. Press the button down to turn it on.

Install YUNMAI App

Search “YUNMAI” in the Android Application Market or the iOS App Store or scan the following QR Code to download, install it and it will get connected to the scale once turned on.

App downloaded via other channels may be unusable. Before pairing with your equipment, make sure your equipment meets the following requirements:

- Compatible with iPhone/iPod touch with iOS 10.0 and above.
- Android 4.3 and above with Bluetooth 4.0.




Scan QR Code


Download

Connect to the scale via Bluetooth

Open YUNMAI App and Bluetooth and it will get connected automatically and your realtime data will be synced to the App

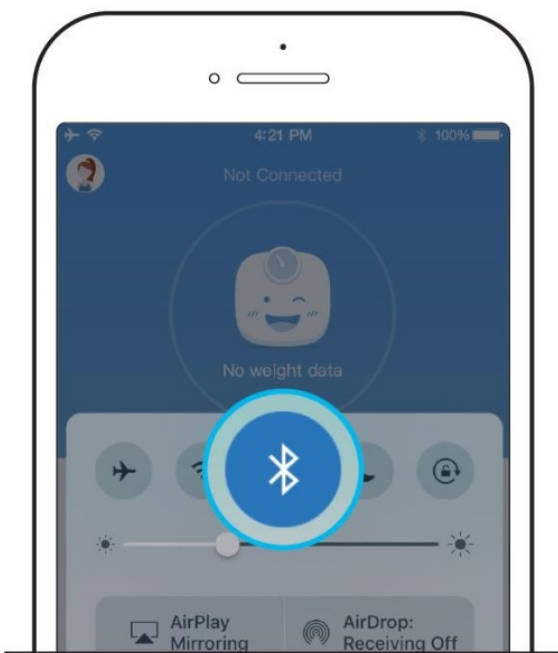

iOS: Slide from the bottom of the screen to open the Control Center and click on (Bluetooth) Logo.

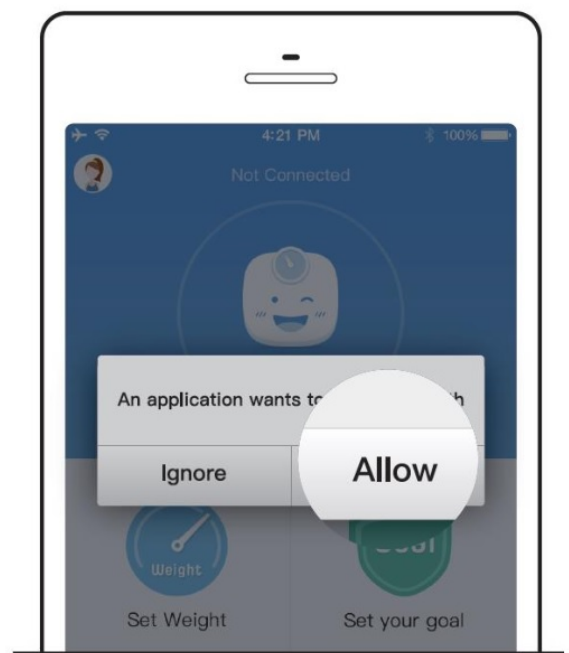


Android: Click on (Allow) to turn on Bluetooth when the system reminds you when to open the YUNMAI App.

Start Using

Correct Standing Posture

Making sure your feet placement is as correct as shown below:

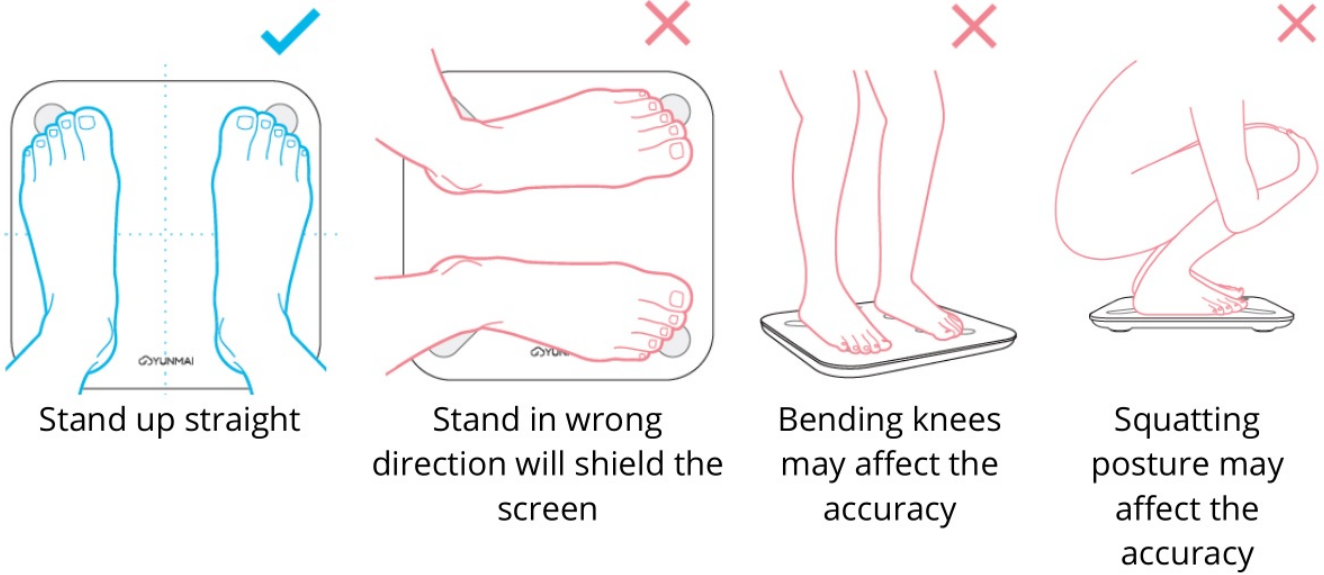

To ensure the accuracy, please stand on the scale barefoot with the right posture.

Display



Body Mass Index

BMI (for men and women)					
Age	Underweight	Normal	Overweight	Obese	Grossly Overweight
18-120	<18.5	18.5-24.9	25-29.9	30-34.9	≥35

Body mass index (BMI) = weight (kg) ÷ height ^ 2 (m)

For example: BMI= 70kg ÷ (1.75m x 1.75m) = 22.86, belongs to the normal range.

Note: BMI Index is a standard that is widely used to measure a person’s health condition.

YUNMAI App

Using Yunmai Smart Scale with the assistance of YUNMAI App, you can manage your personal health data. Yunmai will generate a brief health report and improvement recommendations for you.



Smart Scale

Smart Rope



Connected



50kg

8 hours ago

22.6

BMI

19.2%

Body Fat

Tracking Graphs



Target



Record Weight



Body Report



Home



Me

The picture is for reference only. For more information, please download the latest version of YUNMAI App.

Product Specifications

Product name:	YUNMAI X Smart Scale
Model:	M1825
Measurement unit:	Kg/Lb
Connection:	Bluetooth 4.0
Size:	260 x 260 x 26mm
Net weight:	1.2 Kg
Weighing range:	3 Kg-180 Kg
Accuracy:	± 50g
Working temperature:	23°F-104°F (-5°C-40°C)
Exterior Material:	Eco-friendly ABS
Battery capacity:	1000mAh
Input voltage:	DC 5V
Battery type:	Lithium battery
Charging input:	DC 5V 1A
Input current:	0.5A

Warranty Conditions

A new product purchased in the Alza.cz sales network is guaranteed for 2 years. If you need repair or other services during the warranty period, contact the product seller directly, you must provide the original proof of purchase with the date of purchase.

The following are considered to be a conflict with the warranty conditions, for which the claimed claim may not be recognized:

- Using the product for any purpose other than that for which the product is intended or failing to follow the instructions for maintenance, operation, and service of the product.
- Damage to the product by a natural disaster, the intervention of an unauthorized person or mechanically through the fault of the buyer (e.g., during transport, cleaning by inappropriate means, etc.).
- Natural wear and aging of consumables or components during use (such as batteries, etc.).
- Exposure to adverse external influences, such as sunlight and other radiation or electromagnetic fields, fluid intrusion, object intrusion, mains overvoltage, electrostatic discharge voltage (including lightning), faulty supply or input voltage and inappropriate polarity of this voltage, chemical processes such as used power supplies, etc.
- If anyone has made modifications, modifications, alterations to the design or adaptation to change or extend the functions of the product compared to the purchased design or use of non-original components.

EU Declaration of Conformity

Identification data of the manufacturer's / importer's authorized representative:

Importer: Alza.cz a.s.

Registered office: Jankovcova 1522/53, Holešovice, 170 00 Prague 7
CIN: 27082440

Subject of the declaration:

Title: Personal Scale
Model / Type: YUNMAI-X

The above product has been tested in accordance with the standard(s) used to demonstrate compliance with the essential requirements laid down in the Directive(s):

Directive No. 2014/53/EU
Directive No. 2011/65/EU as amended 2015/863/EU
Prague



WEEE

This product must not be disposed of as normal household waste in accordance with the EU Directive on Waste Electrical and Electronic Equipment (WEEE – 2012/19 / EU). Instead, it shall be returned to the place of purchase or handed over to a public collection point for the recyclable waste. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in fines in accordance with national regulations.



Documents / Resources

 YUNMAI

YUNMAI-X
Mini Smart Body Fat Scale

[YUNMAI M1825 YUNMAI-X Mini 2 Smart Body Fat Scale](#) [pdf] User Manual
YUNMAI-X, Mini 2, Smart Body Fat Scale, YUNMAI-X Mini 2 Smart Body Fat Scale, Fat Scale,
Scale, M1825 YUNMAI-X Mini 2 Smart Body Fat Scale, M1825