
## Pure (Deluxe)

H4101

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DEL USUARIO
- PT MANUAL DO UTILIZADOR
- GR ΕΓΧΕΙΡΙΔΙΟ ΧΡΗΣΤΗ
- TR KULLANIM KILAVUZU
- IT MANUALE PER L'UTENTE
- NO BRUKERHÅNDBOK
- SV ANVÄNDARMANUAL
- DA BRUGERMANUAL

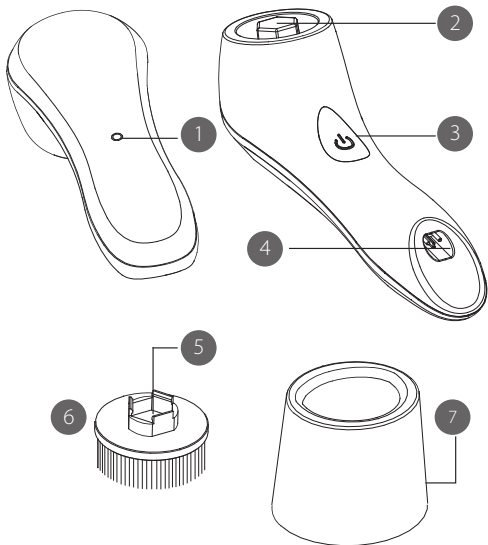---

## Pure Deluxe

# Table of contents

1. Warnings & contraindications.....	2
2. Safety with Silk'n Pure (Deluxe).....	7
3. Get to know Silk'n Pure (Deluxe) .....	7
4. Device controls & settings.....	8
5. Treatment plan.....	9
6. Get started!.....	10
7. Cleaning & maintenance .....	11
8. Troubleshooting .....	12
9. Parts & consumables .....	13
10. Customer service .....	13
11. Disposal .....	13
12. Specifications .....	14
13. Signs & symbols .....	14
14. Warranty .....	15

## 1. Warnings & contraindications

---

EN

This user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Please read all warnings and contraindications before use.

- ⚠ Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if the device is damaged in any way, does not appear to be working properly, or if you see or smell smoke. In such cases, stop using the device and contact Customer Service.
- ⚠ Only authorized personnel are permitted to perform repairs. Do not attempt to repair or open sealed parts of your device. You may expose yourself to dangerous components and cause serious harm to yourself. It will also void your warranty.
- ⚠ The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.
- ⚠ Keep the adapter away from water. The adapter is not waterproof.

- ⚠ Always ensure that your hands are dry before connecting the adapter to the electricity supply.
- ⚠ Do not place or store the adapter where it can fall or be pulled into a tub, sink, water or any other liquid.
- ⚠ Do not reach for the adapter if it has fallen into water. Unplug immediately before retrieving it. Failure to unplug the adapter may result in an electric shock!
- ⚠ Never use the adapter if it has a damaged cable or plug, it isn't working properly, or has been submerged in water. If that happens, unplug the adapter immediately.
- ⚠ Never pull the adapter out of the socket by the cable. Always remove it by taking hold of the adapter.
- ⚠ Place the adapter cable in such a way that it does not pose a tripping hazard.
- ⚠ Do not kink the adapter cable and do not lay it over sharp edges.
- ⚠ Do not leave the device submerged in bath or shower water for more than 20 minutes.

- ⚠ Do not attach wet brush heads when charging the device.
- ⚠ Do not allow face wash to accumulate on the brush heads.
- ⚠ The device is only to be used with the power supply unit provided.
- ⚠ If a battery has leaked, avoid contact with the skin, eyes and mucous membranes. If necessary, rinse the affected areas with water and seek medical assistance.
- ⚠ Do not clean the device or any of its parts in the dishwasher.
- ⚠ Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. For details, see section "11. Disposal".
- ⚠ The battery must be removed from the appliance before it is scrapped.
- ⚠ The appliance must be disconnected from the supply mains when removing the battery.

- ⚠ The battery is to be disposed of safely.
- ⚠ Attention should be drawn to the environmental aspects of battery disposal. Don't throw used batteries in a dustbin. Please contact your retailer in order to protect the environment.

## **Battery Removal**

To remove the battery, the device must be broken apart. The device will no longer function after this procedure, nor can it be repaired. Any warranty will be void when the appliance housing is broken and the battery is removed. It is strongly recommended to wear safety protective goggles and gloves.

Silk'n is not liable for any physical or material damages as a result of performing this procedure.

To remove the battery place the device on a hard, strong surface. Use a hammer to hit the appliance on the seam and continue to do so along the whole seam on both sides of the device.

Once the seam is sufficiently broken, use pliers or a screwdriver to take the housing apart. Unscrew the PCB (printed circuit board) and remove the battery by disconnecting the plug from the PCB.

## Contraindications

- ⚠ This device can be used by children aged from 8 years and above, and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the device in a safe way and understand the hazards involved.
- ⚠ Children must not play with the device or any of its parts.
- ⚠ Cleaning and maintenance must not be done by children without supervision.
- ⚠ For hygienic reasons, do not share the device with others.
- ⚠ Avoid contact with the eyes.
- ⚠ Do not use a dry brush head on the face or neck, because it may cause skin irritation.

### **Do not use this device if you have:**

- ⚠ sore or injured skin on the face or neck.
- ⚠ a pacemaker.

## Consult your physician if you have:

- ⚠ any medical condition or any kind of sensitivity.
- ⚠ acne or any other skin disorder.
- ⚠ a rash, redness or itching during or after use.

## 2. Safety with Silk'n Pure (Deluxe)

---

Silk'n Pure (Deluxe) has been tested for water ingress to the extent covered by IPX7. This means your device is waterproof and suitable for use in the shower or while taking a bath, but do not leave it submerged in water for more than 20 minutes.

The device also has an overcharge protection. This means that you can safely leave the adapter in the device even if it is connected to the electricity supply, without the risk of damage.

## 3. Get to know Silk'n Pure (Deluxe)

---

This user manual refers to both Silk'n Pure and Silk'n Pure Deluxe.

Silk'n Pure (Deluxe) effectively cleanses the skin and removes the outer layer of dry and flaky skin cells. It lifts out impurities from within the pores as well as any build-up of debris from within the folds of fine lines and wrinkles. The device uses vibrations to help smooth the skin's surface, and massages your skin to improve blood circulation.

Silk'n Pure (Deluxe) removes up to 6 times more makeup than manual cleansing, leaving your skin clean, smooth and soft.

Silk'n Pure (Deluxe) is an electrical device and should be used with special attention to safety. **Before use, it is important that you read and understand the instructions on how to use the device including after-use procedures, and follow them strictly.** We recommend you re-familiarise yourself with the instructions in this manual before each treatment.

### 3.1 Package contents

- Your Silk'n Pure (Deluxe) device
- 2 brush heads (Regular, Soft)
- Face brush stand
- Adapter
- User manual

Silk'n Pure Deluxe includes the following extras:

- Travel case
- Face wash travel bottle
- Face Massage brush head
- Eye Massage brush head

## 3.2 Device description

Silk'n Pure (Deluxe) is a small, hand-held, electrically powered skin cleanser. Refer to the diagram at the front of this manual to see the location of each feature.

- Charging light (#1)
- Brush attachment point (#2)
- Control button (#3)
- Adapter socket (#4)
- Brush head center (#5)
- Brush head (#6)
- Face brush stand (#7)
- Travel case (#8) \*
- Face wash travel bottle (#9) \*
- Face Massage brush head (#10) \*
- Eye Massage brush head (#11) \*

\* Silk'n Pure Deluxe

## 3.3 Intended use

Silk'n Pure (Deluxe) is an over the counter device intended for cleansing your face and neck, and massaging your skin. Do not use for any purpose other than as described in this manual.

# 4. Device controls & settings

---

## 4.1 Charging

Before first use, your device needs to be fully charged. For best results, charge your Silk'n Pure (Deluxe) for 6 hours before using it for the first time, even if there is some residual charge left in the device as a result of quality control procedures.

1. For your safety, make sure the device is completely dry before charging.
2. Connect the **adapter** to the device and the electricity supply.
3. When charging, the **charging light** is red and blinks continuously. When fully charged, the **charging light** is constant red.

 If you are charging for the first time but the **charging light** is constant red, you still need to charge the device for 3 hours.

4. We recommend you unplug the device and store the **adapter** after each full charge until you need to charge the device again.

 We recommend you charge the device fully at least every 3 months to sustain battery life.

 Do not attach wet brush heads when charging the device.

## 4.2 Indicator lights and function settings

Your Silk'n Pure (Deluxe) is easy to use and has one **control button** and one **indicator light**. Refer to the diagram at the front of the manual to see the location of each feature.


Button	Action	Function / Mode
Control button (#3) 	Press x 1	Vibrates - low speed
	Press x 2	Vibrates - high speed
	Press x 3	Pulsates - low speed
	Press x 4	Pulsates - high speed
	Press x 5	Device switches off
Charging light (#1) 	Light is red and blinks for 10 seconds	Low battery
	Light is red and blinks continuously	Charging
	Light is constant red (visible only when connected to the electricity supply)	Fully charged. The light will go out when the device is disconnected from the <b>adapter</b> , or when you remove the <b>adapter</b> from the electricity supply.

## 5. Treatment plan

-  For daily use, but not more than twice a day.
-  We recommend cleansing in these areas: forehead, nose, chin, cheeks and neck.
-  The recommended treatment time is 10 seconds per area.
-  Silk'n Pure (Deluxe) is designed to cleanse your skin in just 1 minute. We recommend that you do not to exceed this time.

## 6. Get started!

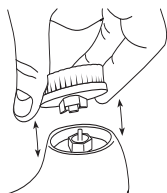1. Select the **brush head** suitable for your skin type. Gently push the **brush head** onto the **brush attachment point**, until you hear a click.
2. Dampen your face with lukewarm water.
3. Moisten the **brush head** with lukewarm water and apply the face wash of your choice onto the brush head.

-  Do not use a dry **brush head** on the face or neck, because it may cause skin irritation.
4. Press the **control button** to turn on the device. The device automatically starts in the low speed vibration mode.
  5. To select the desired mode, keep pressing the control button to cycle through the modes (see section "4.2 Indicator lights and function settings").
  6. Gently press the **brush head** against your skin and move it in an upward circular motion. There's no need to apply any pressure.
-  Try to avoid delicate skin areas such as around the eyes or lips.
7. To finish, rinse your skin with water, and pat dry with a towel.
-  The device is set to switch off automatically after 1 minute. However, if you change modes during a treatment, the timer will reset to zero and will start again.

### 6.1 Silk'n Pure Deluxe accessories

In addition to the handy travel case and face wash travel bottle, Silk'n Pure Deluxe includes two additional treatment accessories. Both are suitable for daily use.

The **Face Massage brush head** improves blood circulation and cleanses the skin mildly, leaving your skin fresher and brighter. Follow the steps in section "6. Get started!"

The **Eye Massage brush head** helps reduce puffiness and dark circles around the eyes. To use this accessory, follow these steps:

1. Cleanse your face and make sure the skin is clean and dry.
2. Place the metal part on the area underneath your eyes and then slide from left to right. Do this for 30 seconds to a 1 minute in one area. You can also use it with your favourite eye serum, lotion or cream.

## 7. Cleaning & maintenance

---

Use warm water up to a maximum of 45°C / 113° F.

EN

### 7.1 Cleaning the brush heads

Clean the Regular and/or Soft **brush head(s)** after each use.

- Make sure the device is switched off. Rinse with warm water and mild soap to remove any build-up that may have accumulated. Then shake off any excess water.
- Dry the **brush head** with a clean cloth. You can also dry the **brush head** by switching the device back on, and rubbing it against a clean cloth for 5-10 seconds.
- ❗ For hygienic reasons, we recommend you replace used **brush heads** every 3 months, or earlier if you see that they are worn.

⚠ Do not allow face wash to accumulate on the **brush heads**.

### 7.2 Cleaning Silk'n Pure Deluxe accessories

- Wash the **Face Massage brush head** with warm water, and remove any build-up that may have accumulated. Then dry it with a clean cloth.
- Wipe the **Eye Massage brush head** with a damp cloth using a mild detergent. Alternatively, wash it with warm water and remove any build-up that may have accumulated. Then dry it with a clean cloth.


### 7.3 Cleaning the device and the face brush stand

We recommend that you clean the device after each use, and the **face brush stand** once a month. Before cleaning:

- Disconnect the device from the **adapter**.
  - Do not use abrasive cleaners or chemicals to clean any parts.
  - Clean the device or **face brush stand** with a damp cloth using a mild detergent, then wipe it dry with a clean cloth.
- ⚠ Do not clean the device or any of its parts in the dishwasher.

## 7.4 Replacing the brush head

1. To remove the **brush head**, simply grip the outer plastic base of the **brush head** with one hand, and the device handle with the other, and pull upwards until it snaps loose.
2. To attach a **brush head**, place the center of the **brush head** on the **brush attachment point**.
3. Gently press the **brush head** and push down until it clicks into place.
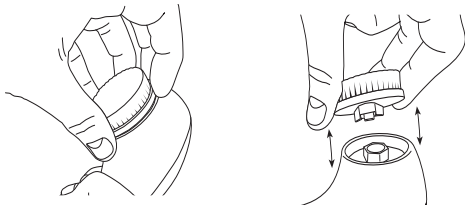
## 7.5 Storage

- You can leave the device in the **face brush stand** and use it as a counter top stand.
- If you don't plan to use the device for some time, we recommend you store it in a dry place.

## 8. Troubleshooting

Do not attempt to repair or open sealed parts of your device. You may expose yourself to dangerous components and cause serious harm to yourself. It will also void your warranty.

- ⚠ Only authorized personnel are permitted to perform repairs. For problems or queries, please contact Silk'n Customer Service.

Problem	Action
The device does not start.	Make sure the device is properly charged. Attach the <b>adapter</b> to the device and the electricity supply and leave the device to charge for at least 6 hours before the next use.
The device does not charge.	Make sure the <b>adapter</b> is connected correctly to the device and the electricity supply. Unplug the <b>adapter</b> and then reconnect it again to completely reset the device. If the problem persists, contact Customer Service.

Problem	Action
The device shuts off all the time.	Make sure the device is properly charged. Attach the <b>adapter</b> to the device and the electricity supply, and charge the device for at least 6 hours before the next use.

## 9. Parts & consumables

You can purchase replacement brushes through the retailer where you bought the device, or you can order directly from our website: [www.silkn.eu](http://www.silkn.eu).

## 10. Customer service

For more information see our website [www.silkn.eu](http://www.silkn.eu). This manual is also available as a PDF download from the website. Contact Silk'n Customer Service to report unexpected operations or events, or any other problem with your device.

Our contact details are:

Country	Service number	Email
Belgium	0800-29316	<a href="mailto:serviceconsommateurbe@silkn.eu">serviceconsommateurbe@silkn.eu</a>
France	0891-655557	<a href="mailto:serviceconsommateurfr@silkn.eu">serviceconsommateurfr@silkn.eu</a>
Germany	089 51 23 44 23	<a href="mailto:kundenservicede@silkn.eu">kundenservicede@silkn.eu</a>
Italy	+31 (0)180-330 550	<a href="mailto:servizioclientiit@silkn.eu">servizioclientiit@silkn.eu</a>
Netherlands	0900-2502217	<a href="mailto:klantenservicenl@silkn.eu">klantenservicenl@silkn.eu</a>
Spain	900 823 302	<a href="mailto:servicioalcliente@silkn.eu">servicioalcliente@silkn.eu</a>
UK	0906-2130009	<a href="mailto:customercareuk@silkn.eu">customercareuk@silkn.eu</a>
Other countries	+31 (0)180-330 550	<a href="mailto:info@silkn.eu">info@silkn.eu</a>

## 11. Disposal

Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. This ensures that the device is recycled professionally, and prevents the release of harmful substances into the environment. Make sure that the device is in an uncharged state before disposal.

## 12. Specifications

Model	H4101	
Adapter	Input:	100-240V ~ 50/60Hz 0.2A
	Output:	5Vd.c.; 500mA
	Model:	RSS1002-025050-W2E (Europe) RSS1002-025050-W3U (UK)
Powered by	Batteries 1x 3.7V; Li-ion; 500mA	
Device classification	Class II	
Protection	IPX7	

## 13. Signs & symbols

The following symbols and signal words are used in this user manual, on the device or on the packaging.

	Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.
IPX7	Protected against water immersion for up to 30 minutes at a depth of 1m.
	Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices.
	This symbol identifies electrical devices that belong to the protection class II.
	The word "GS" stands for tested safety. Products marked with this symbol meet the requirements of the German Product Safety Act (ProdSG) (EU adapter).
	Devices marked with this symbol are only suitable for household use (dry indoor rooms).

	The adapter is equipped with a short-circuit proof safety transformer
IP20	The adapter is protected against solid objects up to 12 mm.
	The adapter has a switching power supply.
	The international efficiency mark indicates that the power supply meets the level VI requirements.
	This symbol indicates the polarity of the pins in the adapter.
	The product fulfils national and international safety standards (GB adapter).

## 14. Warranty

### What Does This Warranty Cover?

Home Skinovations, its applicable affiliates (hereinafter: Home Skinovations) or the authorized dealer from whom you purchased your Silk'n product (hereinafter: the Device) shall repair or replace the Device, if any material defects caused by Home Skinovations are discovered in the Device during the warranty period, whether in the body of the Device or its operation.

### Who is entitled to warranty by virtue of this warranty certificate?

Only you, the original buyer of the Device who purchased the Device from Home Skinovations or its authorized dealers, are entitled to receive the warranty specified in this warranty certificate. In order to be entitled to the warranty service, you need to provide a proof that you are the original buyer of the Device by presenting an original receipt for purchase of the Device. This warranty is not an international warranty. The warranty will only be honored in the country in which the authorized dealer or the Home Skinovations entity from whom you purchased the Device, are located. For a list of Home Skinovations' authorized dealers see: [www.silkn.com/For\\_International\\_Visitors](http://www.silkn.com/For_International_Visitors).

**What period does this warranty certificate cover?**

The warranty covers repair requests received within one year (1) from the Device's original purchase date. Repair requests after one year (1) from the date of purchase are not covered by this certificate of warranty. Repair requests regarding a Device that was repaired or replaced must be submitted no later than: (a) one year (1) from the original purchase date of the Device or (b) ninety (90) days from the date that the repaired or replaced Device was returned to you, the later to occur. If a sticker is attached to your Device stating a warranty period longer than that indicated above, the warranty period mentioned on the sticker will apply, and the Device will be repaired or replaced as part of the warranty until such period expires.

If you purchased the Device in an EU country, you have an additional year (according to the EU-Directive 1999/44) to claim for deficiencies or malfunctioning of your Device, which were present at the time of purchase, but were not apparent then, as determined by Home Skinovations or the authorized dealer. In order to make this extended warranty effective you must contact the authorized dealer from whom you purchased the Device or the local service center and notify them about the production fault immediately, and no later than 60 days after detecting the fault.

**How to get service for the warranty?**

For further instructions how to obtain warranty service and for contact details of the service centers close to you, visit: [www.silkn.com/For\\_International\\_Visitors](http://www.silkn.com/For_International_Visitors), or contact your local authorized dealer, or call our Help Desk at:

North America: Tel: 1-877-DO-SILKN (1-877-367-4556) / Email: [contact@silkn.com](mailto:contact@silkn.com)

European countries: (Including all member states of the European Union, Switzerland, Norway and Russia) Tel: +31 (0)180-330 550 / Email: [silkn@inno-essentials.com](mailto:silkn@inno-essentials.com)

Japan: Tel: 0120-110-469 / Email: [beauty0120110469@sensepil.jp](mailto:beauty0120110469@sensepil.jp)

Israel: Tel: 1-700-700-309, Email: [info@silkn.com](mailto:info@silkn.com)

Other countries: Email: [info@silkn.com](mailto:info@silkn.com)

When asking for customer services you may be asked to fill in a form for returning the product and you shall receive a reference number. In case of a defect that is not covered by the warranty for whatever reason, Home Skinovations can charge you for the shipping and repair costs if you agreed to repair the Device at your expense.

**What is not covered by this certificate of warranty?**

This warranty does not cover any damage resulting from an accident, misuse (including using the Device other than as described in the use instructions accompanying the Device), defective use, lack of reasonable care, irregular and unreasonable use including without limitation commercial or professional use, adding parts not supplied together with the Device, losing parts, using unauthorized spare parts or exposing the Device to any electrical current other than the indicated current. The warranty will be cancelled if the Device was treated by anyone not authorized by Home Skinovations and Home Skinovations will not be liable for any damage resulting from repairs performed by anyone not authorized by Home Skinovations. The warranty does not include replacing used non-durable parts.

**WARRANTY DISCLAIMER**

AR PURPOSE. THIS WARRANTY IS THE SOLE AND EXCLUSIVE REMEDY FOR ANY DEFECT OF THE DEVICE.

**LIMITATIONS**

IN NO EVENT SHALL HOME SKINOVATIONS, ITS AFFILIATES OR AUTHORIZED DEALERS, BE LIABLE FOR ANY INDIRECT, SPECIAL, PUNITIVE, EXEMPLARY, INCIDENTAL OR CONSEQUENTIAL DAMAGES (INCLUDING, WITHOUT LIMITATION, DAMAGES FOR LOSS OF USE OR CHARGES FOR SUBSTITUTE EQUIPMENT), REGARDLESS OF THE LEGAL THEORY ON WHICH THE CLAIM IS BASED EVEN IF THEY HAVE BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. IN THE EVENT THAT THE REMEDIES SET FORTH IN THIS WARRANTY ARE HELD TO FAIL OF THEIR ESSENTIAL PURPOSE OR IF HOME SKINOVATIONS, ITS AFFILIATES OR AUTHORIZED DEALERS ARE HELD TO BE LIABLE FOR DAMAGES ON THE BASIS OF ANY LEGAL THEORY, IN NO EVENT SHALL HOME SKINOVATIONS, ITS AFFILIATES OR AUTHORIZED DEALERS, LIABILITY FOR DAMAGES EXCEED THE WHOLESALE PRICE OF A SINGLE UNIT OF THE DEVICE AS OF THE DATE THAT THE CAUSE OF ACTION ACCRUES. YOU ASSUME ALL RISK AND LIABILITY FOR LOSS, DAMAGE OR INJURY TO YOU AND YOUR PROPERTY AND TO OTHERS AND THEIR PROPERTY ARISING OUT OF THE USE OR MISUSE OF, OR INABILITY TO USE, THE DEVICE EXCEPT TO THE EXTENT CAUSED DIRECTLY BY THE NEGLIGENCE OF HOME SKINOVATIONS.

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS WHICH VARY FROM COUNTRY TO COUNTRY. SOME COUNTRIES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES OR IMPLIED WARRANTIES, SO THE ABOVE DISCLAIMERS AND LIMITATIONS MAY NOT APPLY TO YOU.

## Legal Notice

Home Skinovations Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd, Tavor Building, Shaar Yokneam P.O. Box 533, Yokneam 2069206, ISRAEL.

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Invention Works B.V.