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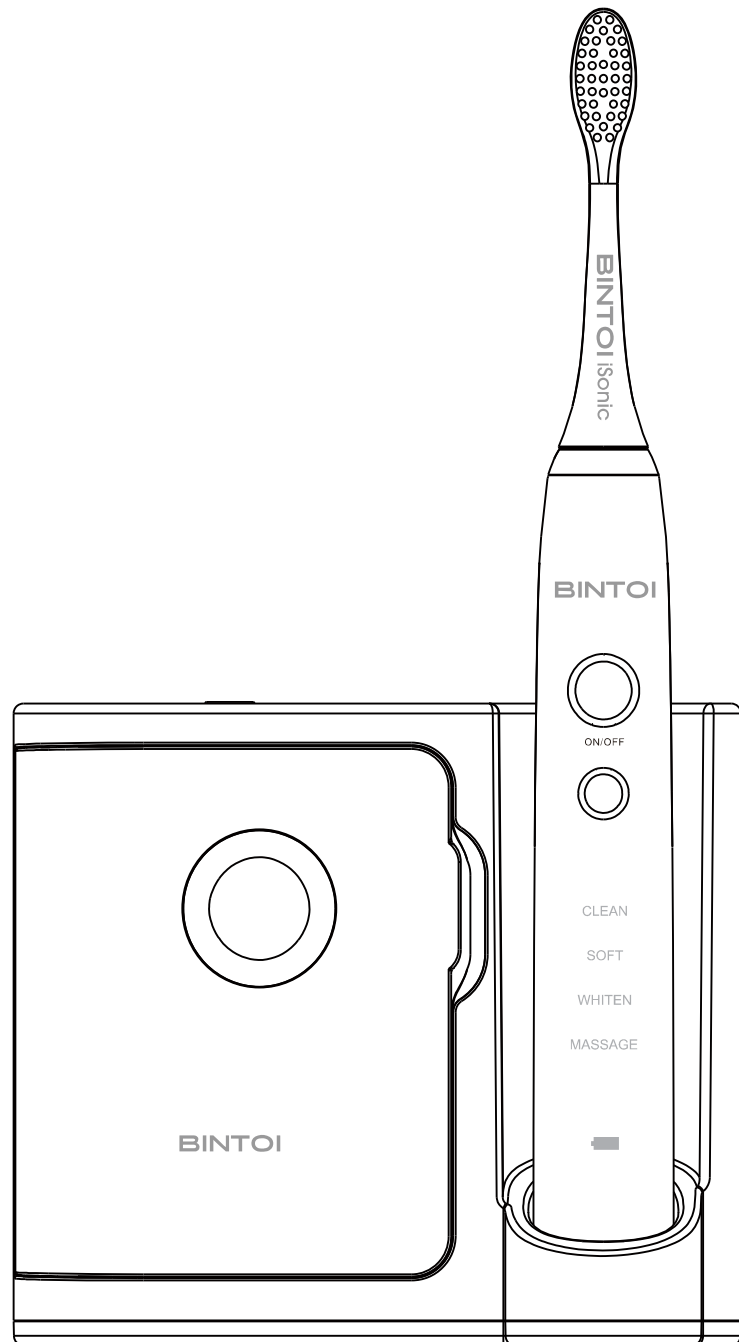
ISO9001 & ISO13485 Certificated Manufacturer

Printed in China

IMPORTANT SAFEGUARDS

The Power toothbrush is safe and effective when used according to directions. Like any electrical appliance, however, basic safety precautions should be followed especially when children are present.

READ ALL INSTRUCTIONS BEFORE USING

DANGER-WATER

To reduce the risk of burns, electrocution, fire, or injury to persons:

- ❶ Do not place or store the charger base with UV sanitizer where it can fall or be pulled into a bathtub, sink or toilet.
- ❷ Do not place or drop the charger into water or other liquid.
- ❸ Do not reach for a product that has fallen into water or other liquid. Unplug immediately.
- ❹ Do not use while bathing.
- ❺ Unplug charger base with UV sanitizer before cleaning and assure that it is dry before plugging it into an electrical outlet.

DANGER-VOLTAGE

To reduce the risk of burns, electrocution, fire, or injury to persons:

- ❶ Do not plug the charger base with UV sanitizer into an outlet with a voltage other than specified on the bottom of the charger.
- ❷ This product is designed to operate within a range of 100-240 volts.
- ❸ Voltage converters do not guarantee compatibility.
- ❹ Do not force the plug into an outlet.
- ❺ If the plug does not fit easily into the outlet, discontinue use.
- ❻ Never wrap the cord tightly around the appliance. Never allow the power supply cord to be pulled, twisted, or severely bent. Untwist the cord before use.

DANGER-WARNINGS

To reduce the risk of burns, electrocution, fire, or injury to persons:

- ❶ Children should be supervised to ensure they do not play with the product.
- ❷ Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only. Do not use any attachments or brush heads other than those recommended by the manufacturer.

- ③ The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
- ④ If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
- ⑤ Never operate this product if it has a damaged cord, plug or charger with UV sanitizer. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- ⑥ Keep the product, including power cord away from heat surfaces.
- ⑦ If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- ⑧ Do not attempt to remove the battery unless the product life is end.
- ⑨ Do not clean the toothbrush in the dishwasher.
- ⑩ Discontinue use of the charger with UV sanitizer if the UV light bulb remains on when the lid is open. UV light can be harmful to human eyes and skin. Keep the charger with UV sanitizer out of the reach of children at all times.
- ⑪ The UV light bulb is hot during or after the sanitizing cycle. Do not touch the UV light bulb when it is hot.
- ⑫ Do not look directly into the UV light.
- ⑬ Do not allow any impact to the UV lamp. If the UV lamp is cracked, chipped or damaged, do not operate. In case the UV lamp breaks, leave the place for 20-30 minutes, and then clean it, do not touch the mercury directly with bare hand especially for pregnant women.
- ⑭ For the best disinfection performance, always keep the UV lamp clean and avoid stains.
- ⑮ In case the UV lamp gets dirty, wipe it with a soft cloth very gently.
- ⑯ Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- ⑰ Do not leave appliance when plugged in. Unplug from outlet when not in use and before servicing.
- ⑱ Do not pull or carry by cord, use cord as a handle, close a door on cord, or pull cord around sharp edges or corners. Do not run appliance over cord.
- ⑲ Do not handle plug or appliance with wet hands.
- ⑳ Turn off all controls before unplugging.

- 21 Replacement of the UV lamp can only be carried out by professional staff, can not be operated by the user.
- 22 Professional maintenance personnel have excessive UV exposure risk.
- 23 This cover is provided with an interlock to reduce the risk of excessive ultraviolet radiation. Do not defeat its purpose or attempt to service without removing.
- 24 This appliance complies with US Federal requirements, including 21 CFR, Chapter 1, Subchapter J, Radiological Health.
- 25 Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- 26 If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- 27 The battery must be removed from the appliance before it is scrapped.
- 28 The appliance must be disconnected from the supply mains when removing the battery.
- 29 The battery is to be disposed of safely.
- 30 Regarding how to remove the battery, please refer to last section <disposal>.

MEDICAL WARNINGS


- 01 If you have had oral surgery or gum surgery within the past two months, consult your dental provider before using the product.
- 02 Consult your dental provider if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- 03 Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
- 04 Consult your physician prior to using the toothbrush if you have any medical concerns.
- 05 The toothbrush has been tested and is certified as per ETL and CE regulations.

- ⑥ Always unplug this product immediately after using.
- ⑦ Never use while sleeping or drowsy.


SAVE THESE INSTRUCTIONS

Configuration

Handle

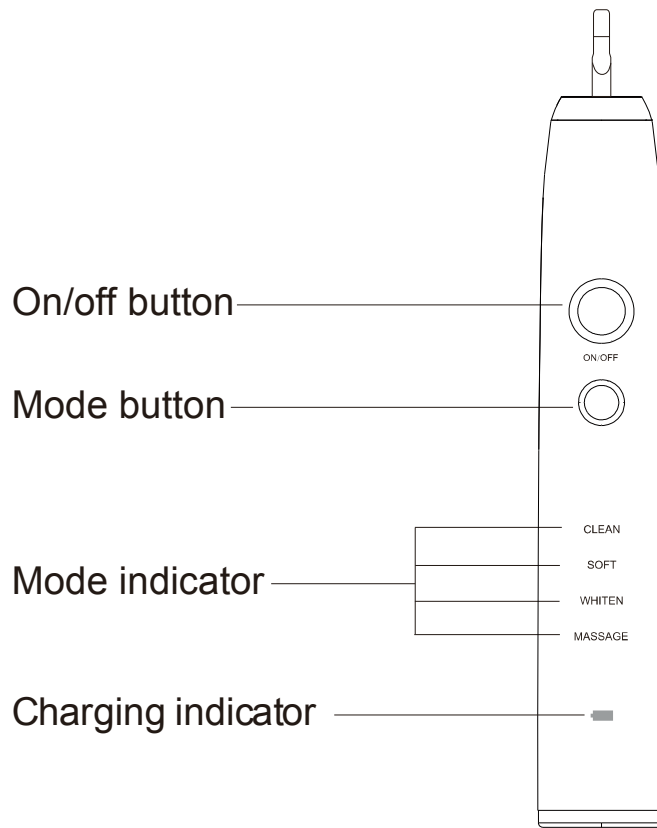Charger

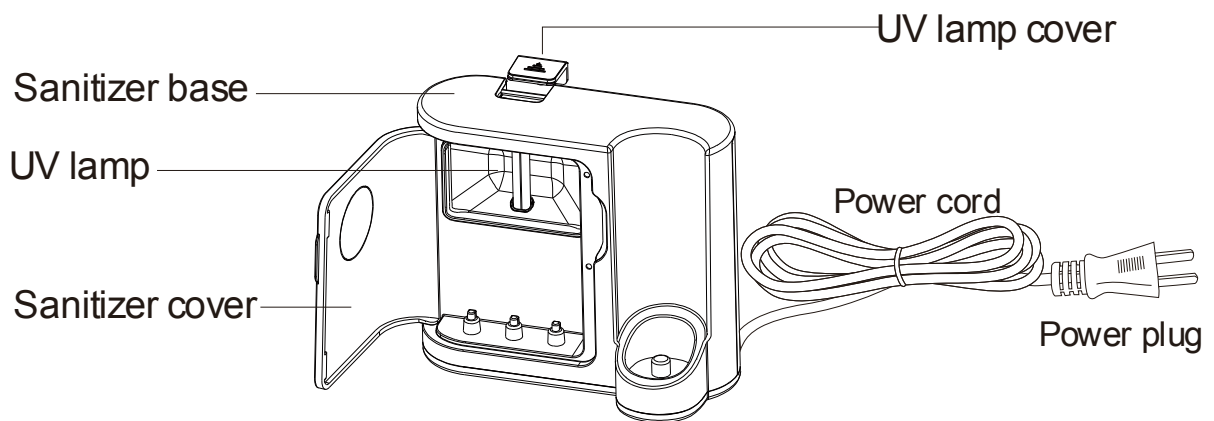
Brush heads

Getting Started


Charging the toothbrush :

Be sure to turn off the toothbrush before charging. Plug the power plug into the outlet. Place the handle in the charger. The charger indicator light will be lit, indicating the handle is charging. After the battery is fully charged, charge indicator will turn solid green, indicating the toothbrush is fully charged.

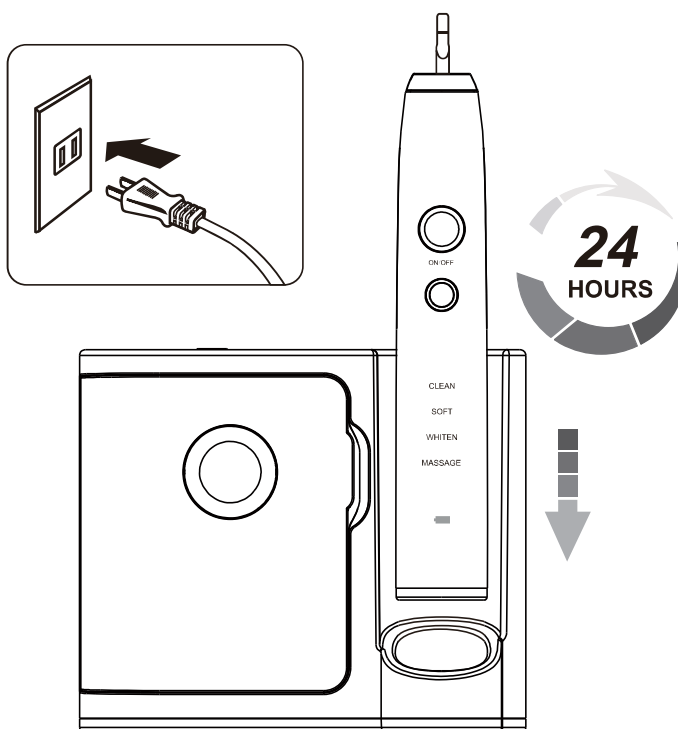- ❶ It takes at least 24 hours to charge the battery fully.
- ❷ The toothbrush must be charged for 24 hours, if:
 - a. this is the first time using the toothbrush;
 - b. the toothbrush has not been used for more than three months;
- ❸ To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.
- ❹ During charging, the toothbrush will get warm, but this is normal.

The operating duration will become shorter as the battery gets older. If the toothbrush operates only for 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact service center for support.



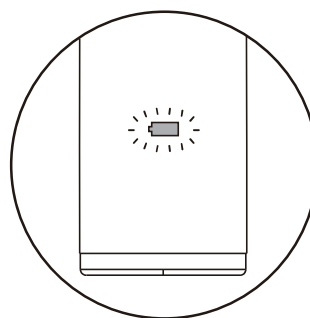
LED charging indicator will be Red during charging; after full charge LED indicator will turn to green

How To Use

- Slide the brush head onto the handle. Push brush head down until it locks into place on handle.
Note: There is small gap between color ring and the handle. (Fig 1)
Press the mode button to select mode. If you know the mode used last time, you need not select mode again. There are 4 modes.
 - CLEAN: Standard mode for superior teeth cleaning.
 - SOFT: Gently caress, carefully clean sensitive mouth.
 - WHITEN: Teeth cleaning with polishing and whitening.
 - MASSAGE: Pulsates for gentle gum stimulation.
- Wet the bristles and apply a small amount of toothpaste.
- Place the bristles against the teeth at a slight angle towards the gum line.(Fig 2)
- Press the on/off button to turn on the toothbrush.

Notes:

- During brushing, if press Mode button, toothbrush will work at next mode and start another 2 minutes cycle and the mode indicator will light up.
- Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.

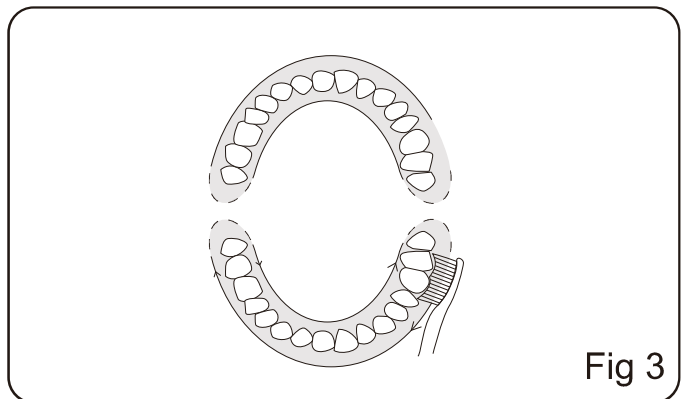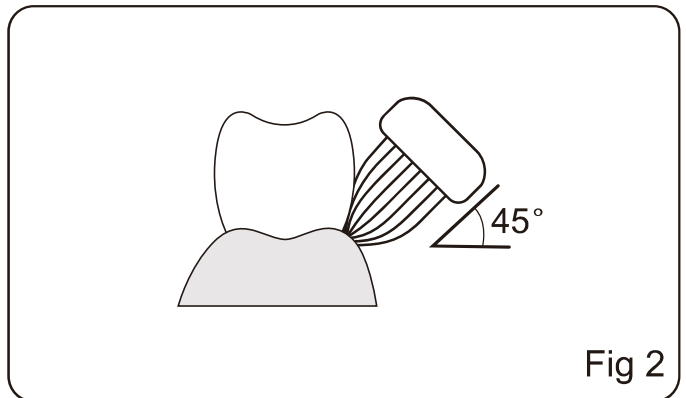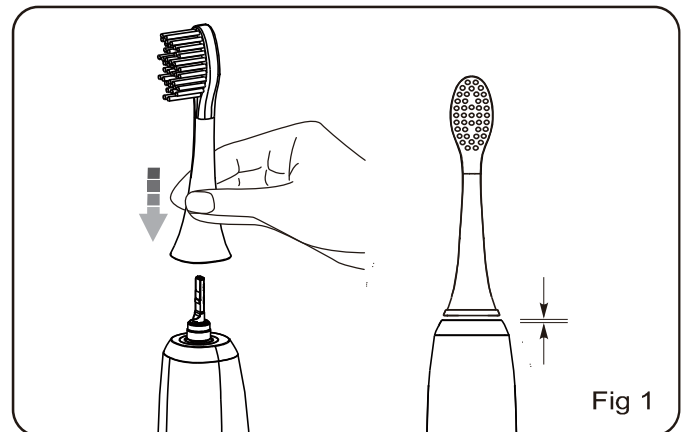- 7 In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds.
- 8 The brushing action will pause between each 30 second interval, indicating when to move on to the next section.
- 9 The toothbrush will automatically stop after 2 minutes of use. (Fig 3)
- 10 Pull brush head out from the toothbrush.(Fig 4)
- 11 Rinse both the brush head and handle under running water wipe them dry. (Fig 5)
- 12 Put the handle back into the charger.

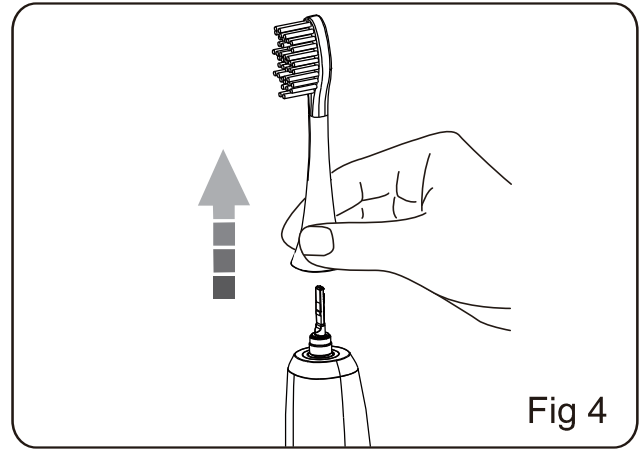
Fig 4


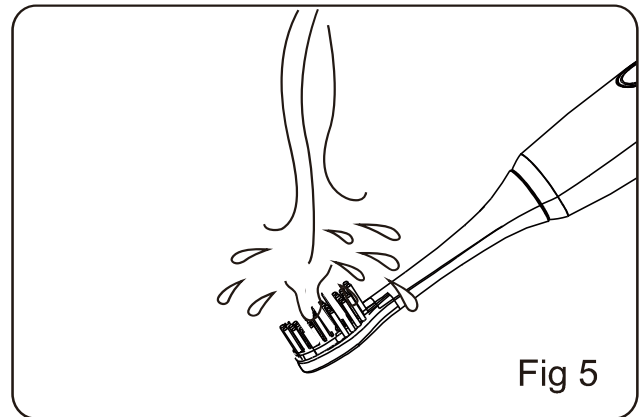
Fig 5

Sanitizing the Brush Head

This product includes a UV sanitizer to sanitize the brush heads between use. The sanitizer has been proven 99,9% effective against E.Coli, Candida Albicans, Salmonella and H. Simplex. It is also 96,8% effective against H1N1.

- 1 Gently shake the toothbrush to remove any excess water after cleaning.
- 2 Open the sanitizer cover.
- 3 Place the brush heads on the stand with the bristles directly facing the UV lamp.
- 4 Make sure the charger with UV sanitizer is plugged into a working outlet.

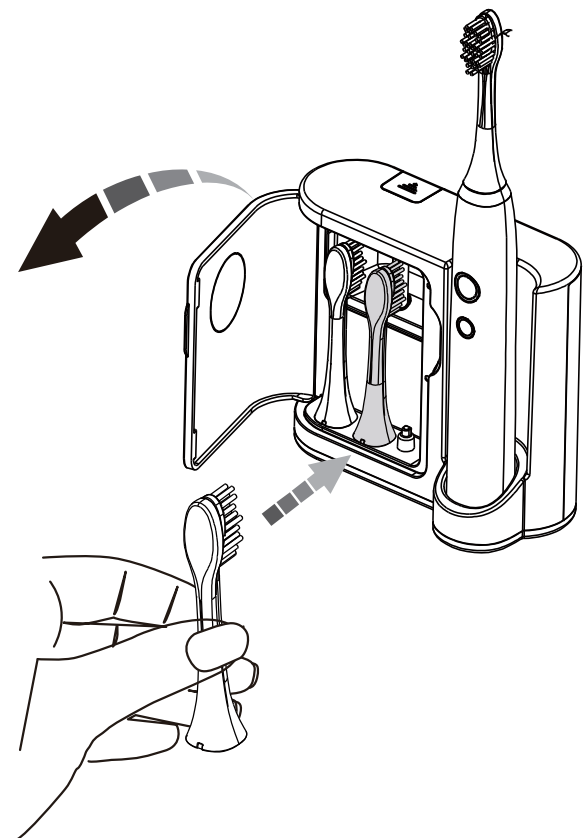Maintenance

- ① Remove the brush head.
- ② Wipe all surfaces areas with a damp cloth and place the handle back into the charger. Do not submerge the handle into water.
- ③ Rinse the brush head and bristles after each use.
- ④ Do not clean in the dishwasher.

Troubleshooting

Problem	Possible Cause	Solution
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the handle for at least 24 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the handle for at least 24 hours.
UV sanitizer light does not turn on.	The sanitizer cover may be open.	Close the sanitizer cover.
	The power plug may not be plugged in.	Plug in the power plug.
	The UV lamp may be defective.	Contact service center for support.

Disposal

The power source of this product is Lithium-ion battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. Do not throw away the appliance or the UV light lamp (lamp contains mercury) with normal household waste at the end its lift. Instead, hand it in at an official waste collection point for recycling. By doing above all, you help to preserve and protect the environment. If you have recycling questions, contact your local waste management office.




Removing and recycling the battery

Do not disassemble the product at any time other than when disposing of it.

Doing

so earlier will destroy the product. Observe basic safety precautions.

- 1 Before removing the battery, be sure to completely drain the unit of power. Please turn on th device by press the on/off button, and let it work until it does not work.
- 2 Use a screwdriver to remove the screw and battery door.
- 3 Hold the handle upside down, press the toothbrush shaft, and remove the two hooks at the same time to release the internal components of the handle.
- 4 Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.
- 5 Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

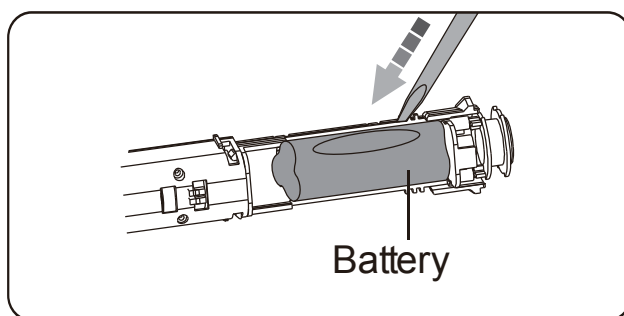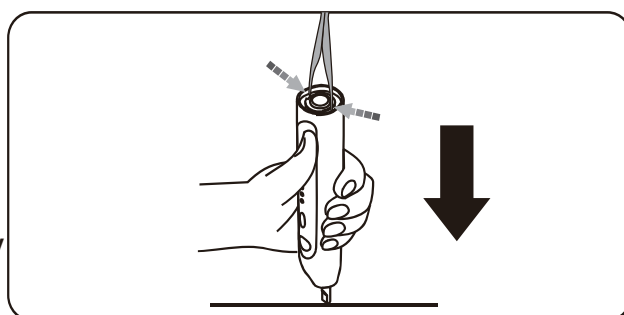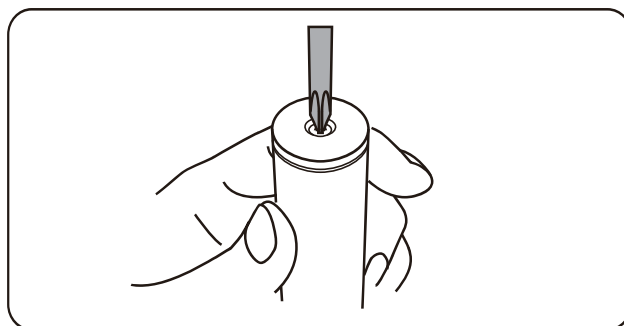Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.

FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Gaurantee

Your toothbrush is guaranteed for 1 year against manufacturers' defects for the original purchaser only, from date of purchase. The warranty does not apply to damage caused by improper handling, accidents, professional use, not following the operating instructions or alterations made to the instrument by third parties.



Warranty only applies to the main device. All other accessories are not covered by warranty.

There are no user serviceable parts inside.

Travel Charging Case with UV Sanitizer

Instructions Manual

Model: D600

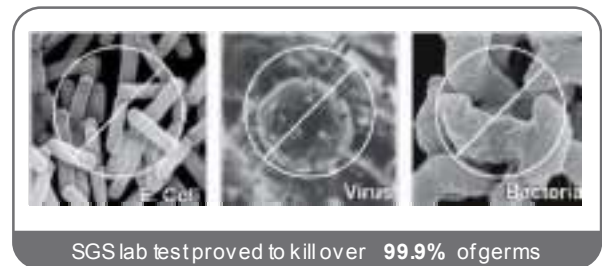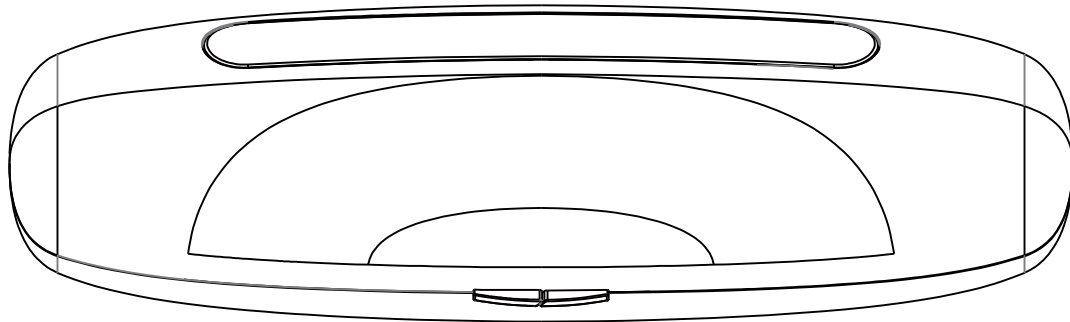IMPORTANT SAFEGUARDS

Please read the safety instructions before using this product. The safety instructions provided here are to help you use this product safely and correctly, and also to prevent you from injuring yourself or others. When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

DANGERS-WATER

- Do not place or store the sanitizer where it can fall or be pulled into a bathtub, sink or toilet.
- Do not place or drop in water or other liquid.
- Do not reach for the sanitizer that has fallen into water or other liquid. Unplug immediately.
- Do not use while bathing.
- Unplug sanitizer before cleaning and assure that it is dry before plugging it into an electrical outlet.

DANGERS-VOLTAGE

- The appliance is only to be used with power supply with provided with the appliance, the adapter with 100-240V/AC input and DC5V 600mA is allowed, and also can be used with DC5V power supply.
- Voltage converters do not guarantee compatibility.
- Do not force the plug into an outlet.
- If the plug does not fit easily into the outlet, discontinue use.
- Always unplug this product immediately after using.
- Never use while sleeping or drowsy.
- The supply cord cannot be replaced. If the cord is damaged the appliance should be scrapped.

WARNINGS

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Warning – This appliance contains a UV emitter. Do not stare at the light source.
- Voltage converters and plug adapters do not guarantee voltage compatibility.
- Never force the plug into an outlet; if it does not easily fit into the outlet, discontinue use.
- Unplug the sanitizer and call customer service if the sanitizer gives off smoke or a burning smell operation.
- Never use the sanitizer if it has a damaged cord or plug. Keep away from heated surface. Do not use sanitizer outdoors or in areas with direct sunlight.
- Discontinue use of the sanitizer if the UV light bulb remains on when the door is open or if the window the door is broken or missing. UV light can be harmful to the human eye and skin.
- Do not operate the sanitizer without the protective screen in place to avoid contact with a hot bulb.
- Do not use attachments other than those recommended by the manufacturer.
- Do not use if dropped into water.
- Do not use this product or any of its components for purposes other described in this manual.
- Do not operate the UV sanitizer where aerosol (spray) products are being used or where oxygen is being administered.


SAVE THESE INSTRUCTIONS

PARTS IDENTIFICATION

Transparent Lens

Charge LED indicator



USB Port

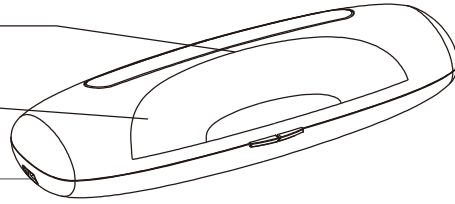
Toothbrush Handle

Removable protective screen for lightbulb


Toothbrush Head

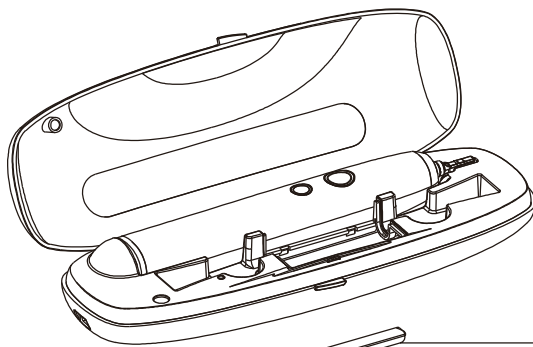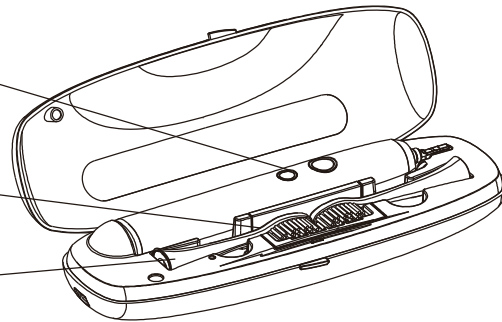UV light bulb

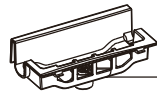
Removable Protective Screen for light Bulb

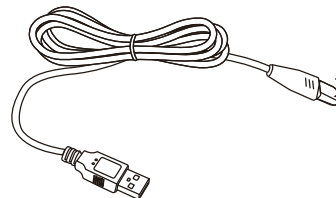USB Cable


The picture is for reference only, please make the object as the standard. This product does not include toothbrush handle and heads.

HOW TO SANITIZER

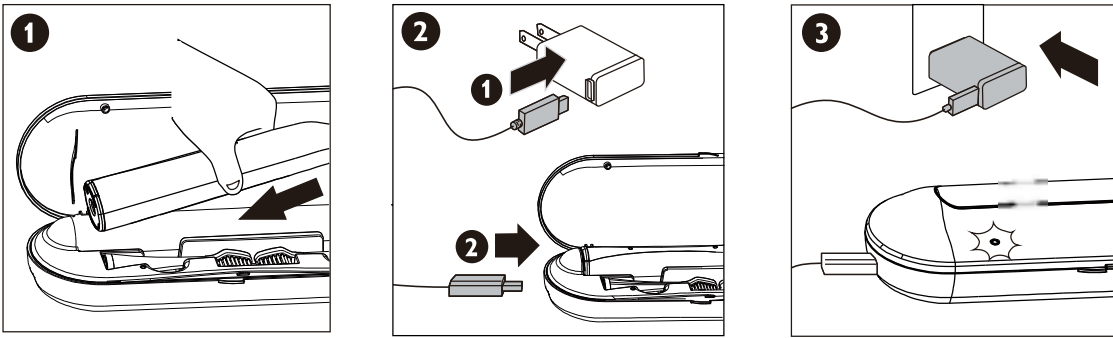
- Place the toothbrush head and handle in the travel case, please ensure the bristles facing the UV light bulb. (Fig.1)
- Put the Micro USB plug into the USB port of case, and another USB plug into the adapter port or other USB power sources (Fig.2)
- Close the top lid and plug the adapter in electrical outlet, UV clean indicator lighten, sanitizing after run 10 minutes UV clean cycle will be automatically shuts off (Fig.3)



HOW TO CHARGING

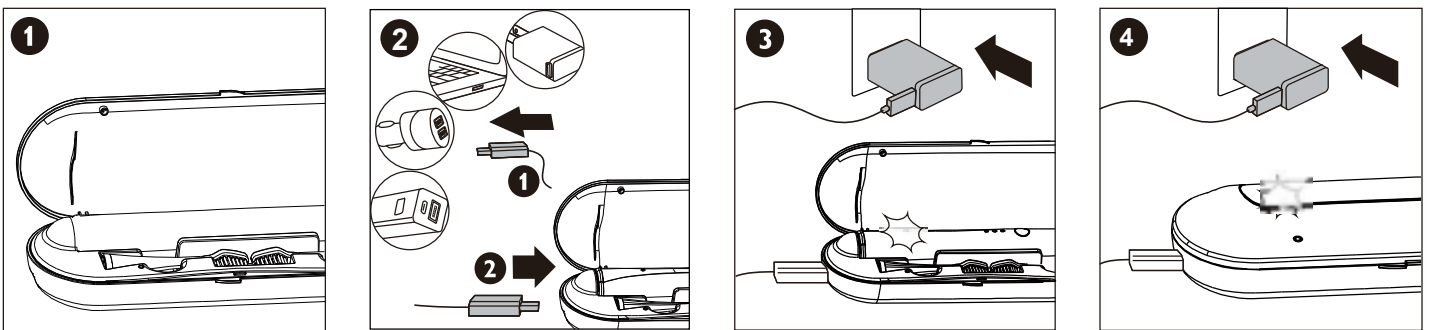- Place the handle in the charging travel case. (Fig.1)
- Plug the USB Cord into the travel case plug and another plug into the adaptor or other USB power sources. (Fig.2)
- Plug the adapter into electrical outlet. The toothbrush charge indicator on the handle will be flashing. (Fig.3)
- Closed the top lid, the charge indicator on the handle flashes red until fully charged. (Fig.4)

REPLACE BULB

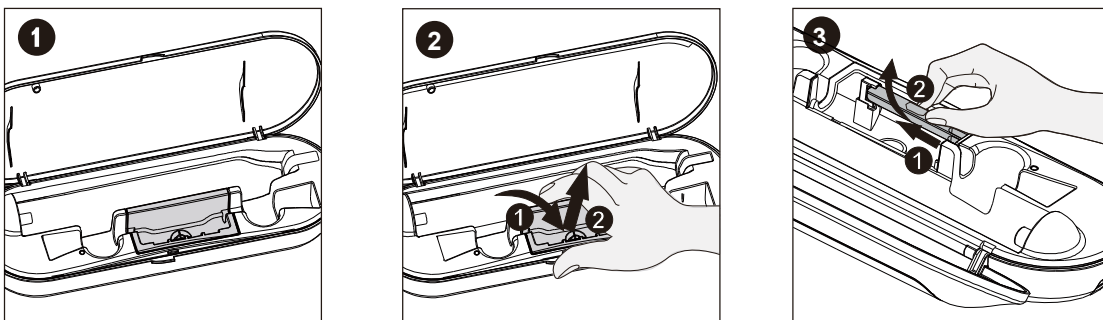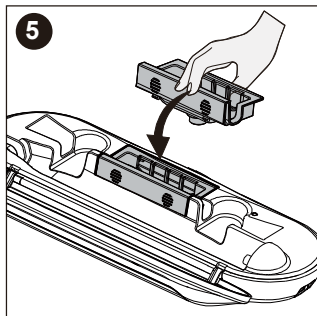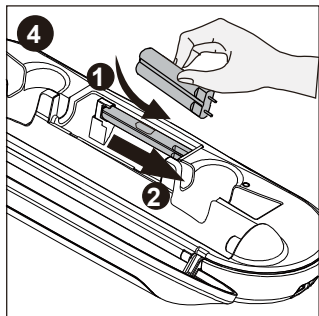- Please take out the toothbrush handle and heads, and place the travel case to best replace direction. (Fig.1)
- Remove the protective screen by gently pulling it forward (1) and lifting it up(2). (Fig.2)
- Remove the UV light bulb by gently pulling it forward(1) and lifting it up (2). (Fig.3)
- **Note:** Do not crush or break the bulb. UV light bulbs contain a small amount of mercury.
- Install the UV light bulb by lifting it down (1)and gently pulling it back(2). (Fig.4)
- **Note:** Please ensure the bulb pins tightness connect the holder base.
- Install the protective screen by lifting it up and lock on the inner tray. (Fig.5)

TROUBLESHOOTING

Problem	Possible Cause	Solution
The charging function doesn't work.	You've just purchased the toothbrush or it has been left unused for more than 3 months.	Charge the case unit for at least 24 hours.
The UV light does not go on.	The lid does not closed.	Make sure the lid is properly closed.
	The UV light bulb has reached the end of its life cycle	Replace the UV bulb.
It only works for a few minutes even after charging.	The toothbrush battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient	Charge the main unit for at least 24 hours.

REPLACEMENT

Replacement parts are available through your customer service enter. Refer to the guarantee and service section of this manual for phone number and website information. The parts available for ordering a UV light bulb & adapter

STORAGE

If you are not going to use the sanitizer for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight. You can use the cord-wrap feature to store the power cord in an organized manner.

ENVIRONMENT

The UV light bulb contains mercury and may be harmful to the environment. Do not crush, break or attempt to disassemble the bulb. Do not throw away the UV light bulb with the normal household waste when you replace it, but hand it in at an official collection point for recycling. Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment. If you have any recycling questions, contact your local waste management office.

FCC STATEMENT

This device complies with Part 15 of the FCC Rules.

Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISO9001 & ISO13485 Certificated Manufacturer

Printed in China