


OPERATION MANUAL

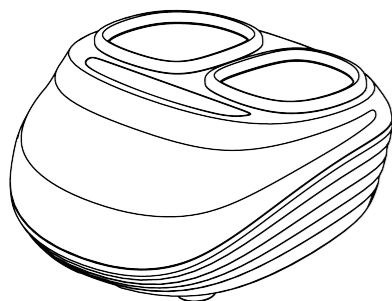

Due to possible design changes, the specific appearance and parameters are mainly in kind

Safety guidelines

1. Description of important matters

- In order to ensure safe and correct use of the product, please be sure to check carefully before use Operation Manual.
- Be sure to use this product in accordance with the product usage in this operation manual.
- Unplug the power plug from the outlet after use or before cleaning the product.
- The device has a heating function. Care must be taken when used by persons who are not sensitive to heat.
- Do not use on sluggish and overly sensitive skin or those with poor circulation part.
- Do not allow children or persons with physical and mental disabilities to use the warming function on their own. Improper handling can be dangerous.
- To avoid fire, electric shock, or injury due to overheating, do not Operate the product under a pillow or blanket.
- Do not cover the use, so as not to cause the motor to overheat and cause malfunction or damage.
- The use of products by people with physical and mental disabilities and disabilities requires the guidance of others and can bear the responsibility responsible for its safety.
- If you feel uncomfortable during the massage, please stop the massage immediately.
- Do not use this product under the following conditions: The power cord or power plug is damaged, malfunctioning, dropped to the ground or damaged, dropped into water, Please contact the service center for inspection and repair.
- Do not use spare parts not recommended by the manufacturer.
- Do not insert or drop objects into any openings in the product.
- If you want to disconnect the power, turn off all operations, and then unplug the power plug from the outlet.
- This product is not intended for self-diagnosis or treatment. If diagnosis or treatment is required, it should be handled by professional medical personnel.

2. The use of the environment

- Do not use the product in a humid or dusty environment to avoid Fault.
- Do not use this product in an environment where the room temperature is higher than 40°C.
- Do not place the product outdoors.
- Do not place the product near electric heaters, near stoves, or in direct sunlight.

3. People who are not suitable to use this product

If you have the following conditions, please consult your doctor before use.

- Those who are being treated by a doctor or who are feeling unwell.
- Heart disease, cerebral arteriosclerosis, cerebral thrombosis.
- Patients with malignant tumor.
- Pregnant and menstruating women.
- Suffer from osteoporosis.
- Those using a pacemaker or implanted medical device.
- Patients with acute illness.
- Minors.

1

4. Pay attention

- Do not expose the control switch to water or liquids.
- Do not operate the control switch with wet hands.
- Do not fall asleep while using the product.
- Do not use this product while intoxicated or unwell.
- Do not use if the cloth cover of the product is damaged.
- Do not use this product for more than 30 minutes continuously.
- Please check the product for damage before use.
- Please do not use a voltage other than the rated voltage.
- Do not use one hour before or after eating to avoid adverse effects.
- When using the product in a low room temperature, do not use the warming function suddenly. It is recommended to gradually adjust the temperature to the normal range before use.
- If the product is moved from cold storage to a warm environment for use, it may Due to the accumulation of water vapor in the parts, the normal operation of the product is affected, and even the hardware failure. It is recommended to place the product at normal room temperature for 1 hour after moving Start using it again.


Warranty Rights

1. This warranty is not transferable.
2. The effective warranty period of this product is one year.
3. When requesting warranty service, please present the original purchase certificate and warranty card.
4. The purchase time is subject to the valid documents issued by the counter.
5. From the date of purchase by the user, you can enjoy a full warranty during the warranty period (the appearance is not covered by the warranty).
6. This warranty does not apply to leases, transactions, commercial-institutions and other residents user.
7. The following situations do not belong to the scope of warranty, and only provide paid maintenance services (charge labor and material costs) user.

- a) All damage caused by human factors, including use under abnormal working environment, do not press Instructions for use and storage, damage caused by improper use;
- b) The user disassembles the machine without permission, repairs or modifies it by himself, or does not have the authorization of the company Three-party dismantling, modification and maintenance;
- c) Damage caused by poor transportation after purchase;
- b) Due to other force majeure factors (such as floods, fires, lightning strikes and earthquakes, etc.) damage caused.

2


Product composition

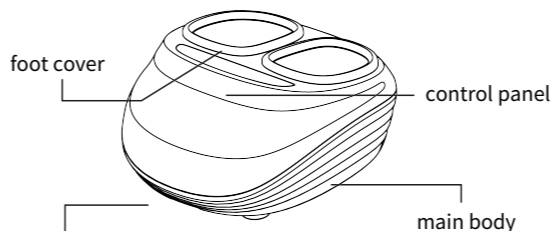

Adapter interface

1. Central surface: 7 independent control buttons - control the massage function of the product.
2. The center of the front bottom: scraping massage head - realize the acupuncture points in the arched area of the foot massage, heated.
3. Front and bottom of the front: rolling massage head - to achieve acupoint massage in the soles of the feet.
4. Front, bottom and rear end: top rubbing massage head - to achieve acupoint massage in the heel area.
5. Inside and outside the surface: airbag - close to the instep, providing air pressure massage.

1. Operating Instructions for Buttons

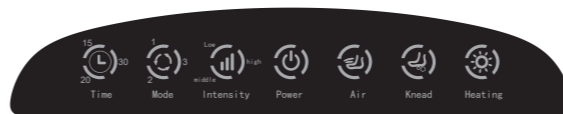

Power:

In standby state, click this key to turn on the power, it will enter automatic mode 1 by default, The working time is 15 minutes, the heating function is turned on; then click the shutdown to turn off all functions.

Time:

In the power-on state, click to cycle to select "15-20-30" minutes of working time.

Mode: In the power-on state, click to cycle through the selection mode "1-2-3"

Mode 1: Forward and reverse scraping + airbag operation.

Mode 2: Forward rotation scraping + airbag operation.

Mode 3: Reverse Gua Sha + Airbag Operation.

Intensity: in the working state of air pressure, click to switch the air pressure intensity "low one medium one high"

Air: In the power-on state, click to enter the pure air pressure working mode, the sole scraping function is turned off in this mode.

Knead: in the power-on state, click to enter the sole scraping mode, the air bag function is turned off in this mode.

Heating: In the power on state, tap the switch heating function, the default high-grade heating The indicator light is always on: then press the low-grade heating prompt light flashes

3

2. Remote control operation instructions

- ⏻ Power : Control the opening and closing of the function of the whole machine.
- 🕒 Time : Set the run time range to 15/20/30 minutes.
- ⌚ Mode : Set 1/2/3 mode.
- 🔊 Intensity : Set the intensity with L/M/H (low to high).
- 🌬️ Air : Press this button to activate the single air pressure function and the scraping stops working.
- 👉 Knead : Pressing this button activates the single foot massage function and the air pressure stops working.
- 🔥 Heating : In the power on state, tap the switch heating function, the default high-grade heating The indicator light is always on: then press the low-grade heating prompt light flashes

Troubleshooting

If you encounter problems when using this product, please refer to the following tips and possible solutions decision method.

If this product is still in a faulty state, please send the product to a service center or inspection and repair by authorized technicians.

Abnormal situation	Possible reason	Approach
Not working	1. Not powered on. 2. The function keys are not turned on.	1. Connect to the power supply or check whether the power supply is well connected. 2. Turn on the on/off button.
Sudden stop of operation	1. Power off. 2. The massage time exceeds 30 minutes. 3. Continuous use time is too long, Product overheating protection. 4. If the force is too large, press ball lock.	1. Check whether the power interface is connected well. 2. Re-open the on/off button. 3. Re-open the on/off button after the product cools down. 4. Re-open the on/off button.

Product keeper

1. Storage

- When not in use, please put the product in the color box and place it in a safe, dry, Cool and ventilated place.
- Do not twist the switch power supply and wires.
- Do not use switching power supplies and wires to hang the product.
- Avoid sharp objects scratching the surface of the product.
- Do not place the product in direct sunlight or high temperature.

2. Cleaning

- Before cleaning, make sure the main power switch is turned off and unplug it.
- Do not use switching power supplies and wires to hang the product.
- Never use aggressive cleaning agents; such as benzene or thinners to clean the product.
- Please use a dry cloth to clean the surface of the product.

4

3. Maintenance

- Keep the product clean and away from high humidity locations.
- If you do not use this product for a long time, please use a dust cover (please buy it yourself) Cover the product to prevent dust from falling into it.

Technical parameter

Product name	Foot massager	Product number	JH-ZL1
Input voltage	110V-220V~	Rated voltage	12V
rated power	24W		
Executive standard	GB4706.1-2005/GB4706.10-2008		

Warranty card

Thank you for purchasing this health product, we will sincerely follow the "Consumer Protection Law of the People's Republic of China" to provide you with services, for your legitimate rights and interests, please be sure to fill in the following content in detail, and send the receipt back to our company.

Purchase Date:	year	month	day	Where to buy:	
Customer name:				gender:	
contact address:				Telephone:	
Product name and model:				Salesperson's signature:	
Seller's phone number:				After-sales phone:	

Certificate

product name: Foot massager

quantity: 1 Person

Date of manufacture: _____

Inspectors: _____



This product has been inspected to meet the specified quality standards, Grant factory.

Warranty Regulations

Dear Customer:
Thanks to Wellcome's products, in order to protect your legitimate rights and interests, Eliminate your worries, and at the same time improve the after-sales service for customers, This warranty is specially formulated, please read it carefully and provide your valuable advice and suggestions.

1. The company's products enjoy a 12-month free warranty from the date of sale. Services, such as the product has a warranty period of more than 12 months from the date of sale. The company will charge the cost of parts for the current product failure as appropriate.
2. Please properly package and transport the repaired product, if it is damaged or lost during transportation, The company is not responsible.
3. During the free maintenance period, if the following situations occur, the company has the right to refuse service or to charge materials and service fees as appropriate.
 1. Product failure or damage caused by improper use by the user;
 2. Impulse, heating, battery, water in the controller or unauthorized installation. When the accident results in burning;
 3. Damage caused by dismantling and repairing by non-maintenance service providers;
 4. Please read this product manual carefully before installing and using the company's products.
 5. The warranty card must be stamped and dated by the sales unit to ensure that your rights.

5