

Fairywill

Sonic toothbrush
P11

User Manual


Thank you for choosing the Fairywill P11 sonic toothbrush. Please read the instructions carefully before using the product and keep them in a safe place for future reference.

1. Product description

- Waterproof and non-slip design.
- The toothbrush has interchangeable brush heads.
- High-frequency vibration provides 62,000 brush strokes per cycle* to effectively remove plaque, prevent gingivitis and deep clean your mouth.
- A USB cable is included for charging the product.
- The toothbrush has a low battery reminder feature. When the charge level is low, the indicator will flash slowly or the toothbrush will not turn on.
- The smart timer will remind you every 30 seconds to change your oral zone cleaning.
- There are 3 modes available; Clean, Soft and White. There is also a mode memory feature built into the toothbrush, which means that when you turn the toothbrush on, the default initial mode remains the same as the mode you set when you last turned the toothbrush off.

2. Product design

- ① Toothbrush head
- ② Power button/operation mode button
- ③ Function (mode)/ Display power/charging
- ④ Handle
- ⑤ Charging port

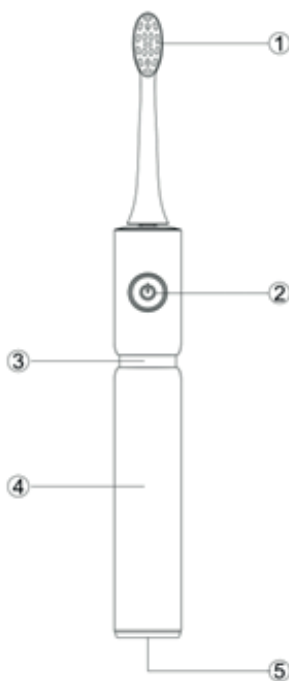

3. Product operation

Usage

- Slide the brush head down into the handle and make sure it is locked in the correct position.
- Turn the toothbrush on by pressing the power button.
- Press the power button repeatedly to switch modes ("Clean", "Soft" and "White").
- Wet the bristles and apply enough toothpaste to the brush.
- Place the toothpaste bristles against your teeth at a slight angle toward your
toward the gumline.
- Gently glide the brush head over the teeth in the front, back, and gums. The brush should glide over the teeth with minimal pressure. Do not apply too much pressure to the teeth/gums or apply the bristles to one area of the mouth for an extended period of time.
- The toothbrush has a 2-minute smart timer. Two minutes is the dentist-recommended optimal time to clean your teeth. There is a pause after every 30 seconds to remind you to move on to another oral zone.
- Move to the next zone when the pause occurs. The toothbrush will automatically turn off after 2 minutes.
- To avoid toothpaste stains, turn the toothbrush on only after you put it to your teeth and remove it from your mouth only after you turn it off.

Replacing the toothbrush head

Assembly: Slide the brush head down into the holder. Make sure it is in the correct position.

Removing: Pull the brush head out of the holder.

Product charging

- If the toothbrush has a low battery when you first start it up, fully charge the unit.
- If the toothbrush has not been used for more than 3 months, fully charge it according to the charging instructions before use.
- If the charging indicator blinks slowly during operation, it means the toothbrush has a low battery.
- To charge the toothbrush, insert the USB cable into the charging port on the bottom of the toothbrush and make sure it is securely connected to the power supply.

- The white charging indicator will flash slowly during the charging process. When the battery is fully charged, the indicator will light up for two minutes and turn off automatically after 2 minutes. Once the toothbrush is fully charged, you can use it twice a day (2 minutes each) for 30 days.
- The toothbrush has a waterproof feature, so you can safely use it in the bathroom.

Operation modes

Clean mode

This mode gives your toothbrush approximately 62,000 brushing strokes per cycle*. The three white indicators light up, the brushing motion is rapid and the vibration intensity level increases. This mode is appropriate for a healthy mouth and offers the most powerful cleaning action.

Soft mode

In this mode, the brush makes about 56,000 brushing strokes per cycle*. Only one white light is illuminated. The brushing motion is gentle and the vibration intensity is the lowest of all modes. Soft mode is suitable for users with sensitive teeth.

White mode

This mode uses 56,000 to 62,000 brush strokes per cycle*. Two white indicators light up, and the brush strokes and vibration intensity change rapidly. This mode is designed for gentle brushing to maintain healthy gums.

*A brushing cycle of 2 minutes is recommended by dentists.

4. Specification

Working temperature: 0-50°C

Operating voltage: DC 3.7 V

Bristle material: Nylon

Input voltage: DC 5.0 V

Rated power: 0.7-1.8 W

Charging time: Approximately 2 hours

Waterproof: IPX7

5. Precautions

- Do not disassemble, repair or alter the product without authorization.
- Persons with serious oral diseases or after oral surgery should consult a physician before use.
- Children, infants and persons incapable of handling the product should be supervised by a parent or adult.
- Do not soak the product in water for an extended period of time.
- Do not place the toothbrush in high temperatures or expose it to sunlight.
- Do not press bristles firmly against teeth or gums.
- This product is not suitable for use by children under 12 years of age.