auronic

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Introduction

- Thank you for choosing our Auronic product! Please always follow the instructions for safe usage. If you have any questions about the product or experience any problems, please contact us by email: service@auronic.nl
- The Heating Pad Vest can be used for short-term heating of your back, neck and shoulders. The Vest is only suitable for household use and not suitable for commercial use.

Warnings and Safety Instructions

 Please ensure that the packaging is kept out of the reach of children! Risk of suffocation!

 Make sure to keep this user manual. Read the manual carefully before using your Auronic Heating Pad Vest!

General

- Be aware of the risks and consequences resulting from misuse of the Heating Pad Vest. Misuse can result in personal injury, or destruction of the vest.
- Use the Heating Pad Vest only for purposes described in this user manual.
- If the Heating Pad Vest malfunctions, or appears to be damaged in any manner, do NOT use the Heating Pad Vest and contact Customer Service.
- The Heating Pad Vest is designed and intended for private, non-commercial use only. The vest is not intended for use in a hospital setting.
- When you are not using the vest, turn the device off and disconnect the device from the power socket.
- Prolonged use of the Heating Pad Vest at a high setting could lead to skin burns.
- Only setting 1 can be used for continuous operation.

Persons and Use

- This device should NOT be used by children under the age of 3 or persons with physical, sensory or mental disabilities, or persons with lack of experience with and knowledge of the device, unless they are supervised or have been given instruction concerning the use by a person responsible for their safety.
- Children from the ages 3 to 8 should only use the vest when they are supervised by a person responsible for operating the controller and who is responsible for their safety.
- Keep the Heating Pad Vest away from children under the age of 3 when using and storing. Make sure children do not play with the vest.
- This vest must not be used by people who are not sensitive to heat or by people who may not be able to react to overheating (diabetics, people with skin conditions or scars, people who are sleeping or unconscious, after drinking alcohol/ using drugs or intoxicatives, or after taking medicine).
- The Heating Pad Vest should only be used on the neck, shoulders and back. Do NOT place the Vest near our eyes or other sensitive areas.
- Never place the vest on or over your head.
- Prolonged use of the vest at a high setting could lead to skin burns.
- Do not use the vest on animals.

Medical remarks

- This device is NOT a substitute for treatments administered by doctors or physiotherapists.
- Do NOT use the Heating Pad Vest/consult your physician if you: Are pregnant; Are fitted with a pacemaker or other medical devices. Have circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins, or other skin conditions; If you are sensitive to heat.
- Before using the vest, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

Electronics

- Do NOT reconstruct, repair or disassemble the Heating Pad Vest. This can damage the device. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the vest to the service centre for repairs.
- Make sure your outlet supplies the correct power voltage to avoid damage to the Heating Pad Vest.
- Always check the power cord and plug before using the Heating Pad Vest. Do NOT operate the vest if the cord or plug are damaged. A damaged cord or plug must be replaced by a qualified person in order to avoid damages and personal injuries.
- Do NOT use the heating pad vest in or near a bath, shower, swimming pool, or over a basin filled with water when it is turned on. Do NOT operate the device with damp or wet hands and do NOT touch the power cord with damp or wet hands. Make sure to prevent any contact with water. Do NOT submerge the heating pad vest in water. Make sure the controller does not come into contact with water. If water does get inside the casing, turn the device off immediately and contact our customer service.
- If the power fails during usage, turn off the vest immediately and pull the plug from the socket.
- The heating pad vest heats up during use. Do NOT use the device if you are sensitive to heat.
- Keep the heating pad vest away from moisture, heat and in or close to areas where the surrounding temperature is relatively hot (e.g. stoves or other sources of heat).
- Do NOT cover the Heating Pad Vest or the controller when it is switched on. Never put blankets or cushions on top of or against the controller. This poses a risk of fire, electric shock or injury.
- Do NOT disconnect the vest from its power supply by pulling the cable.
- Do NOT carry or pull the device by its cable.
- Make sure that you do NOT trip over the cable. The cable should not be kinked, trapped or twisted.

Other

- Do NOT place objects or heat sources (such as hot water bottles) against or on the vest when it is turned on or stored away.
- The Heating Pad Vest can only be operated using the controller.
- Do NOT use the vest when it is folded, creased or wrapped in a cover.
- Neither pierce the vest with, nor attach safety pins or other spiky or sharp objects to the pad.
- Let the vest cool down completely before folding and storing it. Do not fold it too tight.
- The vest is machine washable. Use a gentle program (30°C) and let it air dry. Make sure to ALWAYS disconnect the controller before washing the vest.

Receiving Your Shipment

- When receiving your shipment, make sure that the contents are consistent with the packing list. Notify your distributor of any missing items. If the equipment appears to be damaged, file a claim immediately with the carrier and notify your distributor at once, giving a detailed description of any damage. Save the damaged packing container to substantiate your claim.

Contents of Package

- Auronic Heating Pad Vest - Model AU3133 - EAN 8720195255529

Included Accessories

- 1x User Manual

Product Features

Technical Data


Size	60 x 90 cm
Weight	950 gr
Colour	Grey
Material	Flannel
Cable Length	2,35 cm

Nominal Voltage	220-240V~
Nominal Frequency	50Hz
Nominal Power Consumption	100W
Temperature Range	53°C (+/- 5°C) - 59°C (+/- 5°C) - 65°C (+/-5°C)

Other Features

- Turns off automatically after 90 minutes of continuous use;
- LED Display;
- Machine washable (30°C).

Product Overview

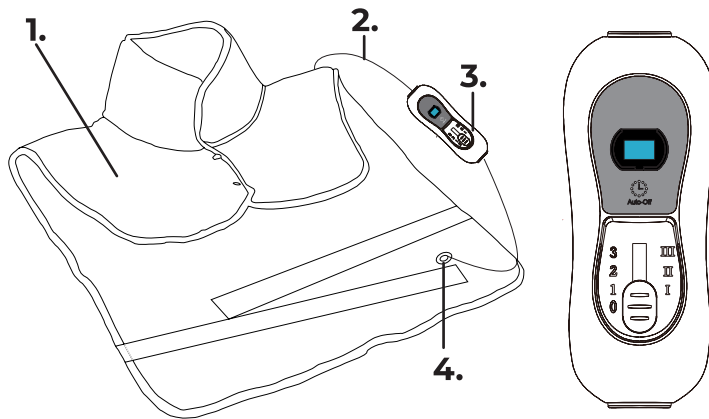1. Heating Pad
2. Power Cord
3. Controller
4. Connector

Temperature Setting

0. - Off
1. - Low Heat - 53°C (+/- 5°C)
2. - Medium Heat - 59°C(+/-5°C)
3. - High Heat - 65°C (+/-5°C)

Operating Instructions

⚠ The Heating Pad Vest should only be used on the neck, shoulders and back. Do NOT place the Vest near our eyes or other sensitive areas. Never place the Vest on or over your head. Do NOT sit on the Vest.

⚠ Prolonged use of the Heating Pad Vest at a high setting could lead to skin burns.

⚠ This vest must not be used by people who are not sensitive to heat or by people who may not be able to react to overheating (diabetics, people with skin conditions or scars, people who are sleeping or unconscious, after drinking alcohol/ using drugs or intoxicatives, or after taking medicine).

Before Use

1. Before use, always inspect the product to ensure there is no damage to the power cord, the controller or the vest itself. The pad connector and control connector are designed to connect only one way.
2. Make sure the pad connector and controller are correctly and securely connected.
3. If there is no power, make sure the AC electrical outlet is functional and the connection between the pad connector and controller is secure. Discontinue use if the power failure continues.

How to Use

1. Connect the plug to a power outlet and slide the slider from 0 to 1. (Low Heat - 53°C +/- 5°C). The Heating Pad will start to heat up. We recommend starting with a lower setting first, then slowly increasing the setting to prevent injuries resulting from the temperature being too high.
2. Slide the slider to setting 2 (Medium Heat - 59°C +/- 5°C) or 3 (High Heat - 65°C +/-5°C) if you prefer a higher temperature. Slide the slider back to a lower setting if the pad becomes too hot. The real-time temperature will be shown on the LED Display.
4. Slide the slider back to position 1 when you are using the pad for an extended period of time.
5. Put the slider back on 0 to turn off the pad. The heating pad turns off automatically after 90 minutes of continuous use.
6. If you are not using the Heating Pad Vest, make sure to turn it off. Slide the slider to 0 and unplug the plug from the power socket.

Maintenance and Cleaning

1. Disconnect the power cord from the electrical outlet and let it cool down for at least 10 minutes.
2. Make sure the pad connector is disconnected from the controller by pressing the tabs on the sides of the connector while pulling the connectors apart.
3. Machine or hand wash the Vest using a mild detergent, cool water and a gentle/delicate wash cycle.
4. Line dry or lay the Vest flat to dry. Do not tumble dry.
5. Make sure the product is completely dry before the next use.

Storage

- Let the Heating Pad Vest cool down completely before folding and storing it. Store the Vest in its original packaging in a clean and dry place. Make sure that you do not damage the power cord. Do not place anything on top of the Vest. Do not store the device in extreme high or extreme low temperatures.

Disposal - Recycle

- Contribute to a cleaner living environment! This product is subject to the European Directive (EU) 2012/19/EU. Do not dispose of the Heating Pad Vest with the household waste. For proper recycling, please contact your local authority or your household waste disposal service for further details on your nearest designated collection point.

Service and Warranty

- If you need service or information regarding your product, please contact the Auronic Customer Support at service@auronic.nl
- Auronic provides a 2 year warranty on its products. To obtain service during the warranty period, the product needs to be returned with a proof of purchase. Product defects have to be reported within 2 years from the date of purchase. This warranty does not cover: damages caused by misuse or incorrect repair; parts subject to wear; defects that the customer was aware of at the time of purchase; damages or defects caused by customer neglect; damages or defects caused by third parties.

Disclaimer

- Changes reserved; specifications are subject to change without stating grounds for doing so.

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© **LifeGoods B.V.**