GESKE

GERMAN BEAUTY TECH

Quickstart Guide

www.geske.com

GESKE Beauty Tech GmbH,
Leipziger Platz 18, 10117 Berlin, Germany
info@geske.com, www.geske.com

SmartAppGuided™
**Skin Scrubber &
Blackhead Remover | 9 in 1**

Short Introduction

Scan for the full instructions
in your language


GESKE German Beauty Tech App
Available in 45 Languages

Apple and the Apple logo are trademarks of Apple Inc.
Google Play and the Google Play logo are trademarks of Google LLC.

1	Exfoliate
2	Moisturize
3	MicroCurrent
4	Warming
5	Warming + Exfoliate
6	Warming + Moisturize
7	Warming + MicroCurrent

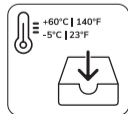Read the safety information before first use.

Keep both instructions for use.

YOUR PRODUCT AT A GLANCE

Your product is a multifunctional MicroCurrent skin scrubber. It is safe to be used up to 1 time per week, at maximum 2 times per week but not more than once per day. Please do not use more frequently. For maximal safety and ideal results, always follow your personalized Skin Care Routine created by the GESKE App.

Download the GESKE App and get your personalized routine.

Intended Use: The Geske & Hello Kitty MicroCurrent Skin Scrubber & Blackhead Remover Devices are over-the-counter home use devices, which are intended for:

- Ultrasonic facial cleansing of blackhead, impurities, & dead skin
- Facial skin stimulation with MicroCurrent technology
- Facilitating the absorption of skincare products

SPECIFICATIONS:

Materials: silicone, ABS, PC, stainless steel | **Storage temperature:** 5 - +60 °C/23 - 140 °F | **Battery:** Li-Ion 3.7V | **Battery Capacity:** 700 mAh | **Maximum Charging Current:** 450 mA | **Charging Voltage:** 5V | **Charging Time:** approx. 3.5 hours | **Operating Time:** Operating time differs depending on activated technologies & intensities. | **Interface:** USB | **Programs:** Vibration: 1 intensity, ION Technology, MicroCurrent: 1 intensity, LED Light Technology: 1 colour

IMPORTANT SAFETY NOTES! READ BEFORE USAGE!

The product is not suitable for use under the following conditions:

- If you have an implanted medical device or other electronic device such as a pacemaker
 - After laser treatment or chemical peeling
 - After cosmetic facial surgery
- In the following cases, a consultation with a physician is necessary as to whether the product is suitable for you to use:**
- If you have a heart disease
 - In case of pregnancy. The safety of use during pregnancy has not been established
 - The product can emit intense light, avoid applying it in the direction of or near an open eye

Further information on use:

- Before using the MicroCurrent Face-Lift Technology, the product must be completely dry

- Use of the MicroCurrent Face-Lift Technology only in combination with the GESKE serum
- Do not use the microcurrent function on the chest area, around the eyes, midline of the neck, genital/grain area, or on red and inflamed skin, over birthmarks and enlarged blood vessels
- Nothing is known about the long-term effects of microcurrent
- As a result of stimulation of the optic nerve, a flashing light may be perceived during use. If this persists when not using the product, seek medical advice
- A slight tingling sensation may occur during application. Decrease the intensity of the Microcurrent to reduce or eliminate tingling
- Electronic monitoring devices may not function properly while the device is in use
- Never use the product while driving or operating heavy machinery
- Do not poke any sharp objects into the product's openings
- The product may only be used with a safety extra-low voltage (SELV) adapter
- We recommend using standard IEC60950 power supplies to charge the device
- Never use the skin scrubbing function in the same place for a longer period of time to avoid skin irritation, swelling or slight bruising
- Do not use the product on irritated or injured skin or open wounds (after any kind of injury), in the case of bacterial or viral infections of the skin (either in the usage area and elsewhere)
- This product is not intended to treat medical skin conditions or diseases. If you have medical skin conditions or diseases, consult your dermatologist before using this product
- Do not use the product if it is 1) untypically warm, 2) mechanically damaged or deformed, 3) discolored
- If pain and/or discomfort occurs during use, discontinue use and seek medical advice
- The eye area is very sensitive: Be extra careful in this area and avoid any contact of the product with the eyelids or the eye itself
- Only use the product as recommended in your personalized routine
- For hygiene reasons, we recommend not sharing your product with others and replacing attachments as needed
- The product is intended for cosmetic use on the face. GESKE assumes no liability for harmful consequences arising from the use on other body areas or other improper use, connection to unsuitable voltage sources, contaminated conducting liquid or balls
- Always make sure the plug and socket are completely dry before charging. Otherwise, electric shocks, a short circuit or fire may result. The product must not be used while connected to the power outlet

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