USER MANUAL

English & Dutch

TWINFLOW DUO

HANDSFREE BREAST PUMP

Table of contents

English manual

Product Description	2
Product Components & Layout	3
Control Panel & LED Indicators	5
Installation	6
Starting Your Pumping Session	8
Battery Status Indicator	12
Recommended Pumping Position	13
Ending a pumping session	14
Cleaning and maintenance	15
Storing breast milk	17
Warranty	18

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English manual

Product Description

The Sensibaby Twinflow Duo wearable breast pump is designed to support mothers during their breastfeeding journey by providing comfortable and effective milk expression.

The pump is designed to mimic a natural sucking rhythm, helping to stimulate milk flow and support regular milk expression when direct breastfeeding is not possible or convenient.

Regular use of the breast pump can help promote effective milk removal, supporting breast comfort and helping to relieve the feeling of fullness that may occur during lactation. Consistent pumping may also support maintaining milk supply as part of a breastfeeding routine.

Using a wearable breast pump can offer flexibility and convenience, helping mothers continue breastfeeding with confidence while balancing daily activities. This can contribute positively to the breastfeeding experience for both mother and baby.


The Sensibaby Twinflow Duo is intended for personal use only.

It is not intended for clinical use, medical treatment, or use as a medical device in healthcare institutions.

Product Components & Layout

The product consists of the following main components:

1. Flange
2. Pump motor
3. Silicone diaphragm
4. USB charging cable
5. Milk collector
6. Bra adjustment buckles
7. valve

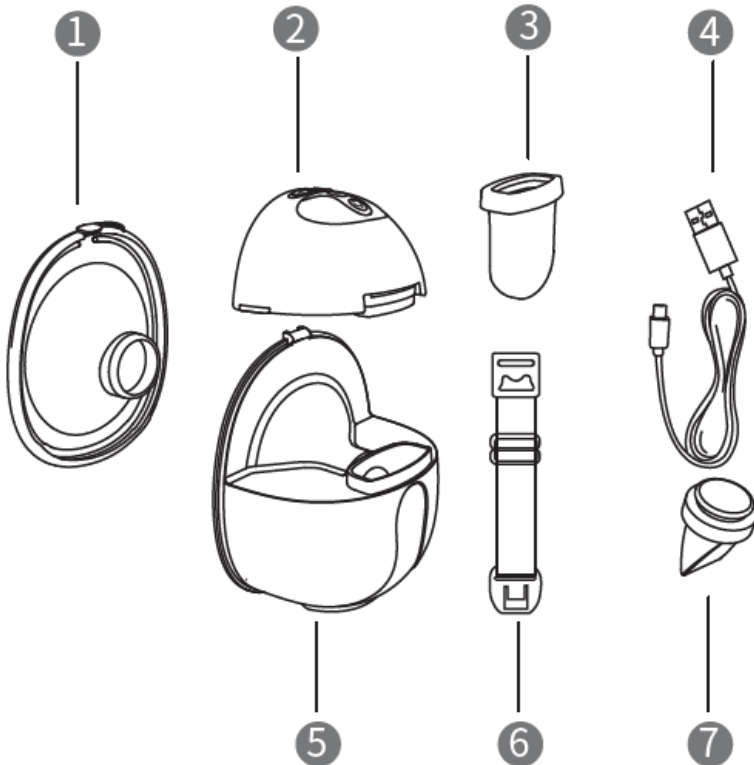Depending on the configuration, the box contains

- 2x Pump motorunits
- 2x Milk collectors
- 2x Silicone breast shields - 24mm
- 2x Duckbill valves
- 2x Diaphragms
- 2x Bra adjustment buckles
- 1x Nipple ruler
- 2x Nipple inserts - 17mm
- 2x Nipple inserts - 19mm
- 2x Nipple inserts - 21mm
- 10x Milk storage bags
- 2x USB charging cables
- 1x Zipper storage bag
- 1x User manual

Functions & Features

Handsfree & wearable design

Dual pumping (left & right)

Multiple suction levels and modes


Quiet operation

Rechargeable battery

Built-in night light function for use in low-light environments

Control Panel & LED Indicators

1. Power display
2. Expression mode
3. Stimulation mode
4. Massage mode
5. Automatic mode
6. Night light power button
7. Level up button
8. Power button/mode switch button
9. Time/minute
10. Level display/ Time display
11. Pause/ Start button
12. Level down button

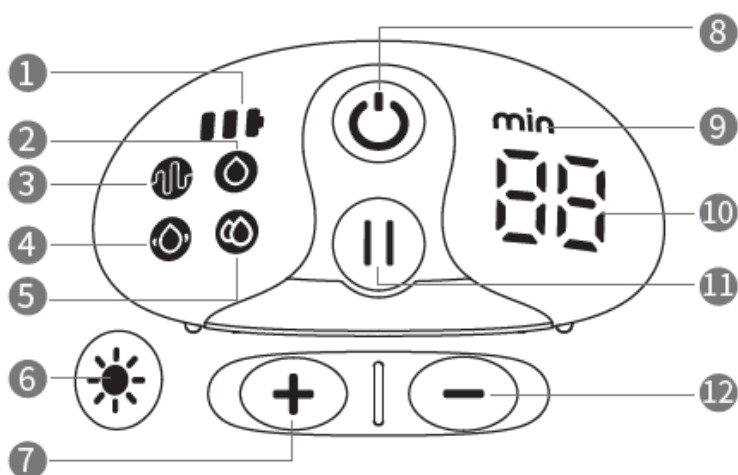Installation


Follow the steps below to assemble the Sensibaby Twinflow Duo correctly.

Important Notes


- Always clean and dry all pump components before assembly.
- Wash your hands thoroughly before handling cleaned parts.
- Make sure all components are correctly positioned and securely sealed to prevent leakage.
- Remove any milk residue from the rims before reassembling the pump.

Assembly Steps

1. Attach the duckbill valve to the bottom of the milk collector.
Ensure it is fitted in the correct direction and fully secured.

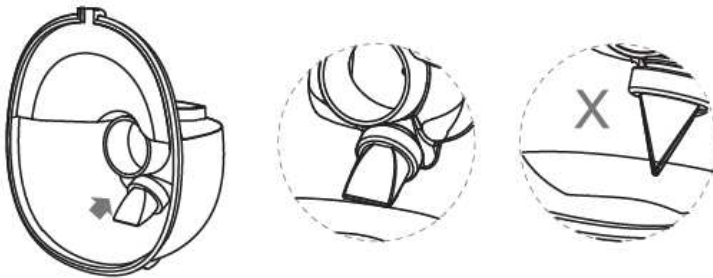
2. Place the silicone diaphragm into the milk collector.
Gently press around the edges with your fingers to ensure it fits tightly and creates an airtight seal.

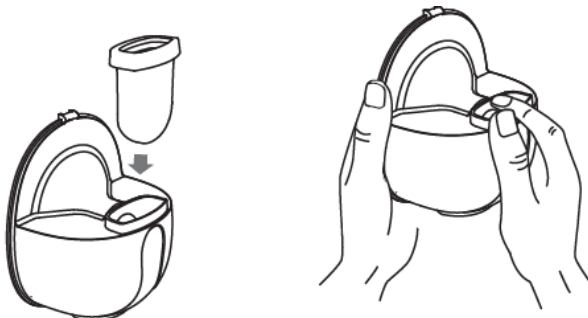
3. Position the breast shield (flange) onto the milk collector.

Align the guide groove on the flange with the corresponding groove on the milk collector, then press evenly around the edges to secure the seal.

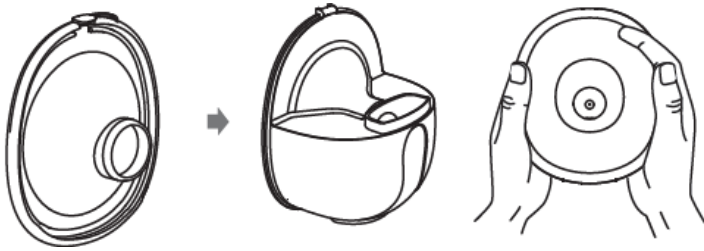
4. Insert the milk outlet plug into the bottom of the flange.

Make sure it is fully aligned and tightly sealed.

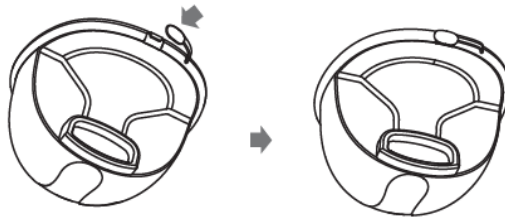5. Place the motor unit onto the top of the milk collector.

Align the notch on the motor unit with the silicone diaphragm and press down gently until it fits securely and forms a complete seal.

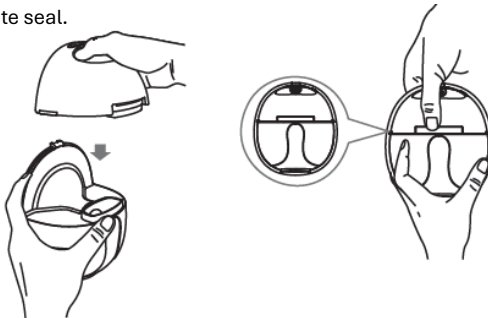Before use, double-check that all parts are properly assembled and firmly connected.

Starting Your Pumping Session

Before You Begin


- Ensure the pump is fully assembled and clean.
- Sit comfortably and keep the pump upright during use.

1. Positioning

Press the breast shield firmly against your breast, ensuring there are no gaps.

Make sure the nipple is centered in the flange opening.

Always keep the pump upright. Incorrect positioning may result in reduced or no suction.

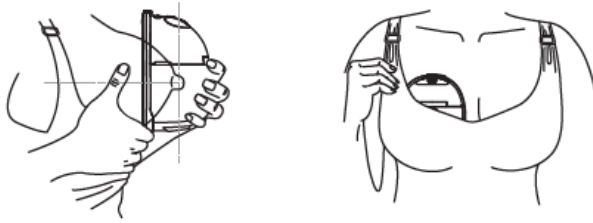2. Power On / Off

Long press the Power button to turn the pump on.

The display will light up.

Long press the Power button again to turn the pump off.


3. Start/ Pause & Mode Selection

- Short press the Start/Pause button  to start or pause pumping.
- While pumping, short press the On/Off button  to switch between modes.

The pump offers four modes:

- Stimulation Mode  – gentle rhythm to stimulate milk flow.
- Expression Mode  – steady suction for milk expression.
- Massage Mode  – gentle pulses for comfort.
- Automatic Mode  – cycles automatically between patterns.

The active mode is indicated on the display.

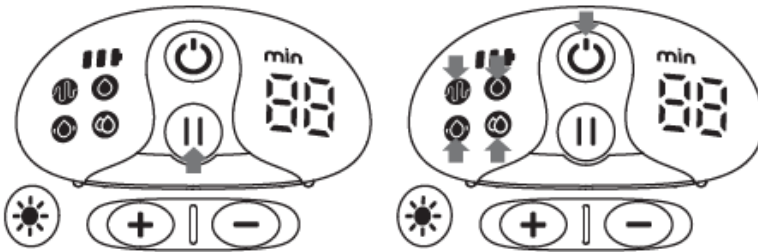
4. Adjusting Suction Levels

Each mode has 9 suction levels.


- Press “+” to increase suction
- Press “-” to decrease suction

Choose the level that feels comfortable for you.

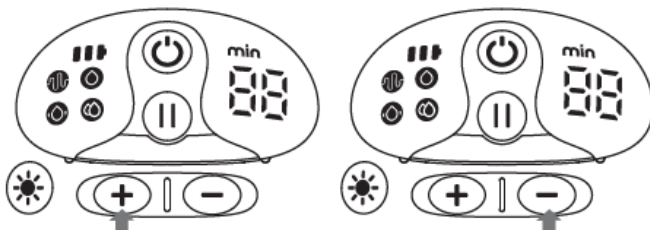5. Display Information

- When a button is pressed, the screen shows the current suction level.
- If no button is pressed for 2 seconds, the display automatically switches to the timer view, showing pumping time in minutes.

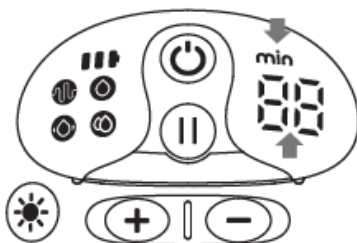6. Night Light Function

- Long press the Light button to turn the night light on or off.
- Short press the Light button to adjust brightness.
- There are three brightness levels available.

The night light operates independently from the pump and can remain on even when the pump is switched off.

Helpful Tips

1. The pump has a memory function. When turned on again, it will start in the stimulation mode and suction level used during the previous session (automatic mode always starts at level 1).
2. In automatic mode, suction levels are not memorized.
If no buttons are pressed, the suction level will automatically increase up to level 7 at 15-second intervals.
Pressing any button will stop the automatic increase.
3. If the pump is paused and no button is pressed for 5 minutes, it will automatically shut down.
4. For safety and comfort, the pump will automatically stop after 20 minutes of continuous operation.

Battery Status Indicator


The breast pump display shows the battery status using three bars. During use, if only one bar is visible and blinking, this indicates that the battery level is low.

Please recharge the pump promptly.

The pump cannot be used while charging.


When connected to a power source, the display will show the battery status only.

When all three bars are visible and no longer flashing, the battery is fully charged.

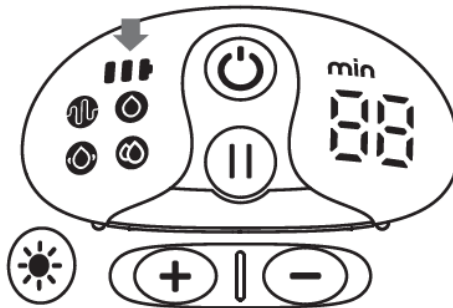After each pumping session, you may briefly press the Power button to check the battery level.

If the battery is low, recharge the pump using a DC 5V / 1A charger.

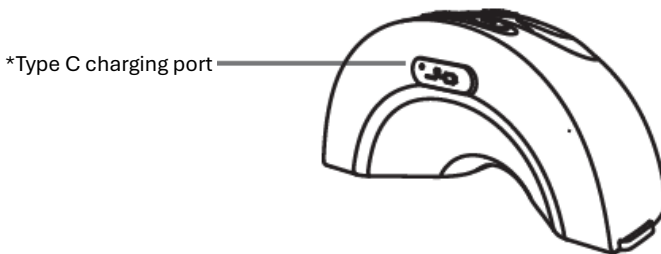Recommended Pumping Position

Sit comfortably in a chair and relax your shoulders.

Lean slightly forward and use a cushion or pillow for back support if needed.

Make sure your body and the breast pump are aligned in a natural, straight position.

This helps ensure proper suction and a more comfortable pumping experience.

Bra Adjustment Instructions

Correct bra adjustment can improve comfort and help keep the pump securely in place.

- Adjust the length of the bra straps so the pump fits snugly without discomfort.
- Attach one end of the bra adjustment buckle to the inner clasp of your nursing bra.
- Adjust the strap on the other end of the buckle to achieve a secure and comfortable fit.

Proper adjustment helps maintain consistent suction during pumping.

Ending a pumping session

Once you have finished pumping, carefully remove the breast pump from your breast.


1. Press the Power button to switch off the pump.

Gently break the seal between your breast and the flange using a finger, then carefully remove the pump while keeping it upright.


2. Detach the motor unit from the milk collector.

Set the motor unit aside and prepare the milk collector for feeding or storage. Clean all used parts according to the cleaning instructions.

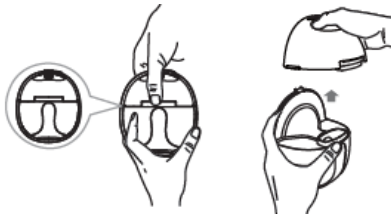3. Pour the expressed breast milk into a storage bottle.

When opening the milk outlet plug, ensure the pouring slot is properly aligned with the bottle opening to prevent spills.

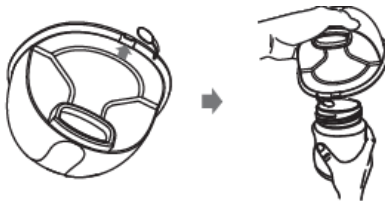Handle expressed milk with care to maintain hygiene and avoid unnecessary loss.

Cleaning and maintenance

Proper cleaning and maintenance are essential for hygiene, safety, and optimal pump performance.

Pump Motor Unit

After each use:

- Wipe the exterior with a clean, damp cloth.
- If needed, disinfect using 70% ethyl alcohol or isopropyl alcohol.
- Do not immerse the motor unit in water or place it under running water.

Flange, Silicone Diaphragm, Valve & Milk Collector

Before First Use

- Boil all washable parts in water for 5 minutes.
- Allow parts to cool completely and air-dry before assembly.
- Do not microwave any pump components.

After Each Use

- Clean all parts using warm water and mild dishwashing soap or rinse thoroughly with clean water.
- Rinse well to remove any soap residue.
- Allow all parts to air-dry completely before storage or reuse.

Valve Cleaning

- Clean the valve gently in warm water with a few drops of mild cleaning liquid.
- Do not insert objects into the valve, as this may cause damage.
- Handle the valve with care during cleaning and disassembly.
- If the valve is damaged, suction performance may be affected.

Recommended Hand-Washing Method

1. Place pump parts in a clean wash basin used only for infant feeding items.
Do not place parts directly in the sink.
2. Add hot water and mild dish soap to the basin.
3. Clean parts using a brush dedicated solely to infant feeding items.
4. Rinse items thoroughly under running water or by submerging them in fresh water in a separate basin.
5. Allow all parts to air-dry completely on a clean, unused dish towel or paper towel in an area protected from dirt and dust.
Do not rub or pat items dry with a towel.
6. Clean the wash basin and brush after each use.
Rinse well and allow them to air-dry. Wash them by hand at least every few days.

Storing breast milk

1. Only store breast milk that has been collected using a clean and sterilized breast pump.
2. Expressed milk should be transferred into a clean, sterilized milk storage bag or bottle with a secure seal.

Place the milk in the refrigerator or freezer as soon as possible after pumping.

Important Storage Guidelines

1. Do not refreeze breast milk once it has been thawed.
2. Do not mix freshly expressed breast milk with already frozen breast milk.
3. Discard any breast milk that remains after a feeding.
4. Expressed breast milk may be:
 - Refrigerated for up to 48 hours (inside the refrigerator, not in the door), or
 - Frozen for up to 3 months.

Clearly label each bottle or storage bag with the date and time of expression.

Thawing Breast Milk

- For best results, thaw frozen breast milk overnight in the refrigerator.
- If needed more quickly, place the sealed bottle or bag in a bowl of warm water.
- Once thawed, breast milk should be used within 24 hours.
- Do not refreeze thawed breast milk.
-

Warming Breast Milk

1. Stand the baby bottle upright in a bowl of warm water or use a baby bottle or food warmer designed for this purpose.
2. Always check the temperature of the milk before feeding your baby.
3. Do not microwave breast milk, as this may destroy valuable nutrients and antibodies and can cause uneven heating.
4. Uneven heating may result in hot spots that could burn your baby.
5. Do not place milk storage bottles in boiling water, as this may cause the bottle to crack or rupture.

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Sensitest B.V.
Sensibaby (brand of Sensitest B.V.)