Instructions for use
Inserts

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
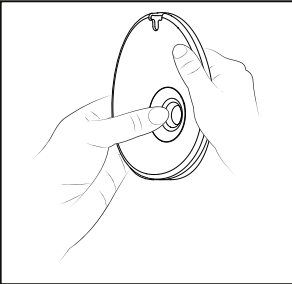
pic. 2

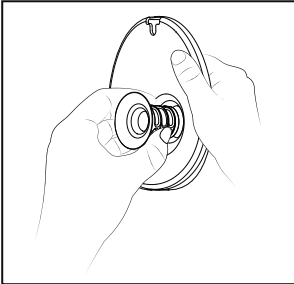pic. 3

Product Description:
The soft silicone inserts are designed to fit within the Medela 24 mm breast shields, enabling precise size adjustment. Compatible with the following Medela breast shields: Motion InBra™, Hands-free, PersonalFit Flex™, PersonalFit™ PLUS, for use with the following breast pumps: Harmony™, Motion InBra™, Swing Maxi™, Solo™, Freestyle™, Pump In Style.
Intended Use:
The Inserts are intended to be used with compatible Medela breast pump systems to reduce the size of 24 mm breast shields for improved fit and comfort.
Sizes Included: 2 × 21 mm, 2 × 18 mm, 2 × 15 mm
Instructions for Use:
Choose the appropriate insert size based on your nipple diameter (see pic. 1). Clean the inserts prior to first use. Place the insert securely and tightly inside the breast shield (see pic. 2). Use the assembled breast shield system with a compatible Medela breast pump as described in the pump's Instructions for Use. Disassemble the insert from the shield before cleaning (see pic. 3)
Cleaning Instructions:
Wash and disinfect the inserts using the same method recommended for your breast shields. For detailed cleaning procedure, refer to the Instructions for Use provided with your compatible Medela breast pump.
Warnings:

- Do not use if the insert is damaged or if you are allergic to silicone.
- Discontinue use if irritation occurs.
- Not for use with non-Medela breast shields / non-compatible pumps.

Additional Info:
Replace after 3 months of use. Refer to the compatible Medela breast pump IFU for additional information.
Material: silicone.
Disposal: via local recycling systems for silicone parts, where available.

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