



URS-2K

Reagent Strips for Urinalysis



HOW TO USE:

1. Immerse the strip into urine for 2 second & take out immediately.
2. Wipe off excess urine against the rim of the cup and hold horizontally for 30 seconds.
3. Match the strip color to the color chart on the bottle label. Make sure to read the correct parameter at the correct time.

IMPORTANT

- Keep away from light and moisture.
- Promptly replace cap after taking out strips.


Do not use in the presence of disinfectants.
Do not touch test areas of reagent strips.
Use it within the expiration date.
Read the label carefully before use.

**Bronstore,
Westmaaspad 101,
2729KH,
Zoetermeer
Nederland**

LOT: 2025102833
EXP: 10-27-2027



IVD

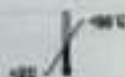

8 721042 450487

Twentyseven® Ketone test (100 strips) - URS-2K LVB