
ROLLATOR | FORTIS

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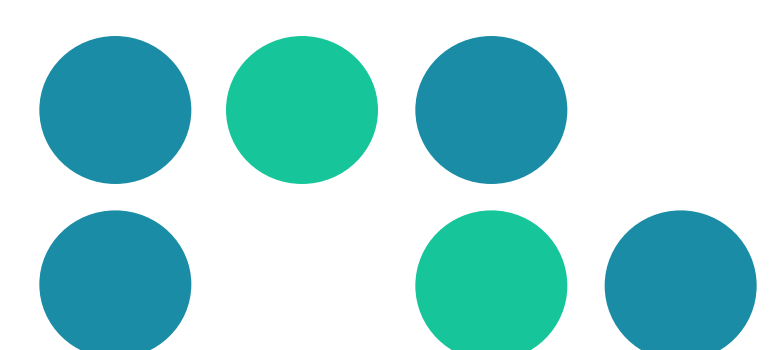Over de FORTIS-rollator

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
ROLLATOR | FORTIS USER MANUAL

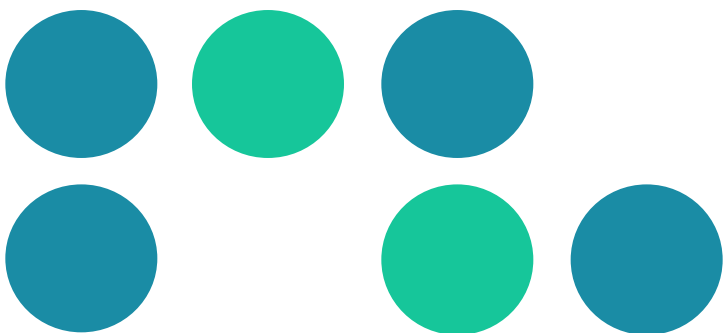About the FORTIS Rollator

The FORTIS rollator is lightweight and made of sturdy aluminum. It features a seat, a shopping bag, and reliable parking brakes. The shock absorption ensures a comfortable ride. The frame is easy to fold for transport and storage. Ideal for daily use by elderly users with reduced strength or mobility.

Use

The rollator is a walking aid with a seat and shopping bag, suitable for people with limited mobility. It is designed for both indoor and outdoor use, up to a maximum user weight (see specifications).

The user must have sufficient strength, vision, and cognitive ability to use the rollator safely. If in doubt, consult your doctor or dealer.

Always read the safety instructions carefully to prevent injury and damage.

Safety Information

The rollator supports a safe and natural walking pattern. You walk with your weight positioned between the hand grips and the rollator in front of you.

Read this manual carefully before use.


Tested according to EN ISO 11199-2:2021.

The Unique Device Identifier (UDI) is stated on the product label and/or packaging.

Safe Use of the Rollator

1. Always engage the parking brakes when the rollator is stationary.
2. Check before use that all wheels rotate freely.
3. Ensure the brakes are functioning properly.
4. Check that the height adjustment is securely locked.
5. Always fold out the rollator completely before use.
6. Pay close attention to your surroundings while walking.
7. Do not place hands in the folding mechanism; this may cause injury.
8. Use the rollator only as intended; it is not a toy.
9. Never sit down without first applying the parking brakes. This prevents rolling away.
10. Do not let yourself fall onto the seat. Sit down calmly and in a controlled manner. Ask for help if needed.
11. Do not hang heavy objects on the hand grips. This affects braking and stability.

By following these safety rules, you will use the rollator safely and reliably.

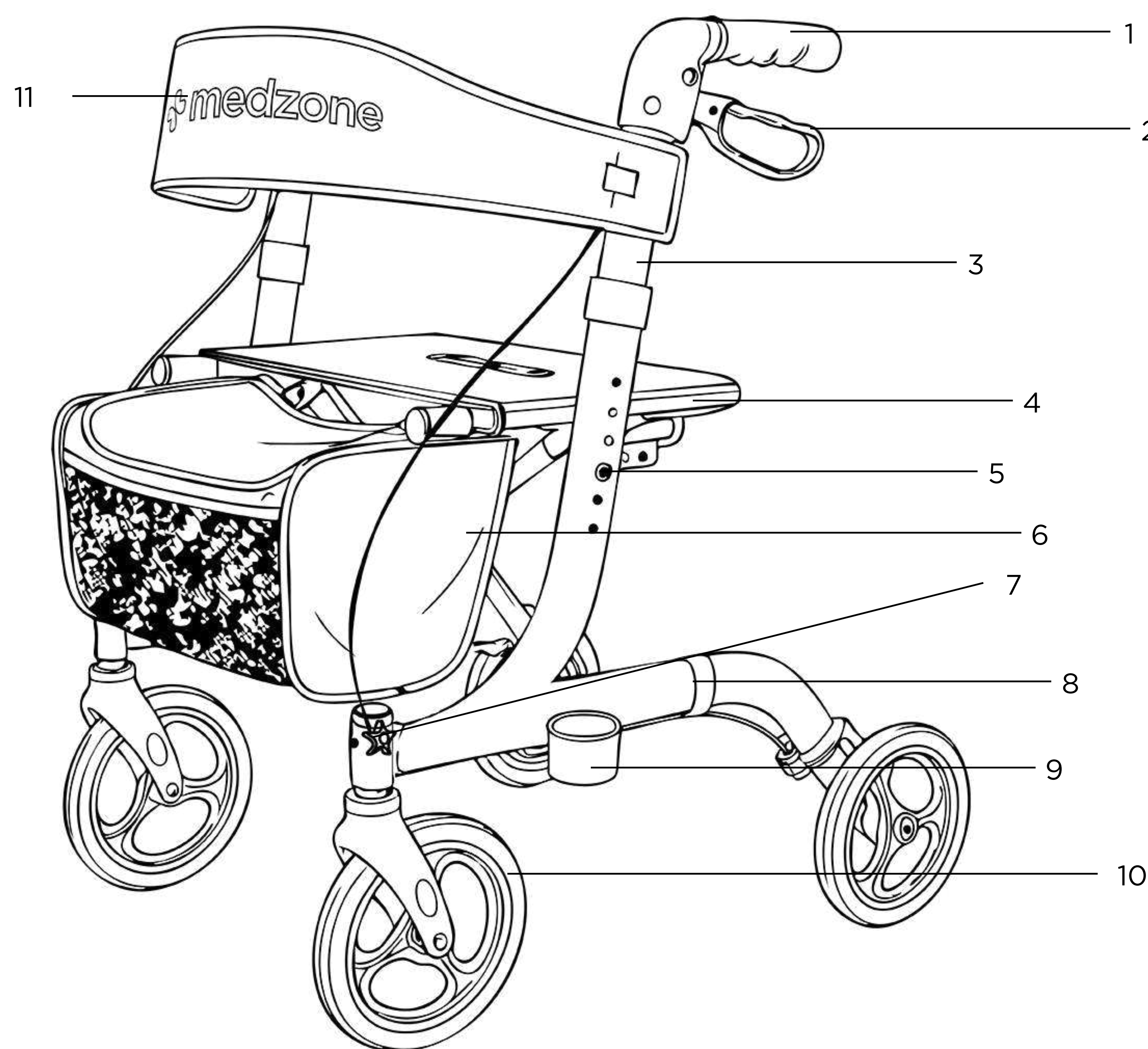Rollator Components

- 1. Hand grip
- 2. Brake lever
- 3. Adjustable backrest
- 4. Seat
- 5. Push button for height adjustment
- 6. Storage bag
- 7. Mounting knob for front wheel installation
- 8. Mounting knob for rear wheel installation
- 9. Push handle
- 10. Wheels (diameter: 20 cm)
- 11. Frame

Disclaimer:

Your rollator consists of several components. Despite the care taken in the design and development of the product, we reserve the right to modify the design and specifications.


Technische specificaties

Width between hand grips		53 cm
Total width		61.5 cm
Total depth		70 cm
Seat width		46 x 25 cm
Seat height		55 cm
Adjustable handle height		82 - 95 cm
Maximum user weight*		136 KG
Product weight		7 KG


**The aluminum frame is designed to support the maximum user weight of 136 kg.*

Assembly Instructions


1. Remove all packaging materials.

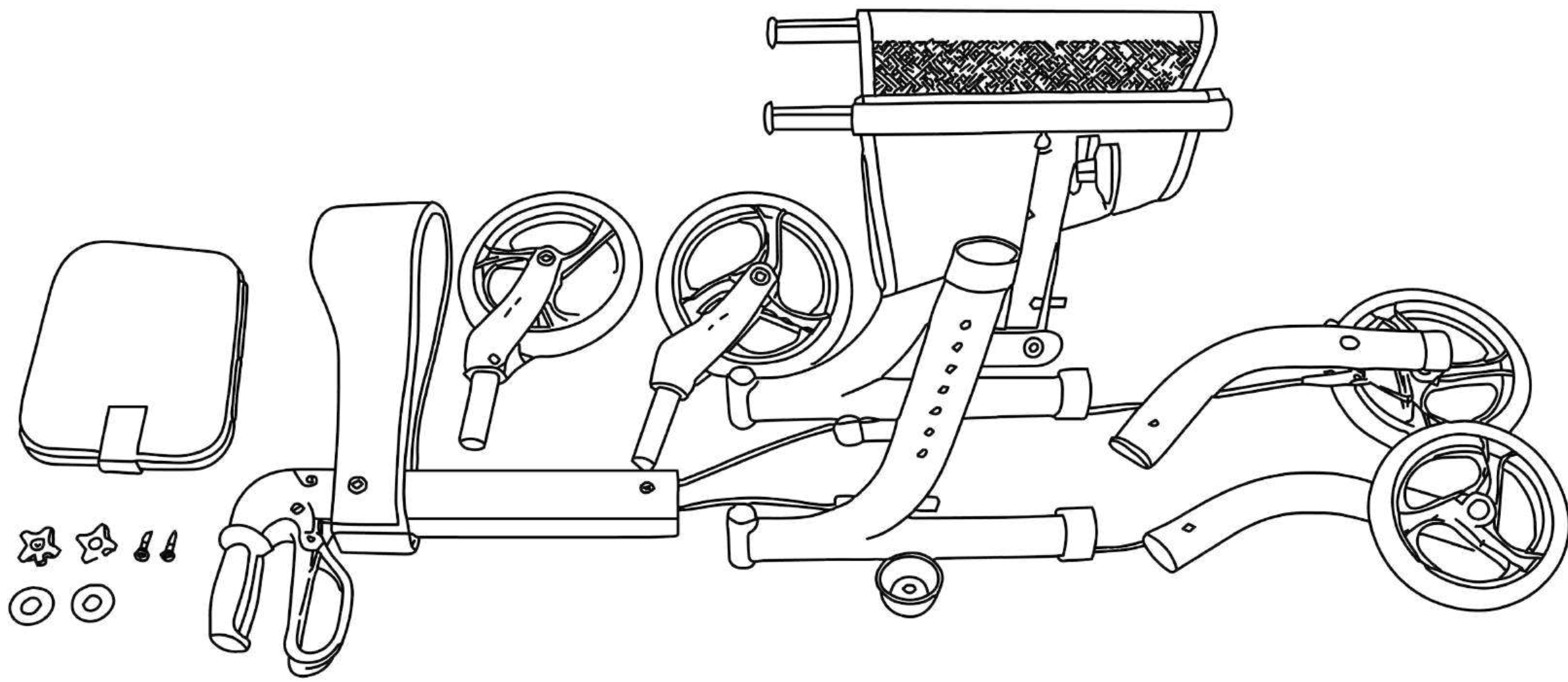
2. Place the frame upright and push both sides downward until you hear a click.

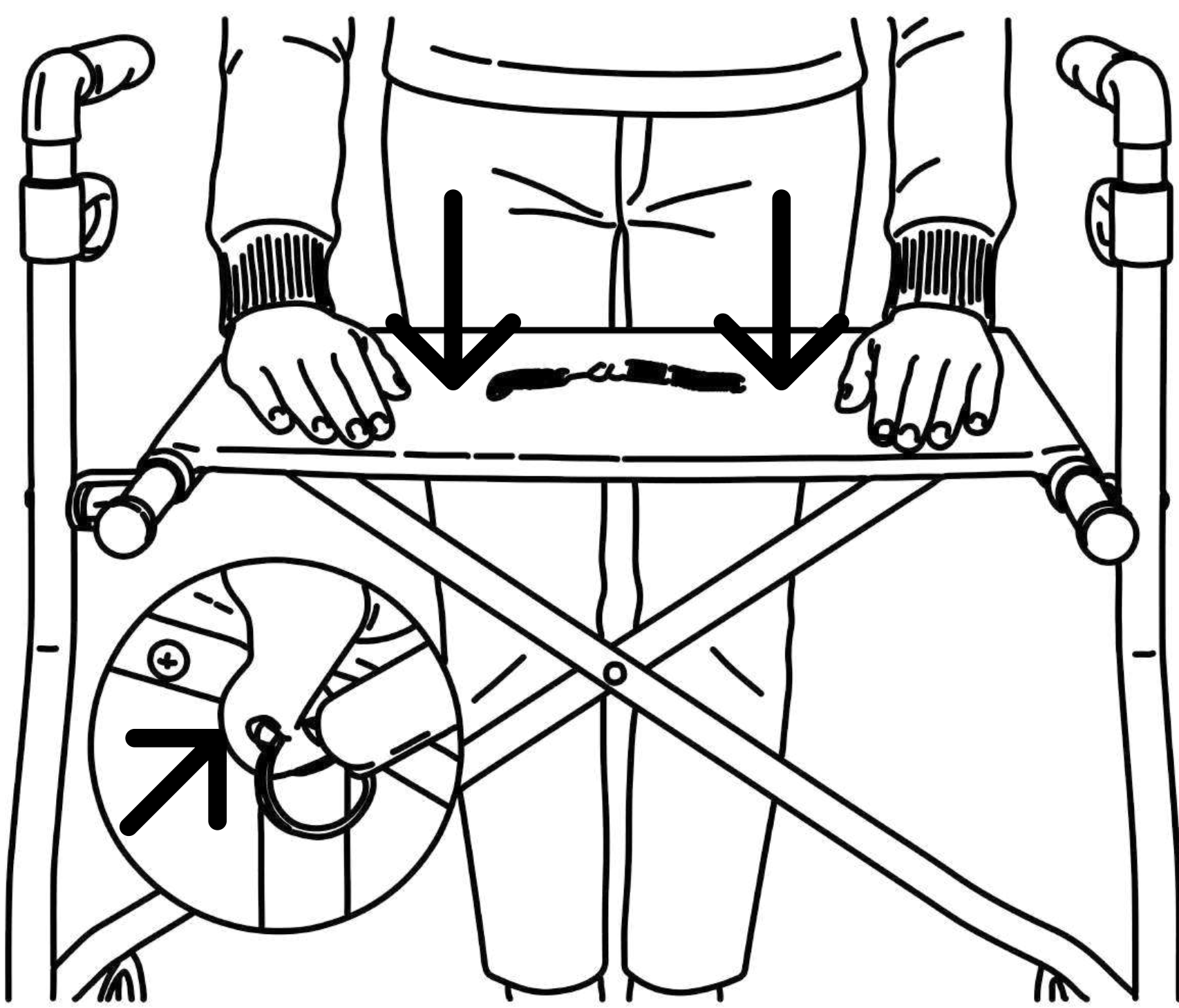3. Install the front wheels. Insert the bolt through the hexagonal hole, add the washer, and tighten the knob securely.

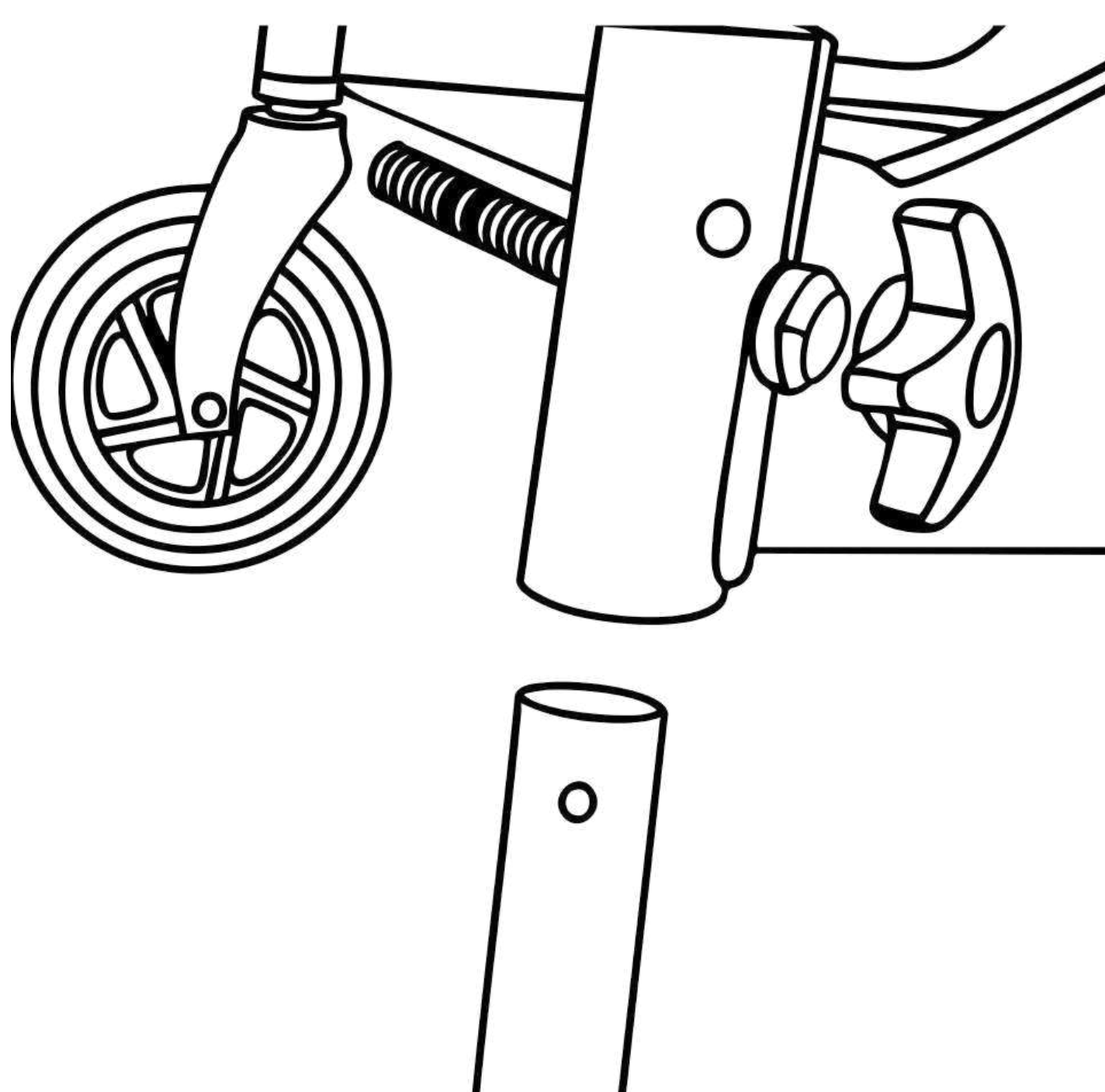
4. Attach the rear wheels. Slide the rear tube into position and press until the locking button clicks into place.

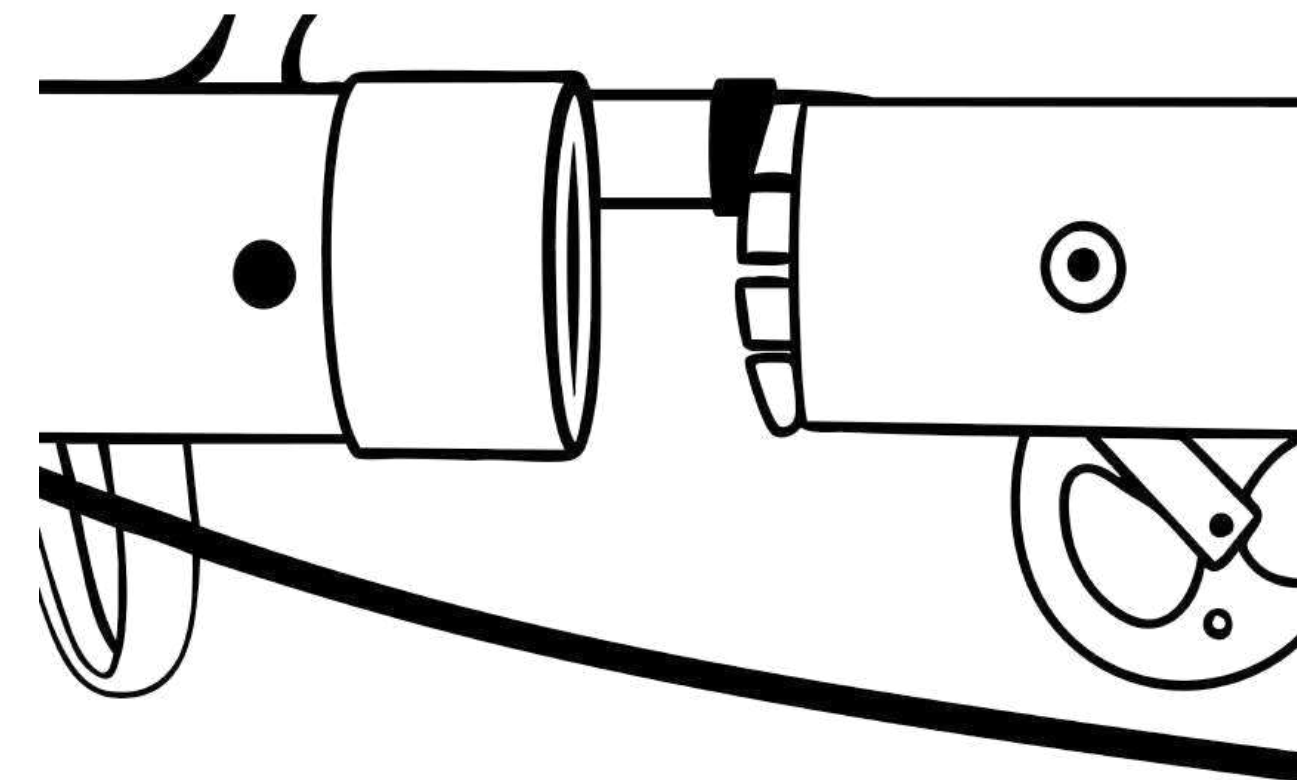Adjusting the Handle Height


1. Press the button (5) on the frame inward.
2. Slide the handle grips up or down to the desired height.
3. Release the button so it locks back into place.
4. Repeat this process on the other side.

Note: Ensure both handle grips are set to the same height and are securely locked before use.

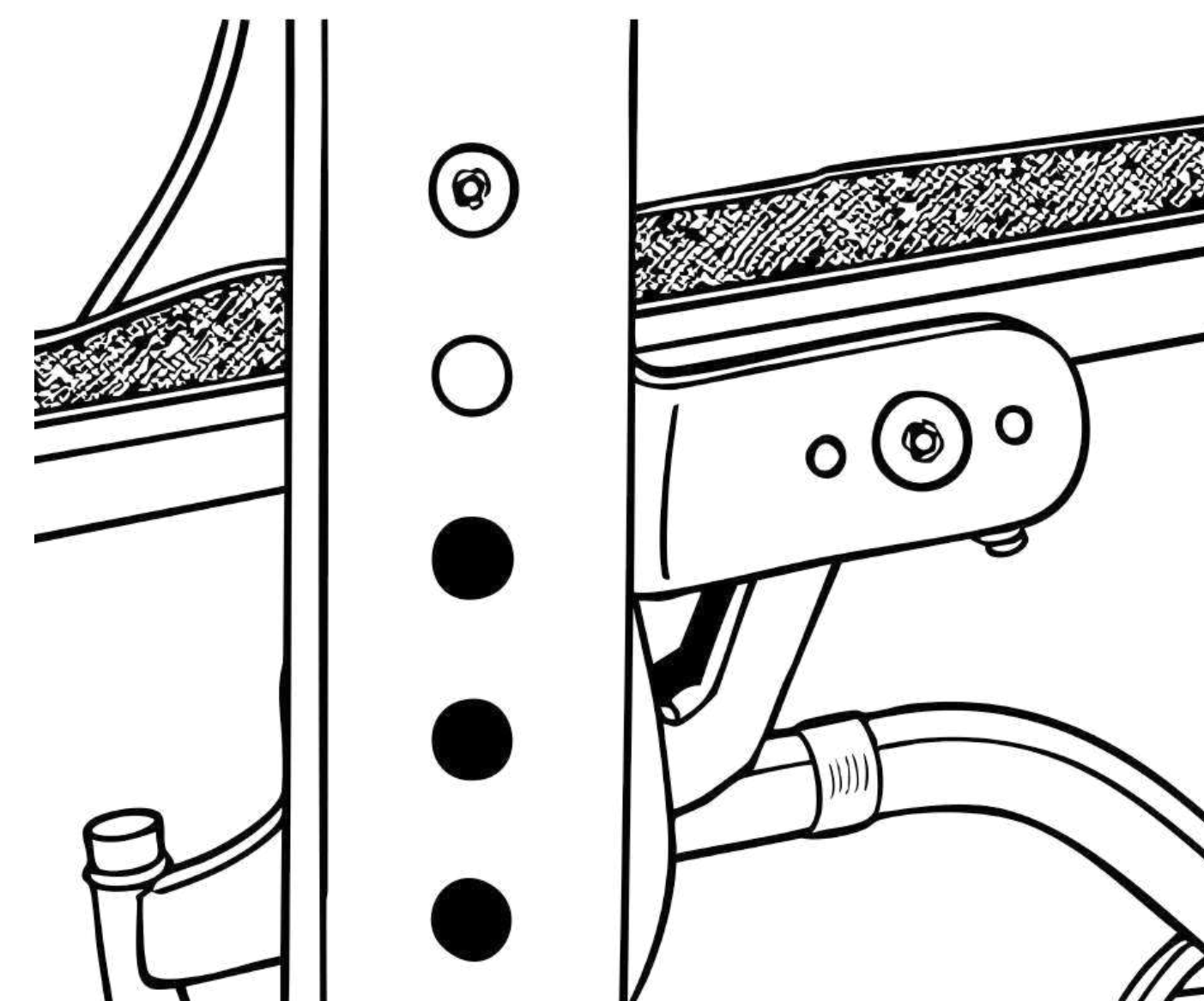
Adjusting the Brakes


- Turn the adjustment screw clockwise to increase brake tension.
- Hold the adjustment screw and tighten the lock nut to secure the setting.

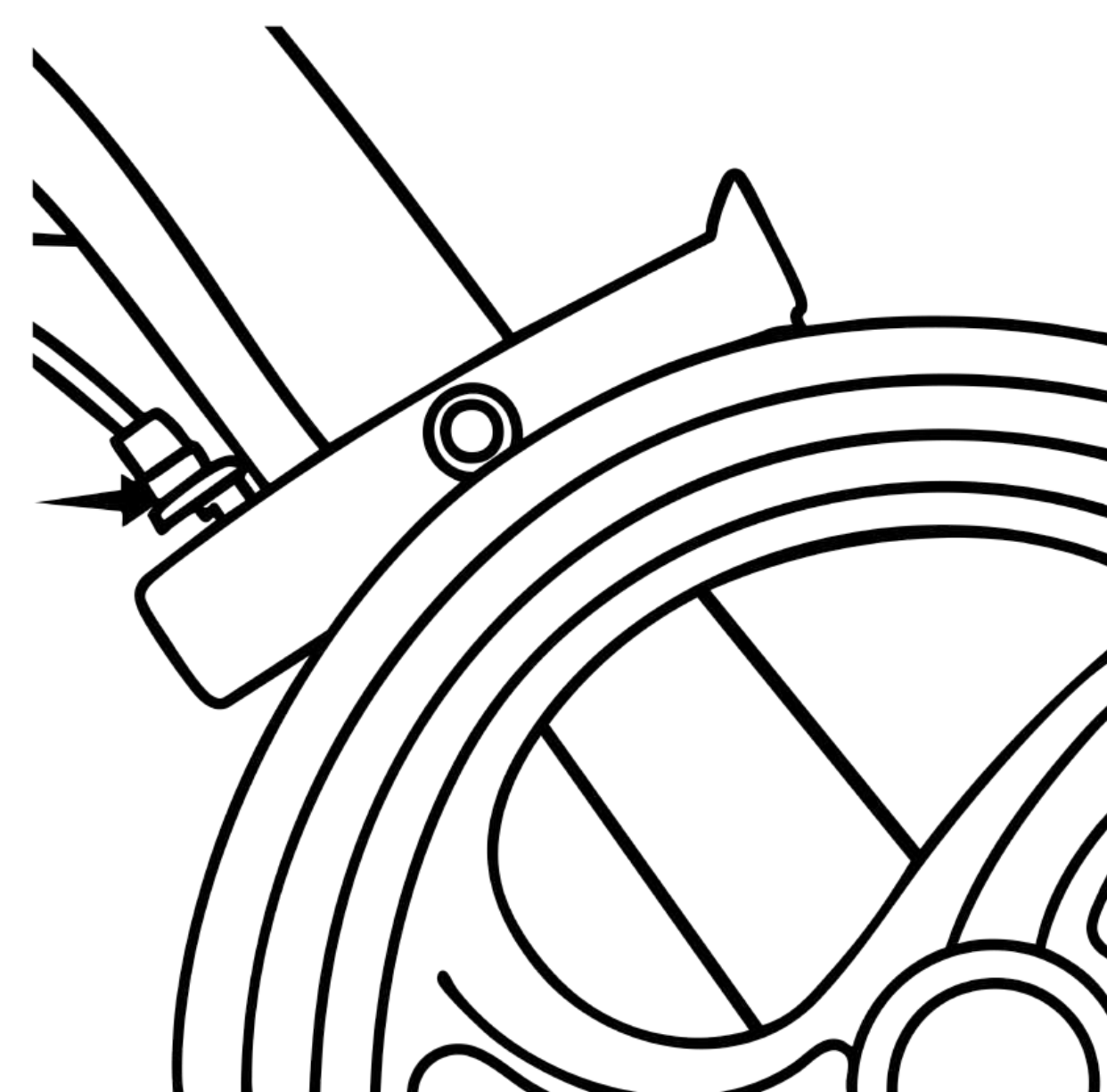
Rollator folded

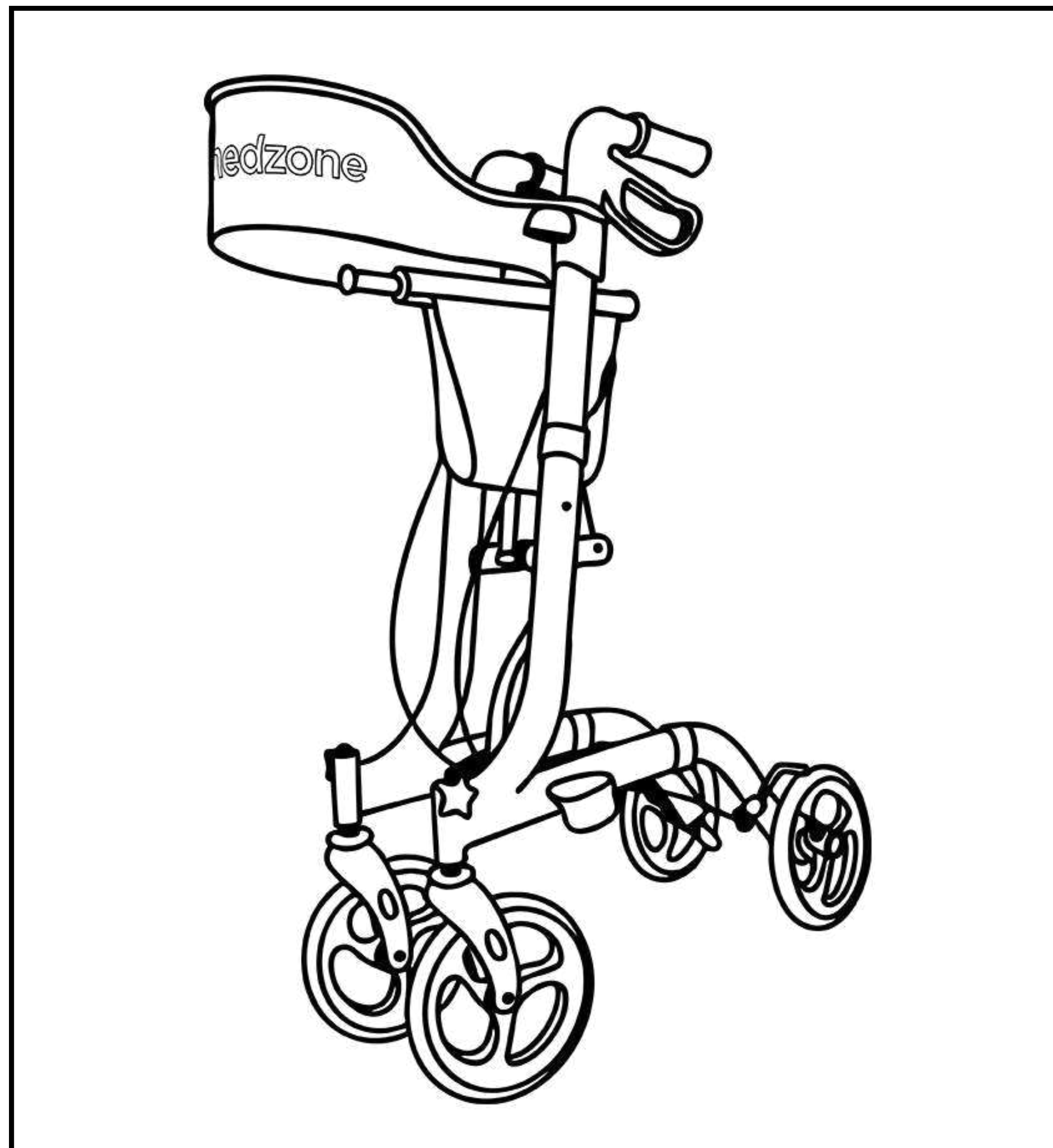
Attaching the shopping bag

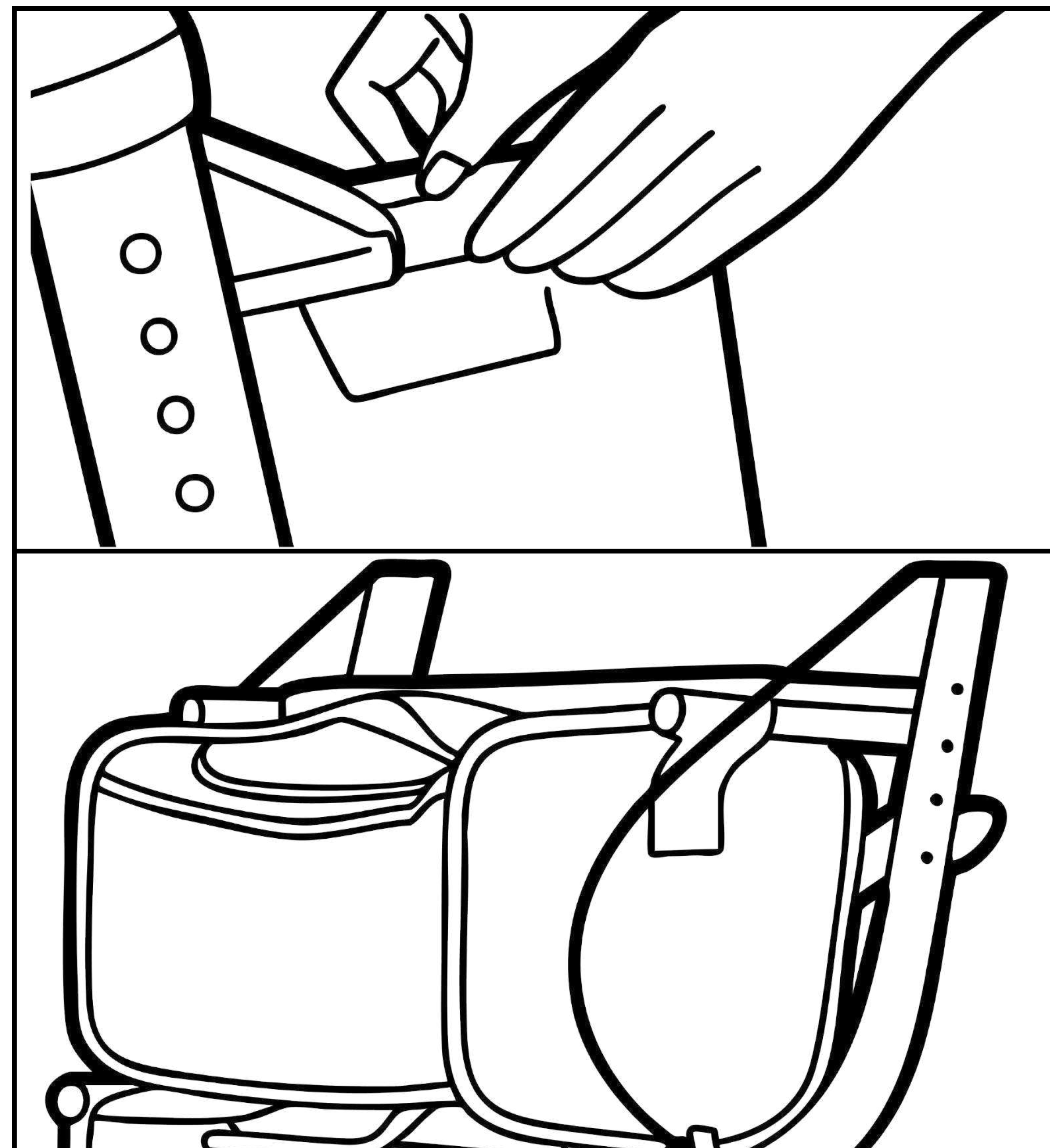Parking brake

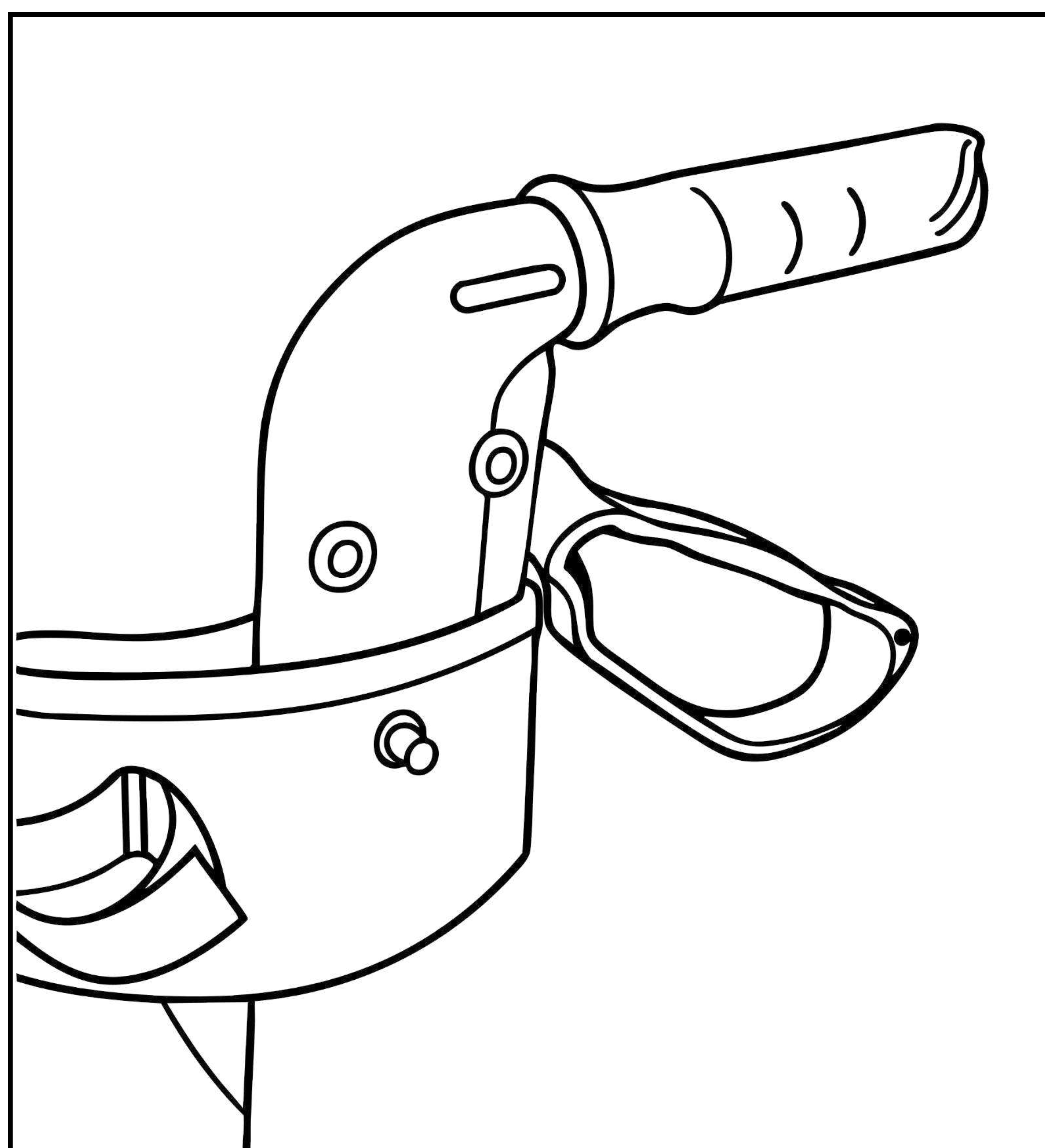Front wheel correctly installed

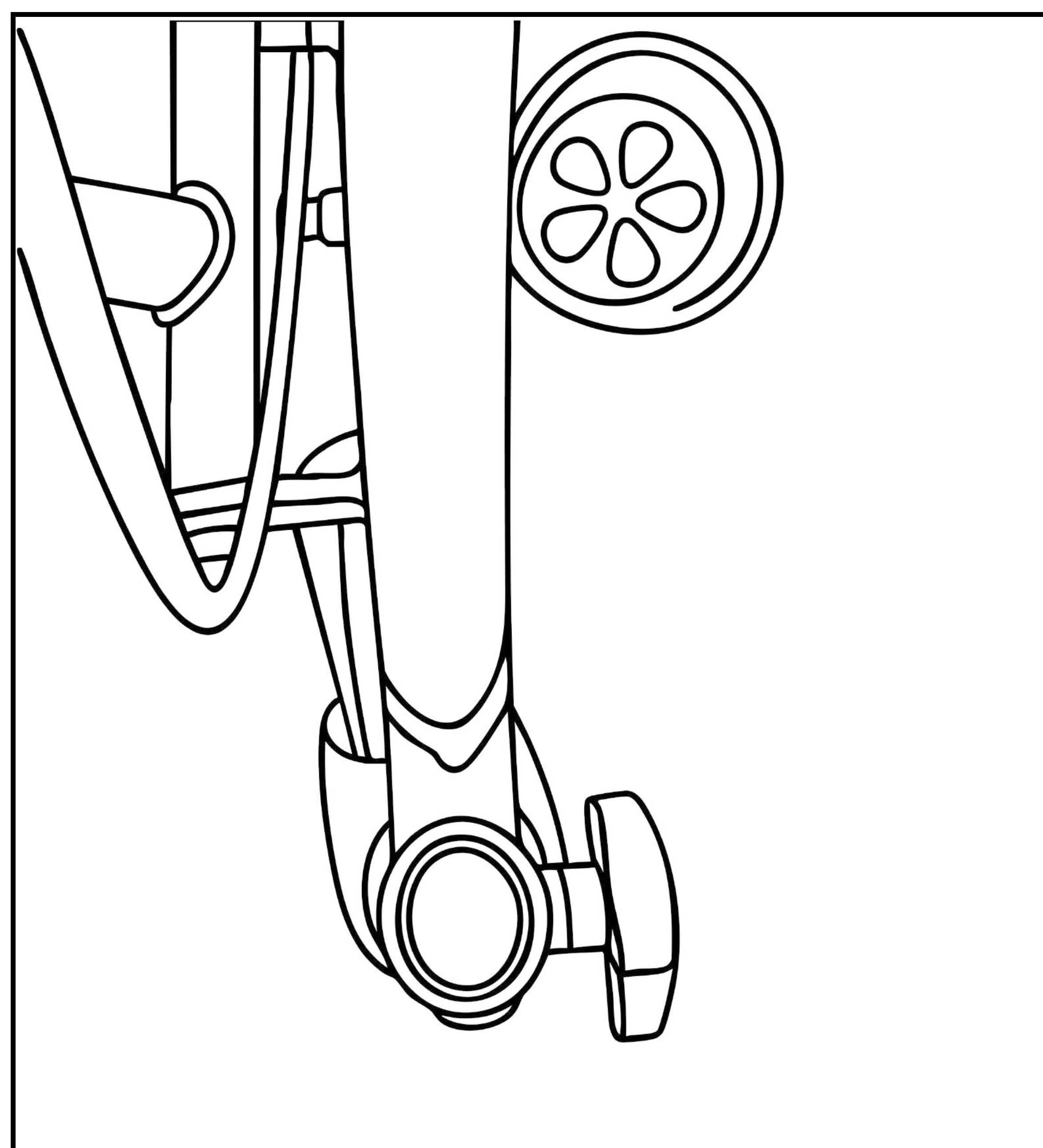Maintenance

Over time, parts may wear out, especially the brake cables and wheels. Have your rollator checked and serviced regularly by an authorized dealer.

Maintenance Schedule:

- **Weekly:**
 - Check for malfunctions
 - Check that wheels and swivel wheels are securely fastened
- **Monthly:**
 - Check hand grips for wear
 - Check that the frame is stable
 - Check brakes and brake levers
 - Check front and rear wheels
 - Check tires for wear or cracks
- **Annually:**

Have the rollator inspected annually by an authorized dealer for inspection, lubrication, and adjustment.

- **Parts Requiring Maintenance:**

- Wheels
- Hand grips
- Brakes
- Brake levers
- Bearings
- Knobs

Cleaning

Clean the rollator with a damp cloth and mild detergent.

- Do not use aggressive cleaners, abrasive agents, or high-pressure cleaners, as these may damage the coating.
- Dry the hand grips and seat thoroughly after cleaning.

Warning

1. Maximum user weight: 136 kg
2. Do not use as a seat while being pushed.
3. Ensure both hand grips are set to the same height.
4. Adjust the height only in consultation with a specialist.
5. Always keep a safe distance between feet and wheels.
6. Loose parts may cause damage or injury.
7. Replace damaged swivel wheels only with original parts.
8. Always fully unfold the rollator before use.
9. Keep fingers away from the folding mechanism when opening or closing.

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FOSHAN HCT MEDICAL EQUIPMENT CO.,LTD, No.501 of 19
SUNGO EUROPE B.V.,
Onebird B.V., Kazernestraat 17, 5928 NL, Venlo,