Zookie

Portable Electric Breast Pump

Instruction manual

Content

Product structure	04
Cleaning and disinfection	05
Operation cautions	05
Assemble	06
Start pumping	08
Charge	10
Instructions For Pouring Breast Milk	12
Storing Your Breast Milk	14

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

Product structure

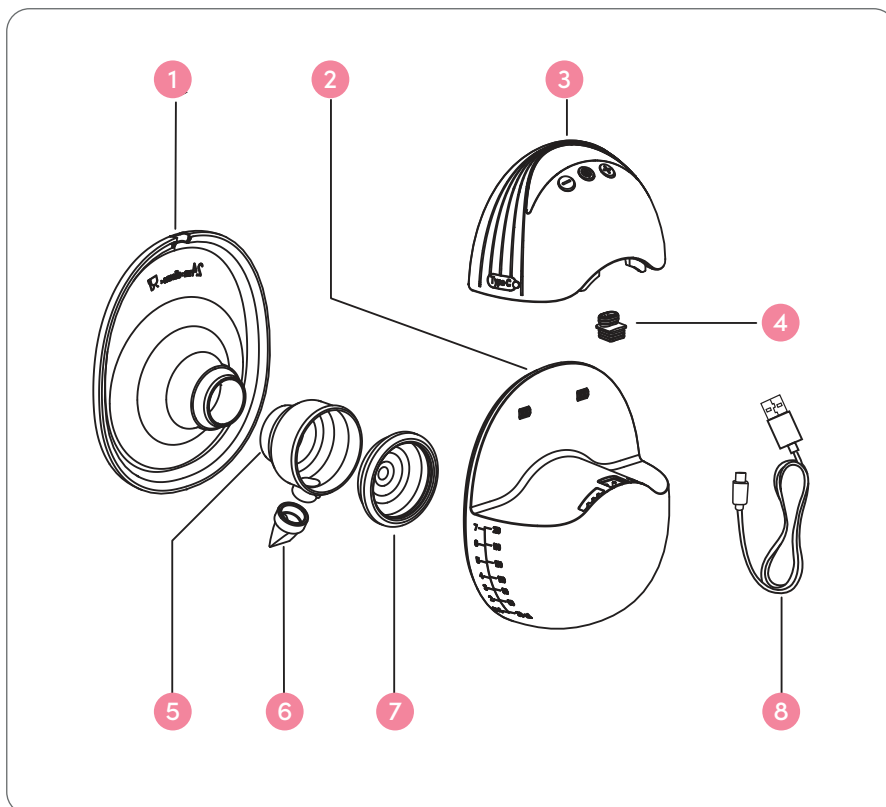Parts list

- | | |
|---------------------------|--------------------------------|
| 1 Siliconen shield | 5 T-junction |
| 2 Milk Collector | 6 Duckbill Valve |
| 3 Host | 7 Airbag |
| 4 Siliconen Plug | 8 Type-c Charging cable |

Control Panel

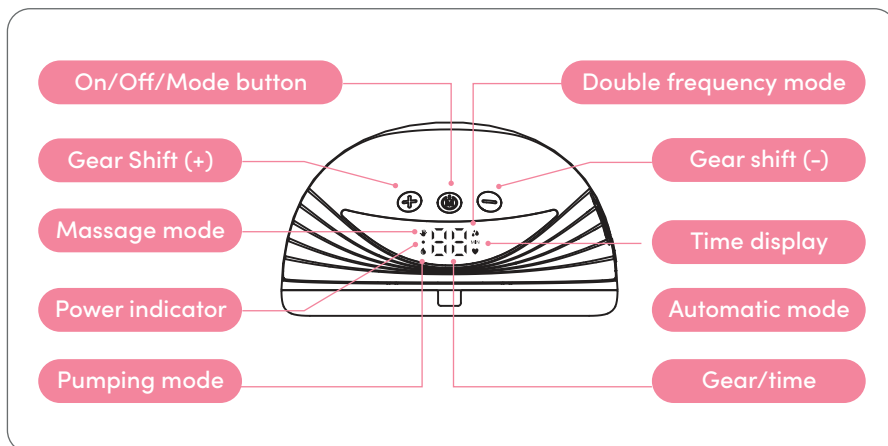Note: This product has a built-in lithium battery. If the battery indicator flashes, we recommend that you use a certified 5V=1A power supply to charge the main unit as quickly as possible.

Cleaning And Disinfection


1. Disassemble, wash and disinfect the parts that come into contact with breast milk (milk collector, silicone shield) before each use.
2. The duckbill valve should not be disassembled and washed separately to avoid loss and affecting normal use of the breast pump.
3. Disassemble and clean the parts that come into contact with breast and milk immediately after using to avoid breast milk residue and bacterial proliferation.
4. It is recommended to use drinking water for cleaning.

Operation Cautions


1. Never wash or disinfect the console, as the breast milk will not enter this component. You can wipe this part clean with a soft cloth.
2. Clean the silicone shield carefully, if it is damaged, your breast pump will not work properly.
3. Clean the milk collector. Do not rub it, do not insert anything. Failure to do so may cause damage. It is recommended to use drinking water for cleaning.
4. The duckbill valve does not need to be disassembled and cleaned separately to avoid loss and affect the normal use of the breast pump.

Assemble


1. Assemble the duckbill valve and airbag in place as shown in the figure, ensure that the duckbill valve is installed in the correct direction, and load the T-junction into the milk collector as shown in the figure. (Make sure the orientation is correct).

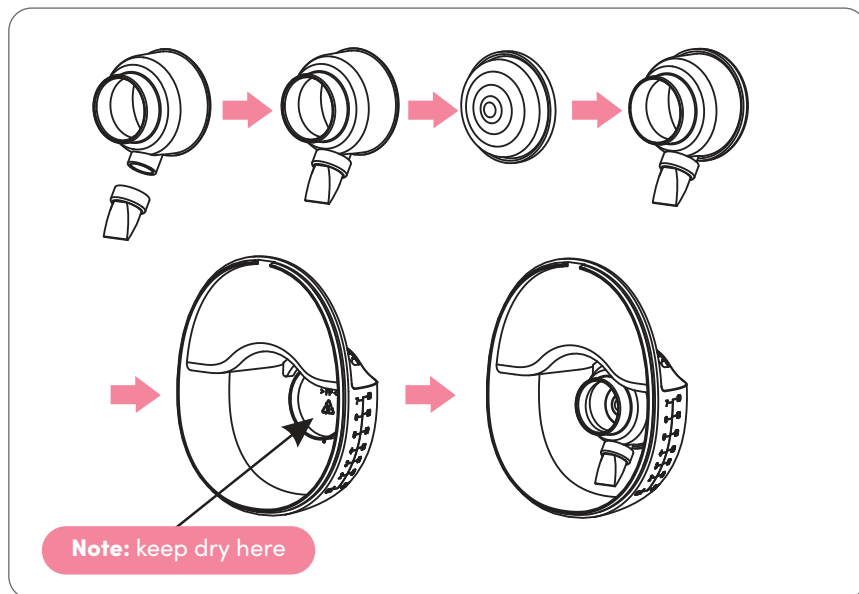
2. To install a silicone shroud on a milk collector, make sure the diversion groove on the silicone shield is aligned with the return groove on the milk collector and press the edges with both hands to ensure the seal is in place.

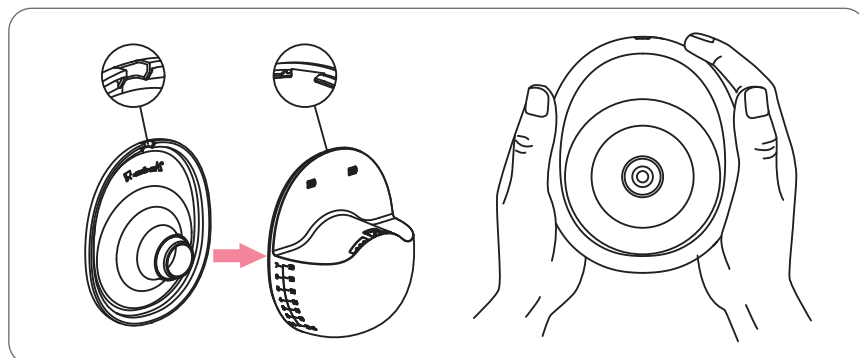3. Insert the silicon plug into the milk collector and press by hand to make the surface flush and ensure a perfect seal.

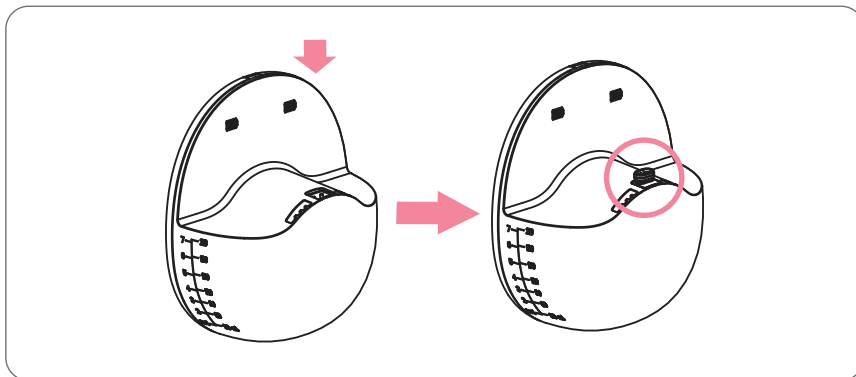
4. Mount the host on the milk collector, align the groove under the host with the milk collector, and fit the milk collector downwards to ensure a complete seal.

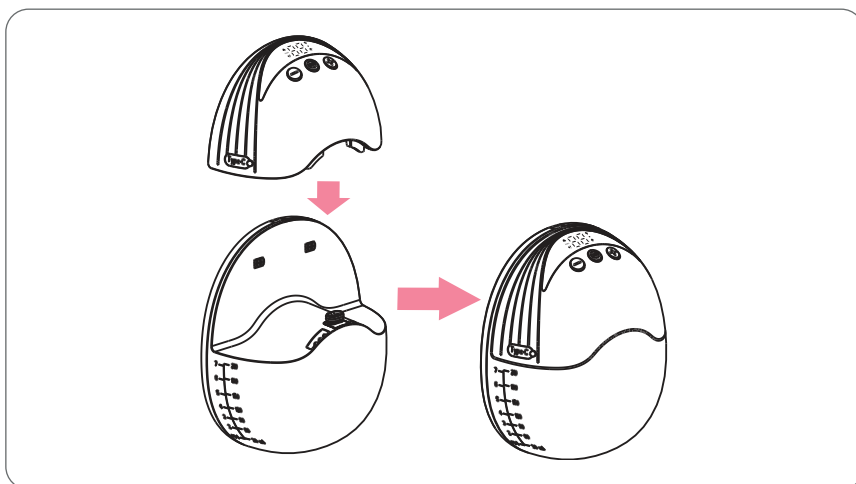
Notice: Power adapters are a common standard accessory on the market. This product is not provided; if necessary, it can be purchased separately. Power adapter specification: 5V==1A.

Start Pumping


1. Press the silicone shield firmly against the breast and make sure your nipple is centered. Stay upright and be careful not to twist the position of the product and chest, otherwise the breast pump will have no suction.

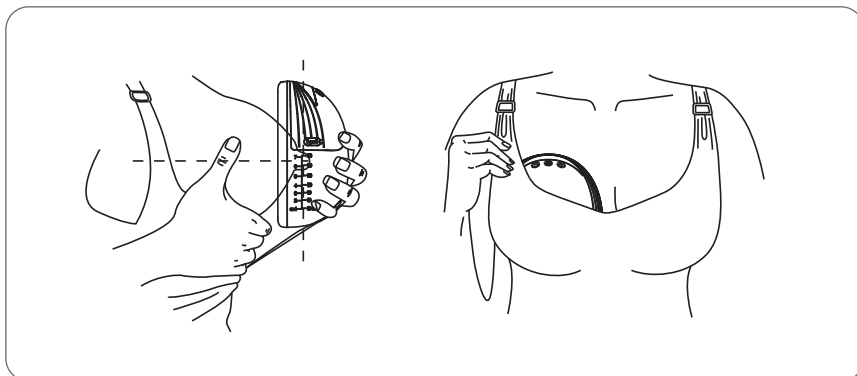
2. Long press  2S to turn on; in the power-on state, long press 2S to shut down. The breast pump will shut down automatically after 20 minutes of operation and the screen will stay off after 30s of operation. Any button can wake it up.
3. You can short press the  button to select the massage mode . This mode is faster than the breast pumping mode and has smaller suction power. During postpartum breastfeeding, the breasts are prone to swelling, pain, lumps and other breast discomfort. Using this mode can alleviate breast problems.

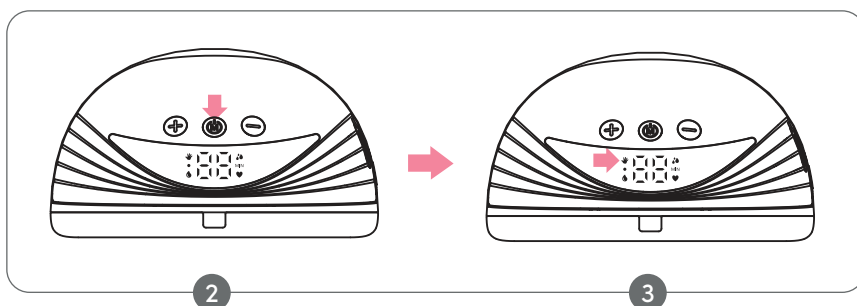
4. Short press the (M) button to select the pumping mode 🍼. This mode is slower than the massage mode, has stronger suction, and truly simulates the baby's sucking, more realistically stimulating the mother's smooth milk production.
5. You can short press the (M) key to select the double frequency mode 🍼. This mode is nine long and one short. In the early postpartum period, the milk ducts are not easy to open. You can stimulate the breasts and promote breast milk secretion by changing the mode.

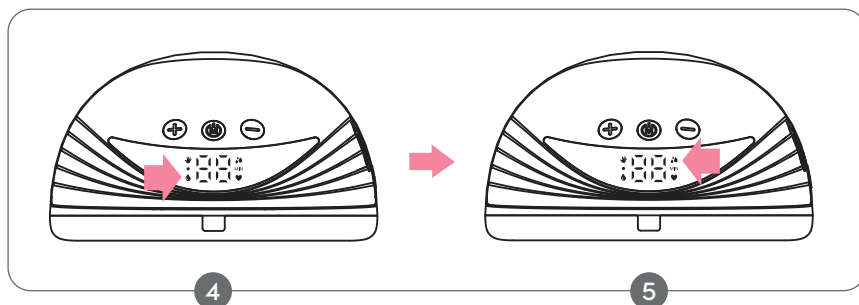
6. Short press the (M) key to select automatic mode ❤️. This mode is between the massage mode and the breast pumping mode. The suction is stable. The automatic mode will increase one gear every 15 seconds and finally stay at the 9th gear. During the upshift process, press the plus and minus keys to exit the upshift.
7. Gear/timing display: the default boot is 1 gear, and the display is "1"; A total of 12 gears. After the selected gear works for 5S without operation, it is transferred to the positive timing mode. The "MIN" timing mode light is on, and the timing is measured in "minutes".

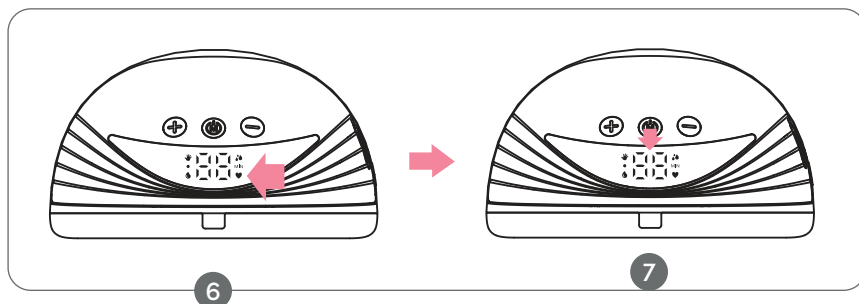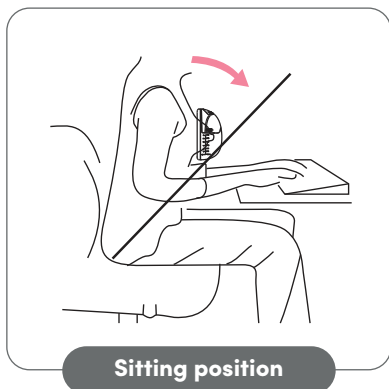Relax and sit in a comfortable chair with a slight lean forward (support your back with cushions). This allows the top of the silicone shield breast pump to tilt slightly upwards, which facilitates the flow of breast milk into the collector and reduces residual milk in the silicone shield.

8. After sucking, press the power button to turn off the breast pump, lean forward slightly, and carefully remove the breast pump from the breast. If there is a small amount of breast milk left over the silicone shield, gently pour the remaining breast milk from the silicone shield by tilting the breast pump.

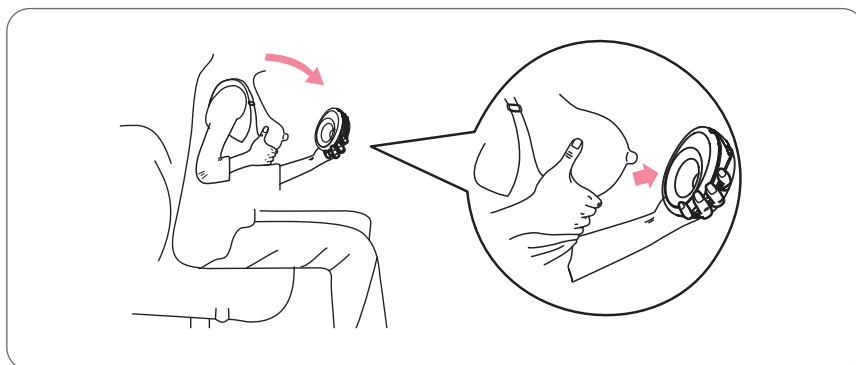
Charge

Battery indicator

1. The blue light flashes (when the voltage is 3.2V), indicating that it is urgent to charge; when the power is too low, the machine automatically shuts down (when the battery voltage is 3V).
2. It turns blue when it's charging and green when it's fully charged.
3. When the battery is charged, the machine will turn on with green light.




Open the USB cover to connect the USB, and connect the adapter at the other end to charge the motor.

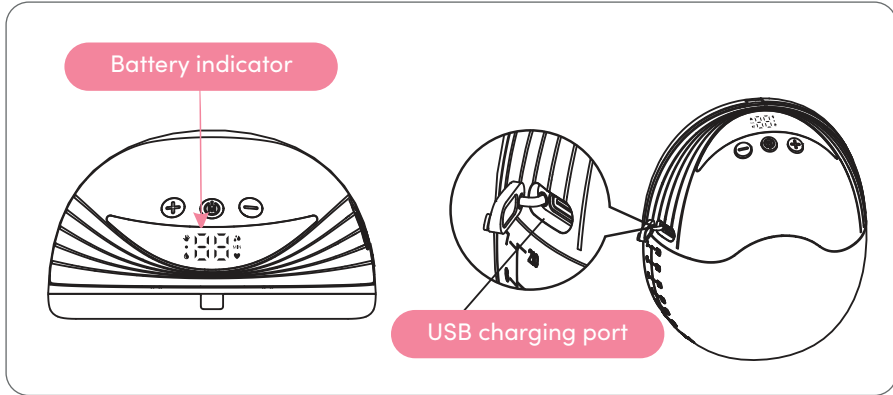**Notice:**

1. The device does not work when charging.
2. This product has a memory function; when you turn on again, the machine will automatically enter the working mode and suction level you used last time.
3. Automatic shutdown after working for 20 minutes after each start-up.
4. Continuous working hours vary depending on the battery or working conditions.
5. Charging time may vary depending on remaining power, temperature, or charging conditions.
6. USB power adapter is a common standard accessory in the market. This product is not available and can be purchased separately if required. Power adapter specification: 5V=1A.

Instructions For Pouring Breast Milk

After pumping, remove the breast pump from the breast to hold the breast pump upright, and pour the milk collector into the bottle, making sure the grooved part of the milk collector is aligned with the mouth of the bottle to prevent milk from spilling.

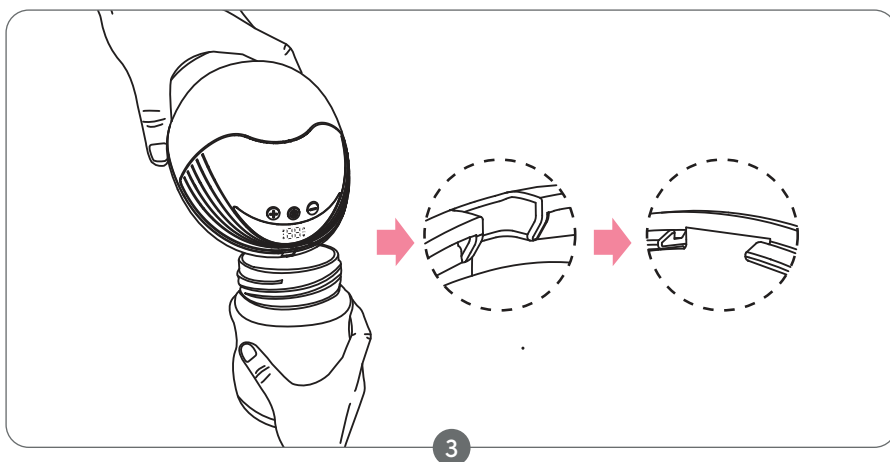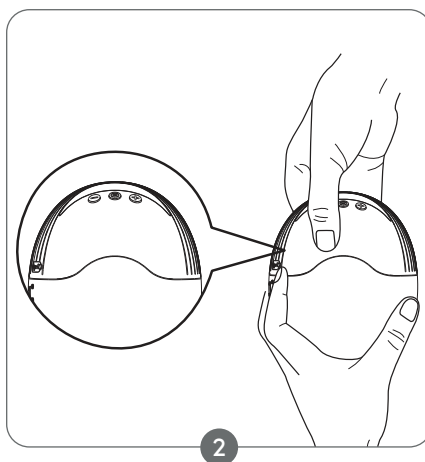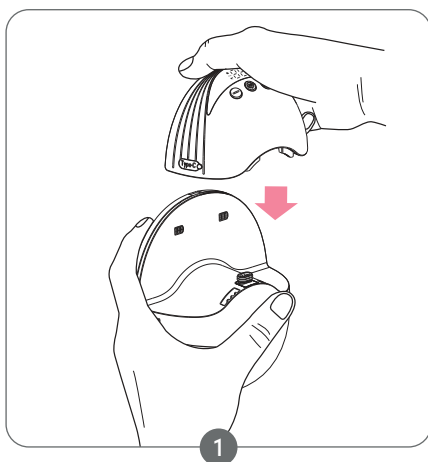Product parameters

Specifications

Product name	Wearable Electric Breast Pump
Product model	MY-375
Charging voltage	DC 5V == 1A
Charging current	730 mA
Battery type	Polymer battery
Battery capacity	1200 mAh
Rated voltage	3.7V
Rated current	1A Max
Working hours	>1.5H
Noise limit	40 mmHg ~ 300 mmHg
Extraction intensity limit	≤55 dB
Working environment	10–30°C, RH 15%–93%, 700–1060 hPa
Storage & transportation environment	-25–50°C, RH 15%–93%, 700–1061 hPa

Storing Your Breast Milk

Caution

1. Do not freeze breast milk repeatedly.
2. Do not mix fresh breast milk with frozen breast milk.
1. Be sure to discard the remaining breast milk after feeding.

Note: Only store breast milk collected from a sterilized breast pump.

Storing your breast milk

Expressed milk can be placed in a clean and sterilized storage bag or a bottle with a seal. Immediately refrigerate or freeze.

Caution

Expressed milk can be refrigerated for up to 48 hours (inside the refrigerator) or in the freezer for up to three months, with the bottle or storage bag clearly labelled with the date and time of expression. The first expressed milk should be consumed first.

Thawing of breast milk


Keep frozen breast milk in the freezer of your refrigerator overnight to defrost it. In case of emergency you can thaw breast milk in a bowl of hot water. After thawing, breast milk should be consumed within 24 hours. To heat breast milk, heat a bottle upright in a bowl of hot water or use a bottle and food warmer to heat it safely and quickly. Safe and fast heating of breast milk with a bottle and food warmer.

Warning

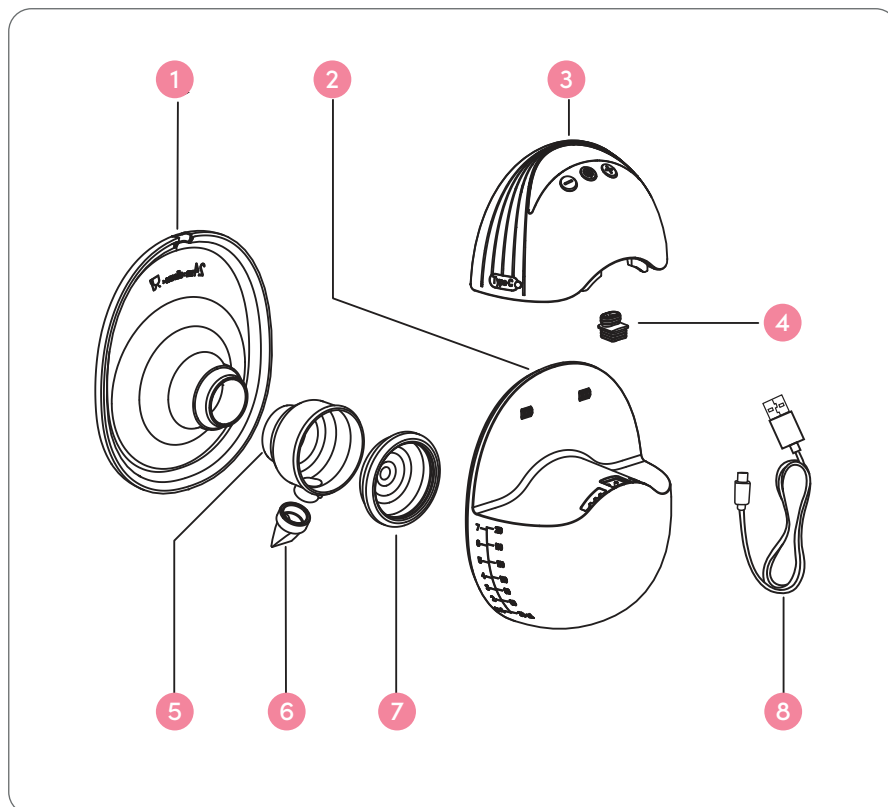1. Remember!! Always check the temperature of your breast milk before feeding your baby.
2. Do not microwave breast milk as this may destroy important nutrients and antibodies.
3. In addition, uneven heating may cause the liquid to overheat locally, which may cause burns.
4. Do not heat the milk storage bottle by immersing it in boiling water to prevent it from cracking.

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