




# Sonic Facial Device

User' s Manual

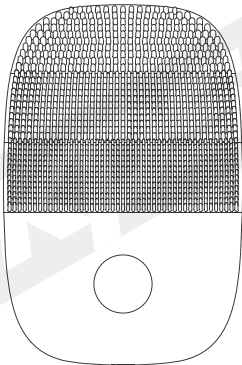

Born for love Living for Beauty

Welcome to use inFace

inFace Sonic Facial Device is independently developed and manufactured by Shenzhen inFace Technology Co. Ltd., which is affiliated to Shenzhen Migoo Smart Co. Ltd. Since the company was founded, it has been designing the core concepts of practicality, functionality and safety, aiming to create products with technology and intellectual, and striving to explore and innovate.

## **Elfin Cleanser** inFace Sonic Facial Device

### FEATURES

Open up a new cleansing frontier of your skin.

Cleaning — the first step of skin care.


Every face wash is a professional enjoyment of deep cleaning, gentle care and precise cleaning, laying a solid foundation for subsequent skin care.

1. The centrifugal motor controls the rotating speed (1200rpm to 4650rpm, and the cushioning design forms the tapping actions at the moment of buffering ), and the rotating effects skin directly through sonic vibration technology. Imitating the manual operation helps to absorb dirt from deep layer to skin surface, at the same time, avoid thick pores while cleaning.
2. Five-level cleansing massage mode to meet a variety of skin needs.
3. Silicone certified by the FDA (US Food and Drug Administration), and thickened to 4.5mm. The slope design of brush direction, 15 degree to 0 degree, fits the facial pore growth direction, which is targeted to clean, no dead angle, well care and do no harm to your skin surface .
4. International certified IPX7 waterproof (not waterproof when charging), it can withstand column watering and could be washed in the water directly.

## ACCESSORIES AND PRODUCT APPEARANCE

The following items should be found in the package:

- ① inFace Sonic Facial Device
- ② Charging Cable
- ③ User's Manual and Warranty Card

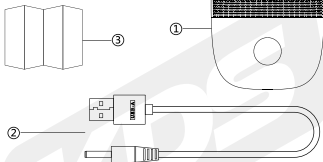

## PRODUCT SPECIFICATIONS

Product Model: MS2000 Pro

Product Size: 80mm X 53mm X 30.4mm

Net Weight: About 83g

Battery Capacity: 400mAh

Charging Voltage: 5V  $\overline{\text{---}}$

Charging Power: 1.5W

Charging Time: About 2.5h

Continuous Working Duration: About 240 minutes

Working Temperature Range: -10°C~40°C

Texture/Ingredient: FDA Food Grade Dow Corning Silicone

## OPERATION METHOD

1. In the off state, long press the button for 1.5 seconds to turn on the device. The product will automatically shut down if there is no operation within 90 seconds.
2. Short press the button to switch the 1-5 level.
3. In the on state, long press the key for 1.5 seconds to turn off the device.



## USING METHOD

1. Moisturize the face with foam cleanser and moisturize the silicone brush.
2. Long press the device for 1.5 seconds to turn on; short press the button to set the appropriate intensity.
3. Touch your face with the brush and move it slowly, stay 2-3 seconds along the contours of your face.
4. Wash off residual foam on your face and continue your daily facial care.




## CLEANING METHOD

1. Long press the button to turn off the device and clean it with clean water.




2. Shake off the water on the product, put it on the base (do not cover the transparent cover) and place in a ventilated dry place.

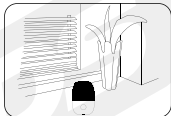

Notice: Do not use organic solvents, such as alcohol, benzine or thinners, to clean the product.

## CHARGING METHOD

1. When the battery is low, the white LED above the button will continue to flash.
2. Connect the USB adaptor to the charging interface at the bottom of the device to charge.
3. When charging, the white LED breathing; when full, the white LED is always on and no longer automatically extinguished; unplug the charging cable, the white LED is off. The device cannot be turned on during the whole charging process.

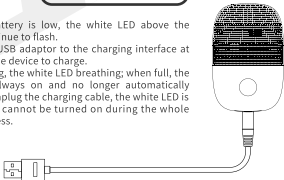

## PRECAUTIONS

1. Please read the instruction carefully before using the product.
2. The product vibrates while working. If you feel uncomfortable, please short press the button to adjust the intensity or stop using it.
3. Follow your doctor's advice for acne, wounds, or surgery, use with caution.
4. Sensitive skin is recommended to use the device with the weakest level.
5. Pregnant women should consult a physician before use.
6. The product is not recommended for long-term use in showers or baths.
7. The product is not recommended for long-term soaking in water to avoid affecting some properties.
8. Do not use while charging.
9. Avoid sun exposure so as not to affect the performance of the silicone and reduce the service life.
10. Do not use organic solvents, such as alcohol, benzene or thinners, to clean the product.
11. Turn off the power when cleaning the product.
12. In case of danger, do not expose children to the product.
13. Users who had the facial surgery/shaping are recommended to use the product under the guidance of a physician.
14. Clean the product in time after use.
15. For the first time to use the product, it is recommended to use up the battery, then fully charge it.
16. The final interpretation right of the instrument belongs to Shenzhen inFace Technology Co., Ltd..

## COMMON QUESTIONS

If you find any issues during use, please check the abnormality according to the following table. If you cannot solve the problems, contact the after-sales service.

| PHENOMENON   | SOLUTIONS   |
|--|---|
| During use, the product will automatically shut down if there is no operation within 90 seconds. | Software settings, Non-product quality issues, users can turn it on again.  |
| During use, the irritation caused by sonic vibration may make the skin uncomfortable.            | Please adjust the intensity appropriately to use it.  |
| What is the best frequency of using the product?   | Use once a day or every other day (depending on the skin condition to adjust the frequency); if you make up, it is recommended to use it every day. |

| PHENOMENON  | SOLUTIONS  |
|---|--|
| <p>What facial cleanser products are recommended to use with the product.</p> | <p>Foam-rich cleansing product is recommended. If using a non-squeezed foam product, use the non-activated device to gently rub the cleansing cream/powder/glue out of the rich foam, then turn on the device for use.</p> |

### PRESERVATION METHOD

- 1.If it is not used for a long time, please clean the product and store it in a dry and ventilated place.
- 2.Keep out of reach of children.

### ABOUT WARRANTY

- 1.The main body of the product (excluding accessories such as storage boxes and charging cables) enjoys one-year warranty from the date of purchase.
- 2.The return and change policy of inFace is in accordance with the service standard of "New Three Guarantees National Policy (China)" .
- 3.In case of the damage to the silicone material, please avoid the sun exposure.
- 4.If there is any quality problem after opening the box, please contact the after-sales customer service in time.

### LIMITATION OF WARRANTY

1. Failure to use or maintain the product according to the product manual;
2. Machines repaired by non-inFace authorized repair center;
3. Free gifts and outer packing are not subject to warranty service;
4. Damage caused by natural disasters and other force majeure;
5. Product failure or damage caused by accident or human factors (including operation error, incorrect insertion and extraction, scratch, disassembly, collision, improper voltage input, etc.).

If there is any change to the above list, the company's new provisions or the country's updated relevant laws and regulations shall prevail. Please consult the national service hotline 400-007-8388 for the specific warranty terms.



## WARRANTY CARD

Product Number: \_\_\_\_\_

Product Name: \_\_\_\_\_

User Name: \_\_\_\_\_

Phone Number: \_\_\_\_\_

User Address: \_\_\_\_\_

Dealer Information: \_\_\_\_\_

Date of Purchase: \_\_\_\_\_

Order Number: \_\_\_\_\_

Trouble Description: \_\_\_\_\_

\_\_\_\_\_

\_\_\_\_\_

(Fill in the warranty card details when applying for the repairing service. The warranty period is one year for the whole product.)



inFace Sonic Facial  
Device

Product Name: \_\_\_\_\_

Product Model: MS2000 Pro

Production Date: \_\_\_\_\_

Quality Inspection: \_\_\_\_\_

**CERTIFICATION**

Shenzhen Migoo Smart Co., Ltd.





Please Scan the QR code to  
re-view the User Manual

### AFTER-SALES CONTACT

Phone: 0086-4000078388

Service WeChat: migoosmart

Working Hours: 9:00-18:00, Monday to Saturday  
(China time, Day off in Holidays and  
weekend)

Company: Shenzhen Migoo Smart Co., Ltd.

Address: Rm.807, Building U6, Junxiang U8 Intelligent  
Manufacturing Industrial Park, Gushu  
Community, Xixiang Sub-District, Bao'an  
District, Shenzhen, China

E-mail: [service@migoosmart.com](mailto:service@migoosmart.com)

Post Code: 518000



WeChat  
Official Account



Service WeChat

※ The final interpretation right belongs to  
Shenzhen Migoo Smart Co., Ltd.